

SCOTTSDALE POLICE DEPARTMENT CRIME LABORATORY BLOOD ALCOHOL FACE SHEET

ANALYSIS DATE Aug 6, 2026

SEQUENCE NAME 06Aug26

EQUIPMENT

Pipettor ☐ Hamilton ML600GA23163 ☒ Hamilton ML600EA20851
Gas Chromatograph ☐ Agilent US14173023 ☒ Agilent US2320A024

INSTRUMENT CALIBRATION

Vial 1 0.02 calibrator Lot FN03062417 Coefficient of determination (r^2)
Vial 2 0.10 calibrator Lot FN04042408 0.99999
Vial 3 0.20 calibrator Lot FN07012401
Vial 4 0.40 calibrator Lot FN04092401

CALIBRATION VERIFICATION AND RESOLUTION TEST


Vial	Sample	Expected result	Measured result	Manufacturer/lot
5	Blank	Not detected	Not detected	SPD lab 061026WB
6	Volatiles mixture	6 compounds	6 compounds	SPD lab 031726MIX
7	Aqueous control	0.400 g/dL	0.403 g/dL	Lipomed C1475825
8	Aqueous control	0.040 g/dL	0.040 g/dL	Lipomed 529862
9	Blood control	0.203 d/dL	0.203 g/dL	ACQ 4190425129/3
20	Aqueous control	0.100 g/dL	0.101 g/dL	Lipomed C1455120
31	Blood control	0.203 g/dL	0.203 g/dL	ACQ 4190425129/3
42	Aqueous control	0.100 g/dL	0.101 g/dL	Lipomed C1455120
53	Aqueous control	0.400 g/dL	0.406 g/dL	Lipomed C1475825
54	Aqueous control	0.040 g/dL	0.041 g/dL	Lipomed 529862
55	Blood control	0.203 g/dL	0.204 g/dL	ACQ 4190425129/3
56	Blank	Not detected	Not detected	SPD lab 060326AQ

SUBJECT SAMPLES

Subjects in the sequence 20 Subjects requiring reanalysis 0

ADDITIONAL NOTES: All testing proceeded as expected.

Run valid ☒
Run invalid ☐

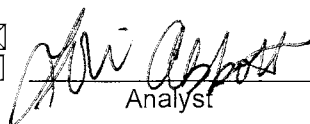
Analyst

Run valid ☒
Run invalid ☐


Technical Reviewer

Scottsdale Police Department Crime Laboratory


Sequence Quality Assurance Summary

Sequence: 06Aug26

Analyst: Abbott

Sample Name	Vial	Measured Value (g/dL)	Expected Value (g/dL)	Percent Difference	Absolute Difference (g/dL)
blank 061026WB	5	negative	negative	-	-
0.400 C1475825	7	0.403	0.400	0.75	0.003
0.040 529862	8	0.040	0.040	0.00	0.000
0.203 4190425129/3	9	0.203	0.203	0.00	0.000
0.100 C1455120	20	0.101	0.100	1.00	0.001
0.203 4190425129/3	31	0.203	0.203	0.00	0.000
0.100 C1455120	42	0.101	0.100	1.00	0.001
0.400 C1475825	53	0.406	0.400	1.50	0.006
0.040 529862	54	0.041	0.040	2.50	0.001
0.203 4190425129/3	55	0.204	0.203	0.49	0.001
blank 060326AQ	56	negative	negative	-	-

Calibrator	Ethanol Area	N-propanol Area	Ratio
0.020	13.731	196.677	0.070
0.100	70.803	195.367	0.362
0.200	141.368	193.355	0.731
0.400	284.148	192.824	1.474

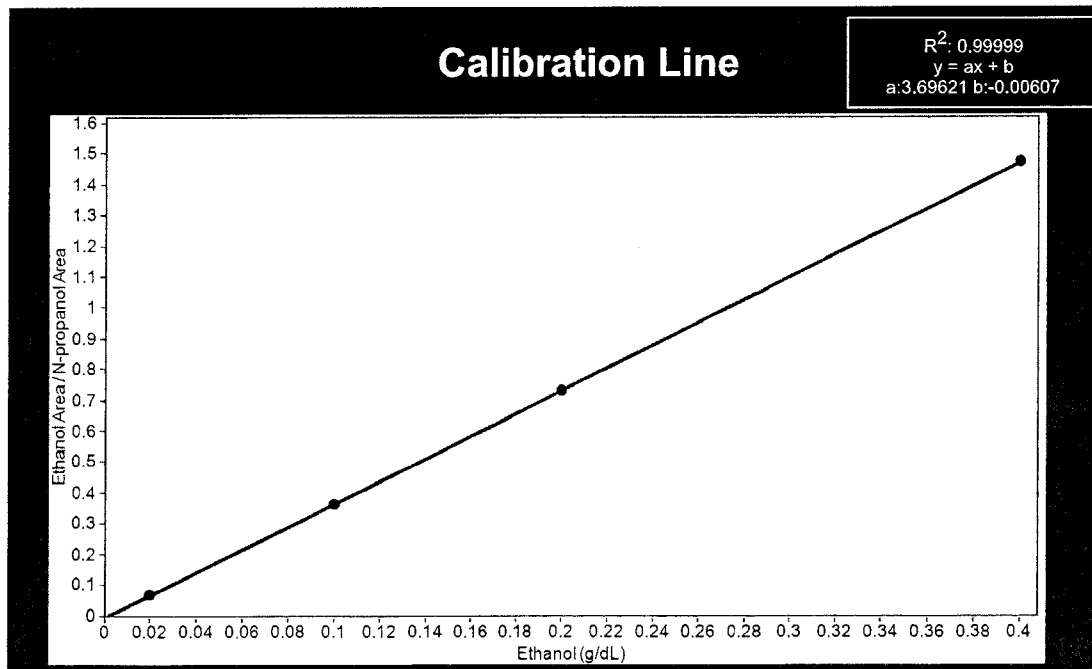

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	FN03062417	Vial:	1
Description:	0.020 Calibrator	Sequence:	06Aug26
Method:	ethanol quant.M	Injection date:	8/6/2026 2:33:08 PM
Instrument:	US2320A024 US23180007	Analyst:	Abbott 
Data file:	C:\Chem32\1\Data\06Aug26\1.D		

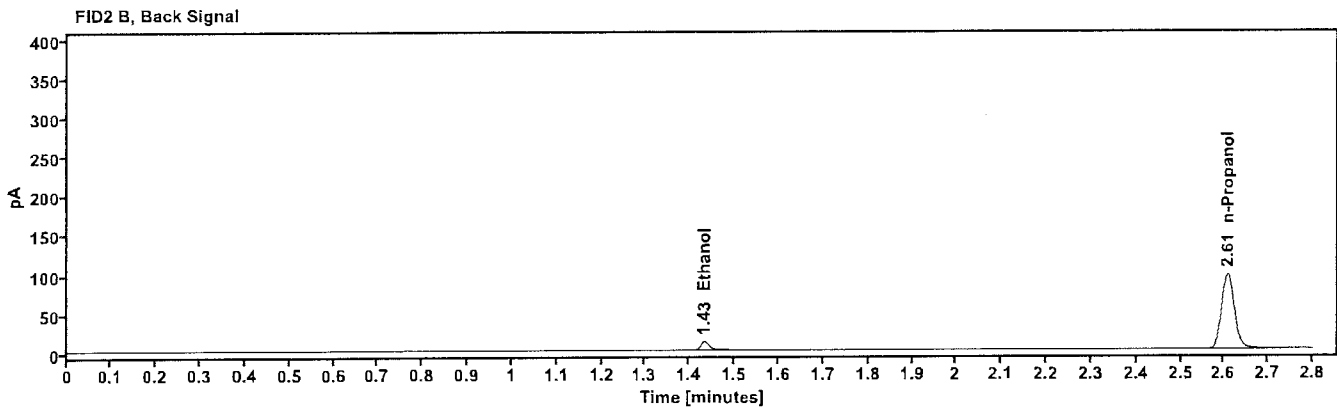
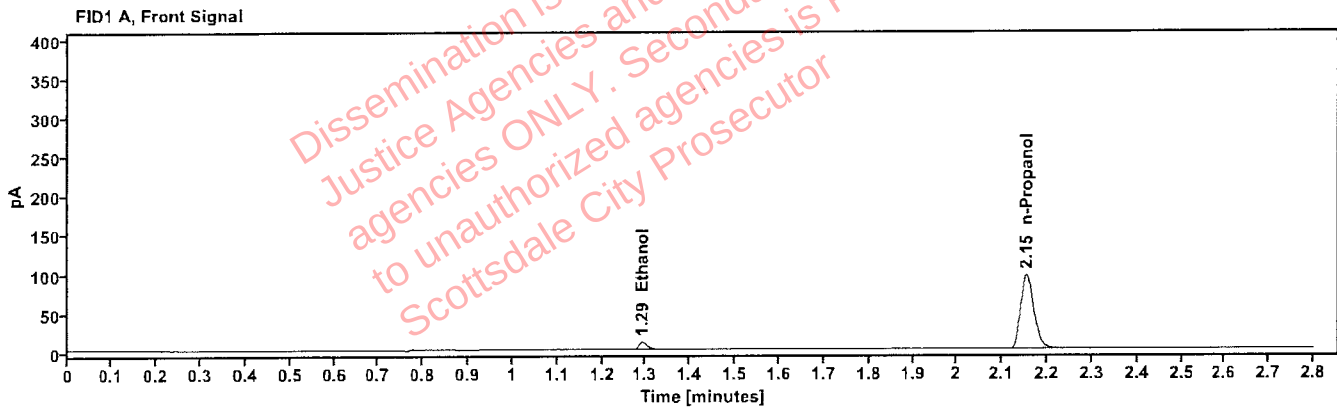


Table 1: FID 1 A (column DB-ALC1)

Compound	Time (min)	Peak Area
Ethanol	1.292	13.731
n-Propanol	2.154	196.677

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.433	14.328
n-Propanol	2.609	201.113

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	FN04042408	Vial:	2
Description:	0.100 Calibrator	Sequence:	06Aug26
Method:	ethanol quant.M	Injection date:	8/6/2026 2:36:55 PM
Instrument:	US2320A024 US23180007	Analyst:	Abbott 
Data file:	C:\Chem32\1\Data\06Aug26\2.D		

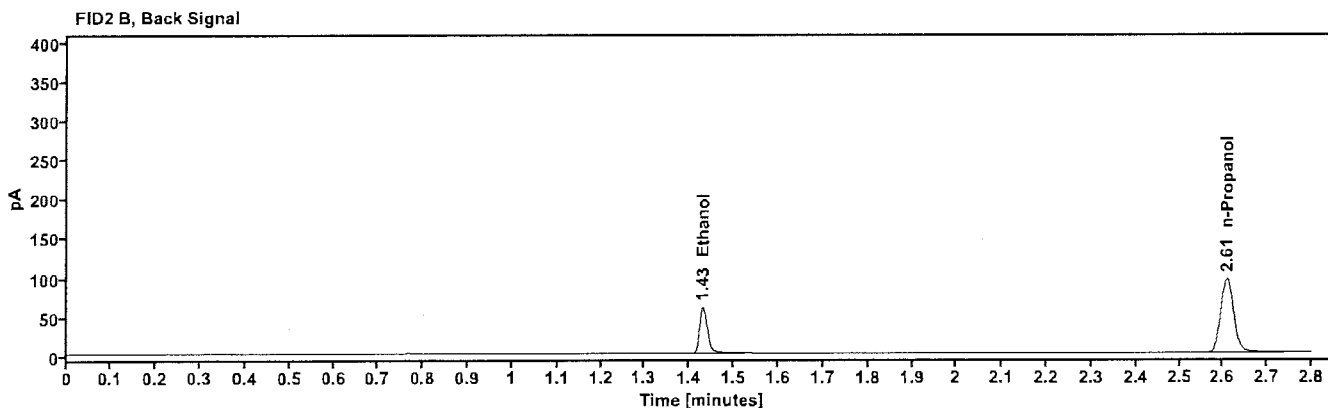
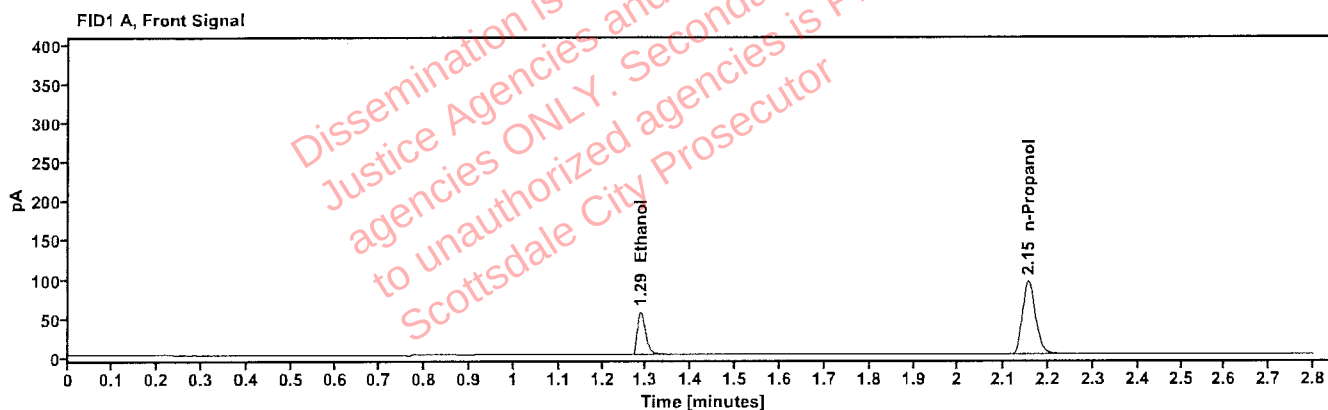


Table 1: FID 1 A (column DB-ALC1)

Compound	Time (min)	Peak Area
Ethanol	1.287	70.803
n-Propanol	2.154	195.367

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.430	73.274
n-Propanol	2.609	199.470

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	FN07012401	Vial:	3
Description:	0.200 Calibrator	Sequence:	06Aug26
Method:	ethanol quant.M	Injection date:	8/6/2026 2:40:43 PM
Instrument:	US2320A024 US23180007	Analyst:	Abbott <i>ph</i>
Data file:	C:\Chem32\1\Data\06Aug26\3.D		

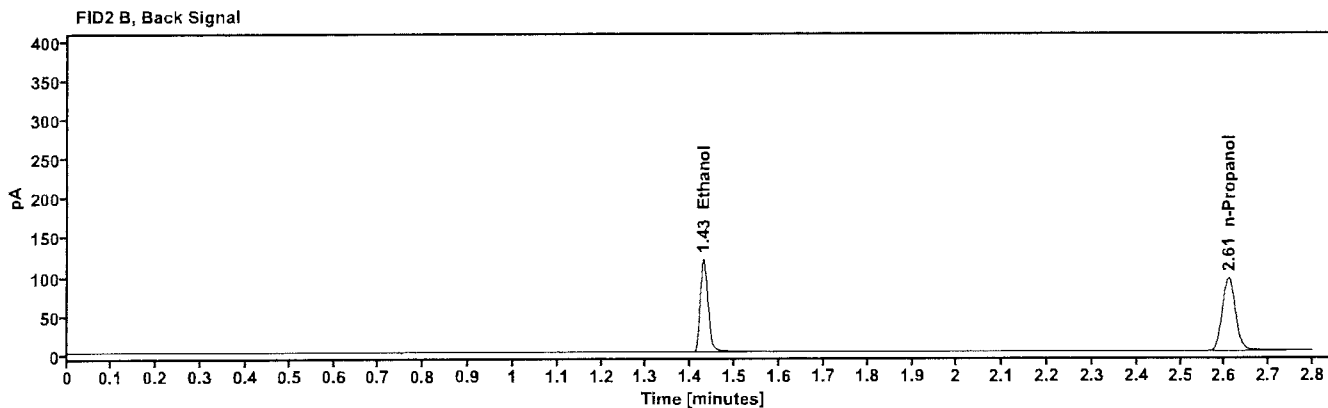
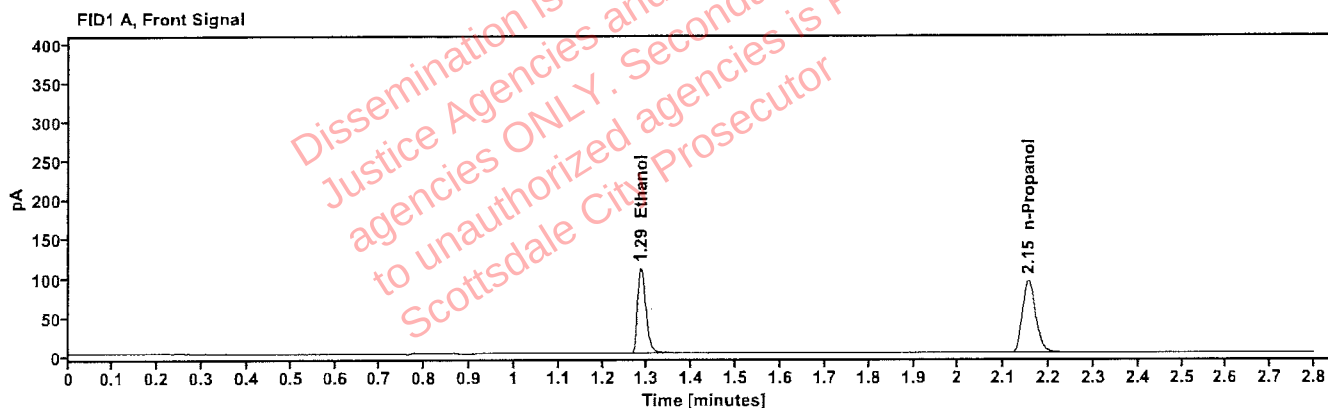


Table 1: FID 1 A (column DB-ALC1)

Compound	Time (min)	Peak Area
Ethanol	1.286	141.368
n-Propanol	2.154	193.355

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.429	146.075
n-Propanol	2.608	197.307

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	FN04092401	Vial:	4
Description:	0.400 Calibrator	Sequence:	06Aug26
Method:	ethanol quant.M	Injection date:	8/6/2026 2:44:32 PM
Instrument:	US2320A024 US23180007	Analyst:	Abbott <i>du</i>
Data file:	C:\Chem32\1\Data\06Aug26\4.D		

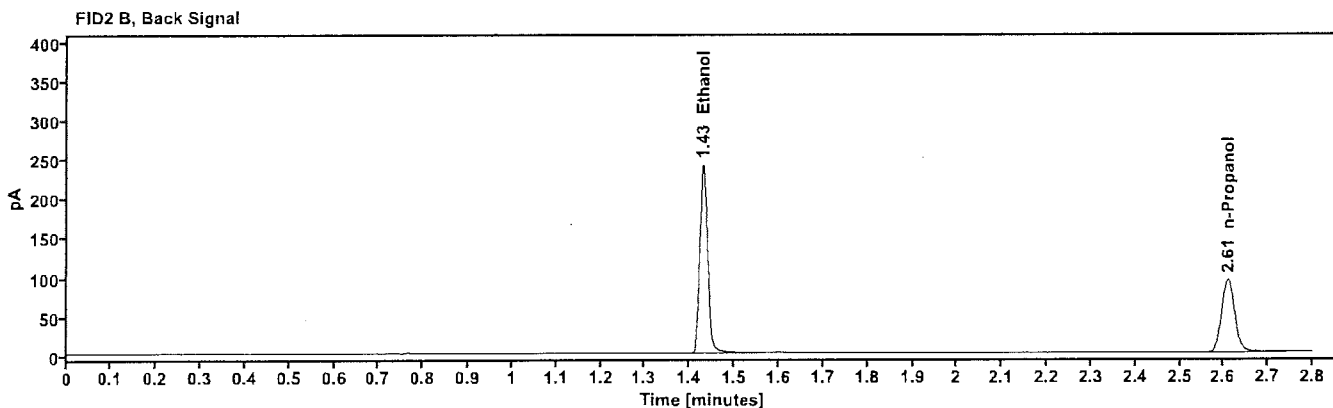
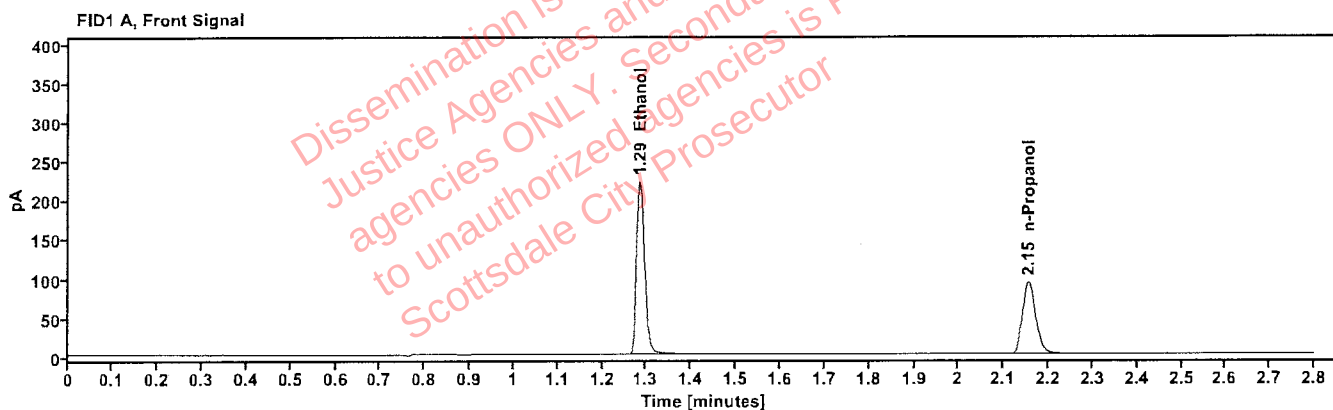


Table 1: FID 1 A (column DB-ALC1)

Compound	Time (min)	Peak Area
Ethanol	1.285	284.148
n-Propanol	2.154	192.824

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.429	293.081
n-Propanol	2.609	196.706

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	061026WB	Vial:	5
Description:	negative	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 2:48:19 PM
Data file:	C:\Chem32\1\Data\06Aug26\5.D	Analyst:	Abbott 

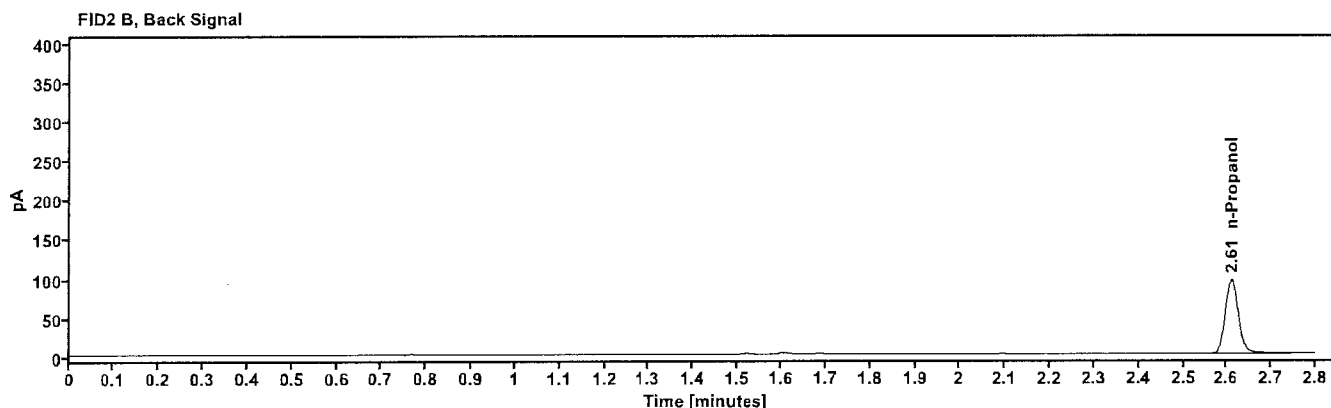
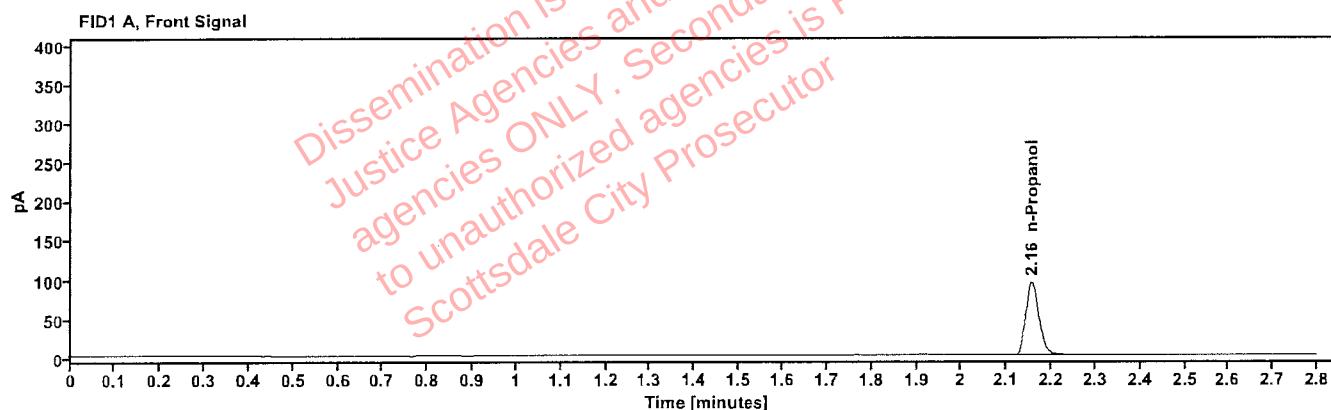


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
n-Propanol	-----	2.155	194.842

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
n-Propanol	2.610	199.390

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	031726MIX	Vial:	6
Description:	Volatile Mix	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 2:52:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\6.D	Analyst:	Abbott <i>du</i>

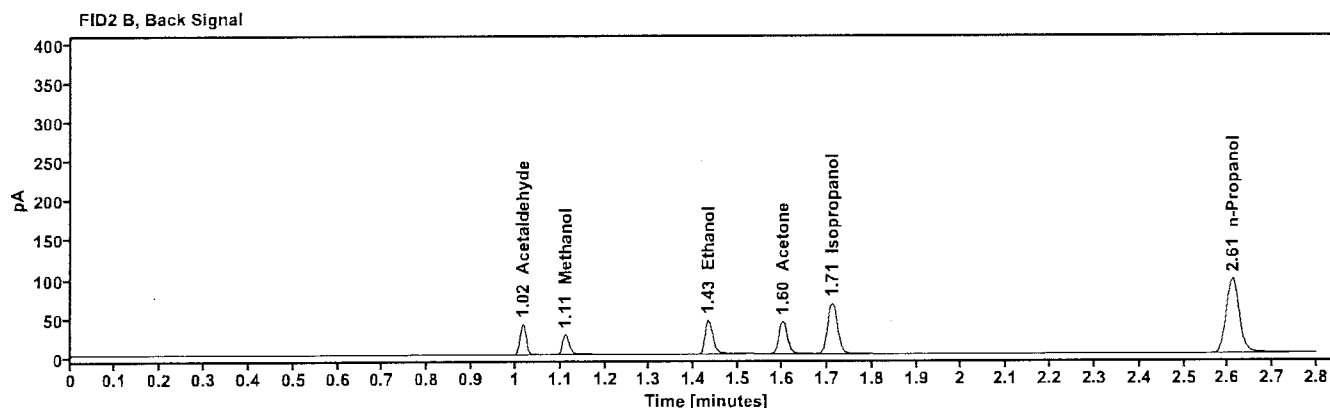
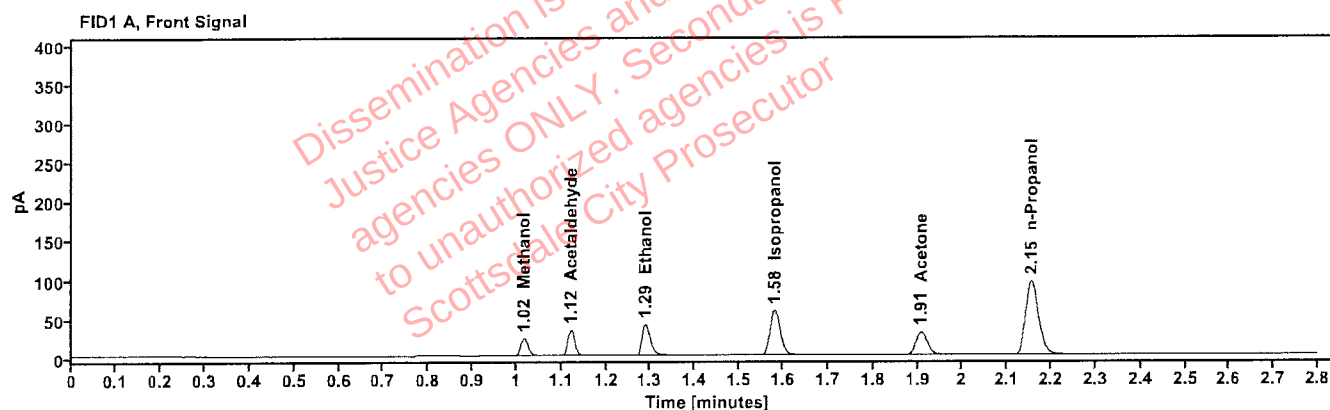


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
Methanol	-----	1.017	26.689
Acetaldehyde	-----	1.122	36.171
quantified Ethanol	0.0741	1.289	52.772
Isopropanol	-----	1.581	95.962
Acetone	-----	1.908	54.126
n-Propanol	-----	2.155	196.992

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Acetaldehyde	1.017	37.596
Methanol	1.111	27.368
Ethanol	1.432	54.574
Acetone	1.601	55.970
Isopropanol	1.711	98.885
n-Propanol	2.609	200.743

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	C1475825	Vial:	7
Description:	0.400	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 2:57:07 PM
Data file:	C:\Chem32\1\Data\06Aug26\7.D	Analyst:	Abbott <i>dh</i>

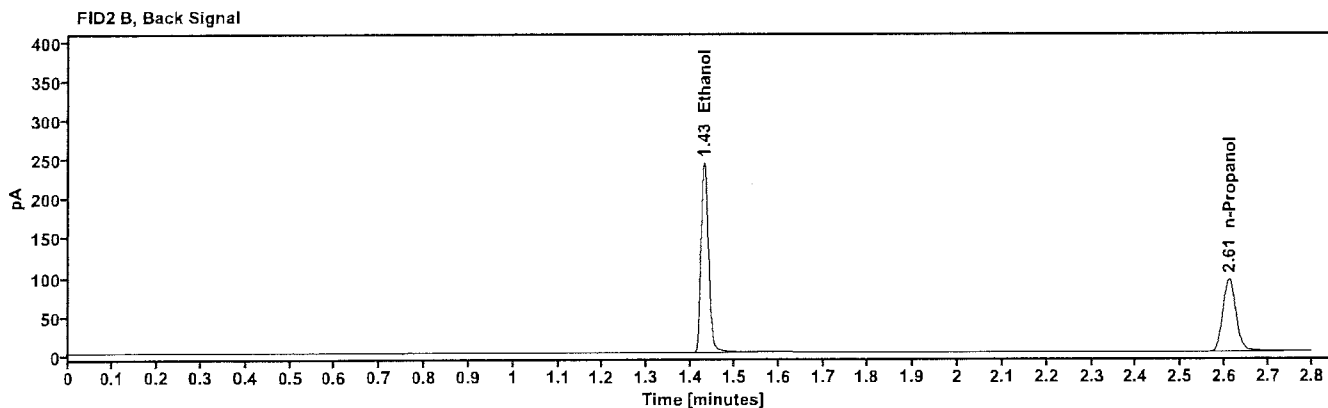
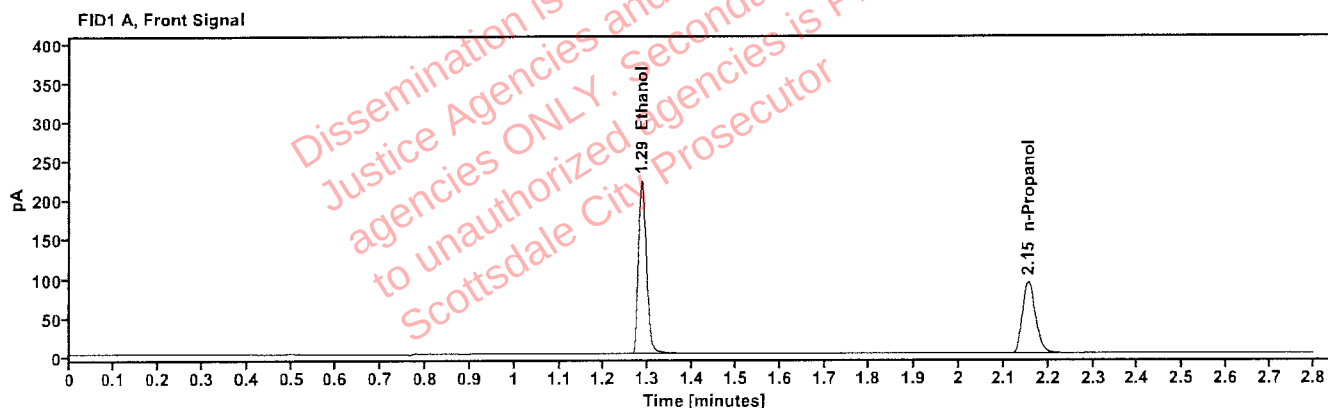


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.4039	1.285	284.711
n-Propanol	-----	2.155	191.477

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.430	293.489
n-Propanol	2.610	195.239

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	529862	Vial:	8
Description:	0.040	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:00:56 PM
Data file:	C:\Chem32\1\Data\06Aug26\8.D	Analyst:	Abbott <i>DA</i>

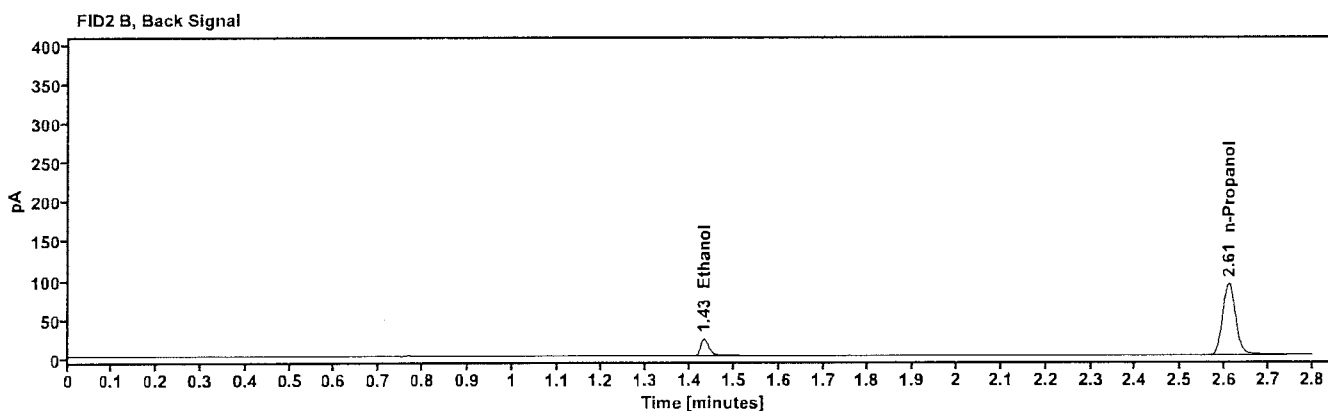
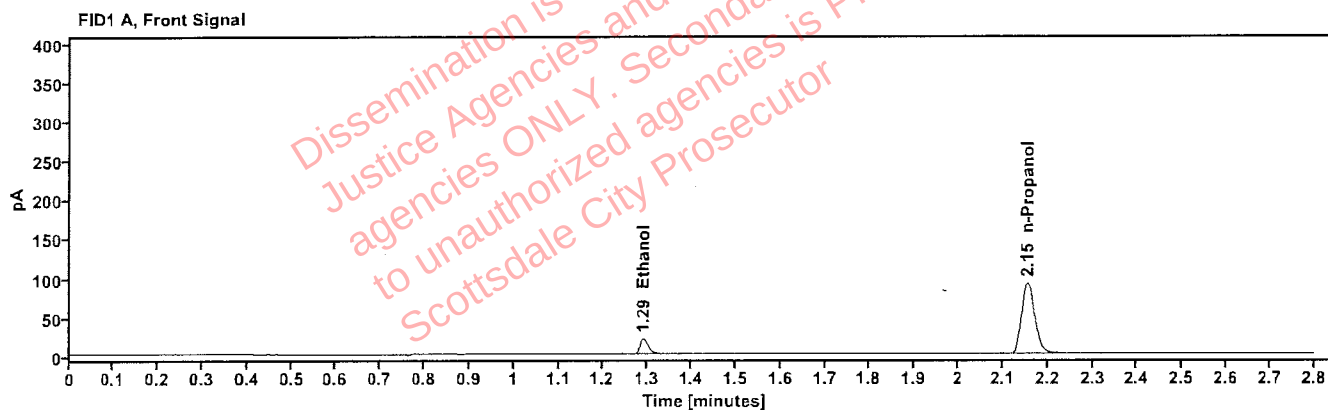


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0401	1.290	26.936
n-Propanol	—	2.155	189.302

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	27.935
n-Propanol	2.609	193.306

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	4190425129/3	Vial:	9
Description:	0.203	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:04:44 PM
Data file:	C:\Chem32\1\Data\06Aug26\9.D	Analyst:	Abbott

Abbott

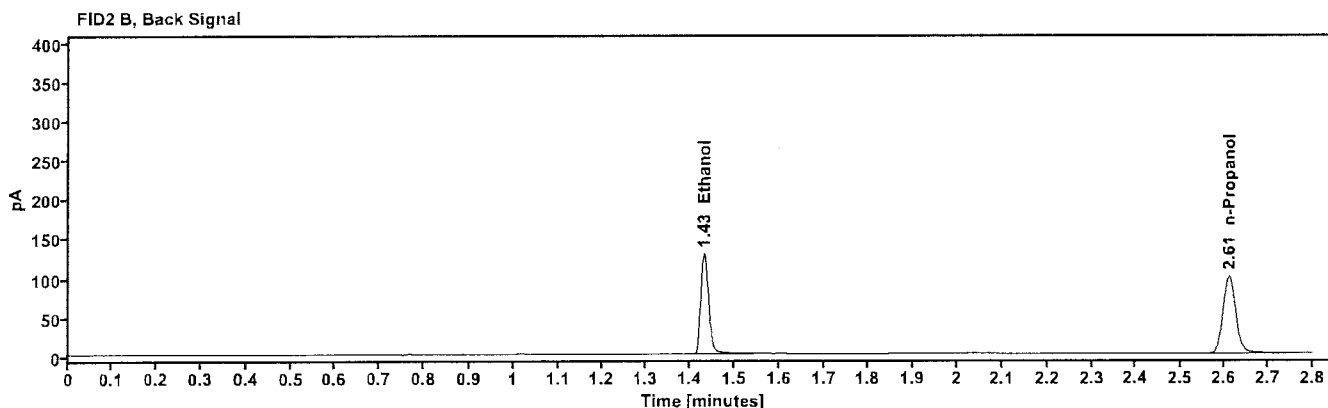
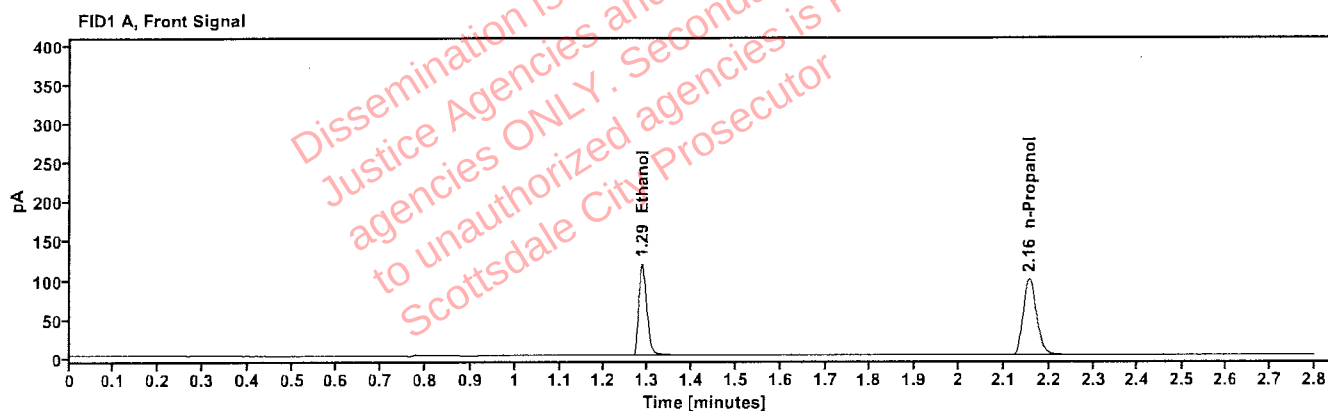


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2032	1.287	151.850
n-Propanol	----	2.156	203.825

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	156.820
n-Propanol	2.611	207.854

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	C1455120	Vial:	20
Description:	0.100	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:48:56 PM
Data file:	C:\Chem32\1\Data\06Aug26\20.D	Analyst:	Abbott <i>ah</i>

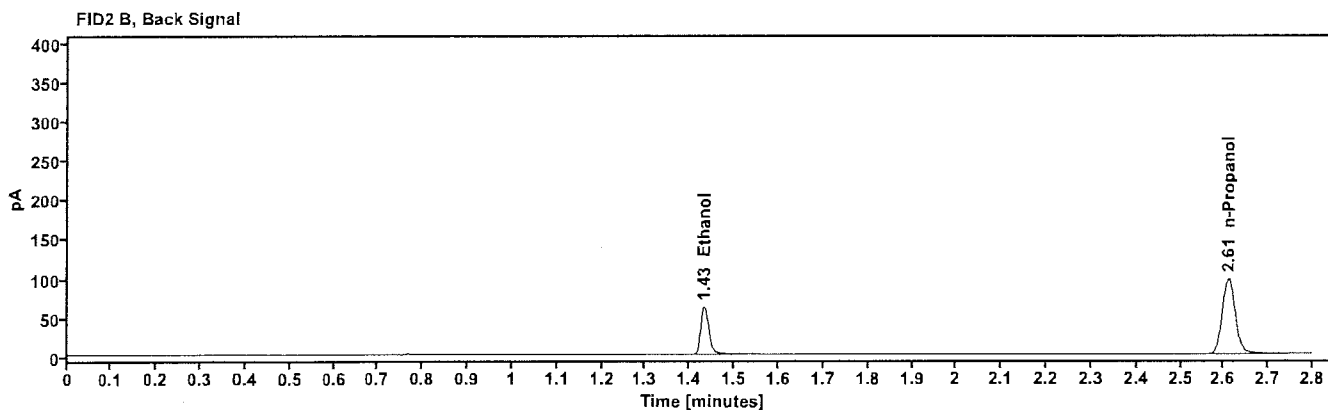
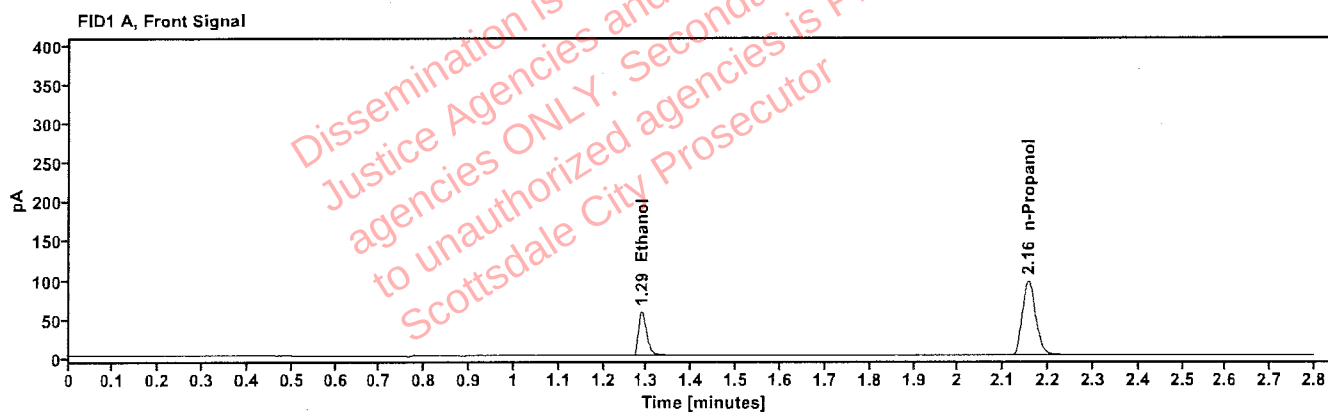


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1016	1.288	73.673
n-Propanol	-----	2.155	199.335

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	75.946
n-Propanol	2.610	203.300

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	4190425129/3	Vial:	31
Description:	0.203	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:33:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\31.D	Analyst:	Abbott

sh

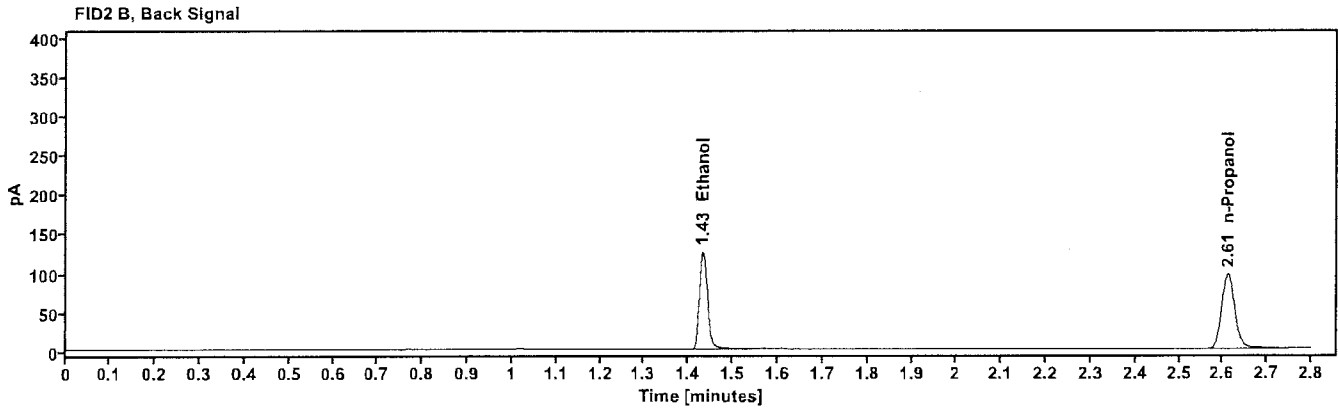
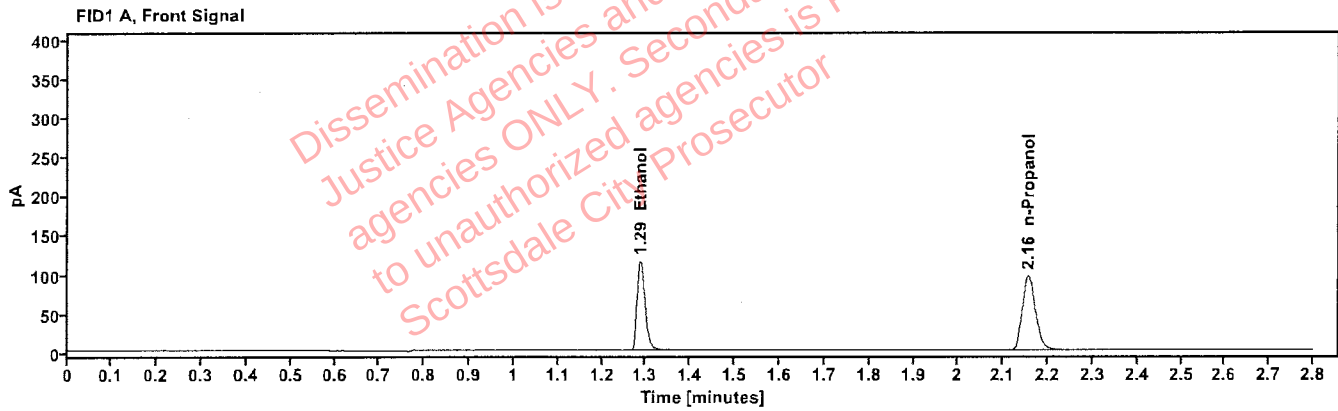


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2036	1.287	147.929
n-Propanol	-----	2.156	198.197

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	152.783
n-Propanol	2.611	202.046

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	C1455120	Vial:	42
Description:	0.100	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:16:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\42.D	Analyst:	Abbott

[Signature]

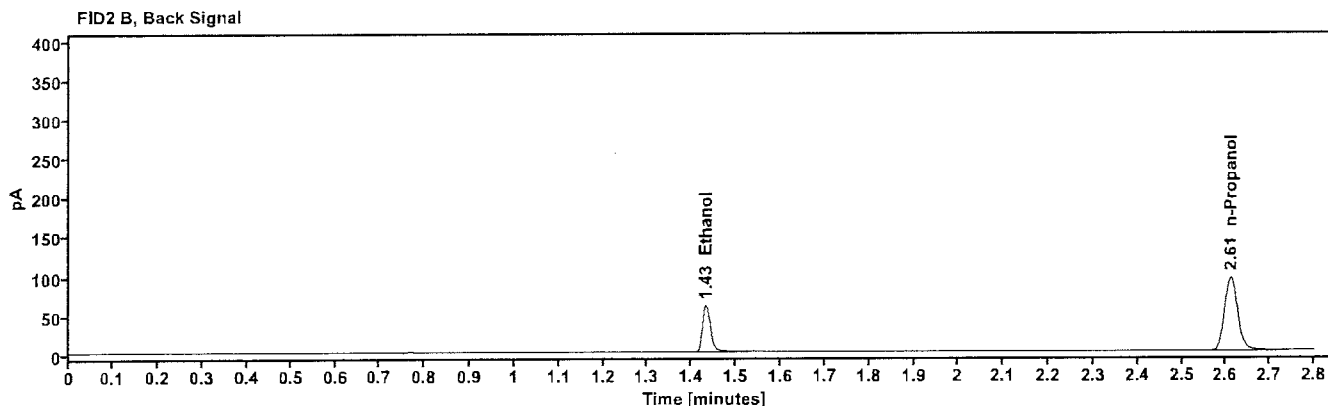
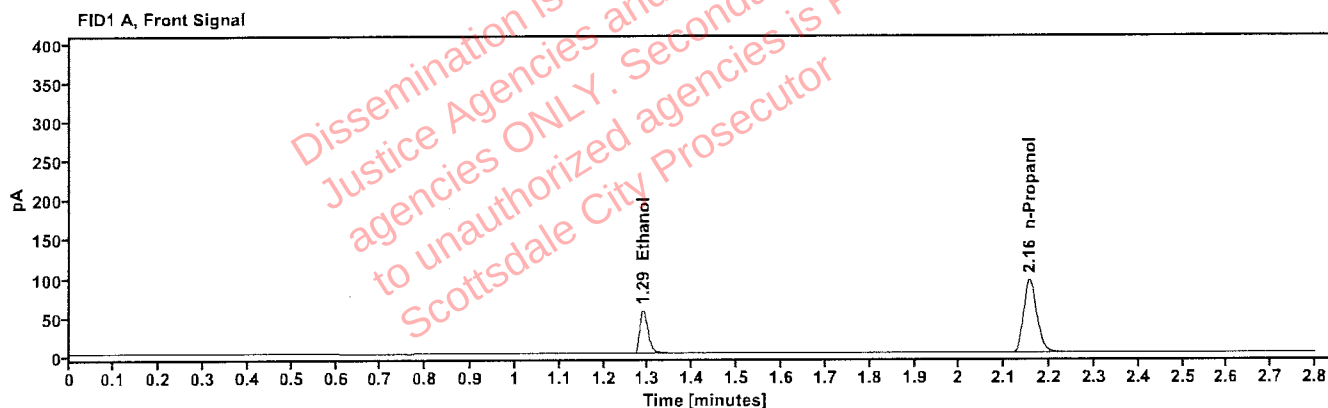


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1017	1.288	72.261
n-Propanol	----	2.156	195.404

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	74.491
n-Propanol	2.611	199.030

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	C1475825	Vial:	53
Description:	0.400	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 6:00:19 PM
Data file:	C:\Chem32\1\Data\06Aug26\53.D	Analyst:	Abbott

Abbott

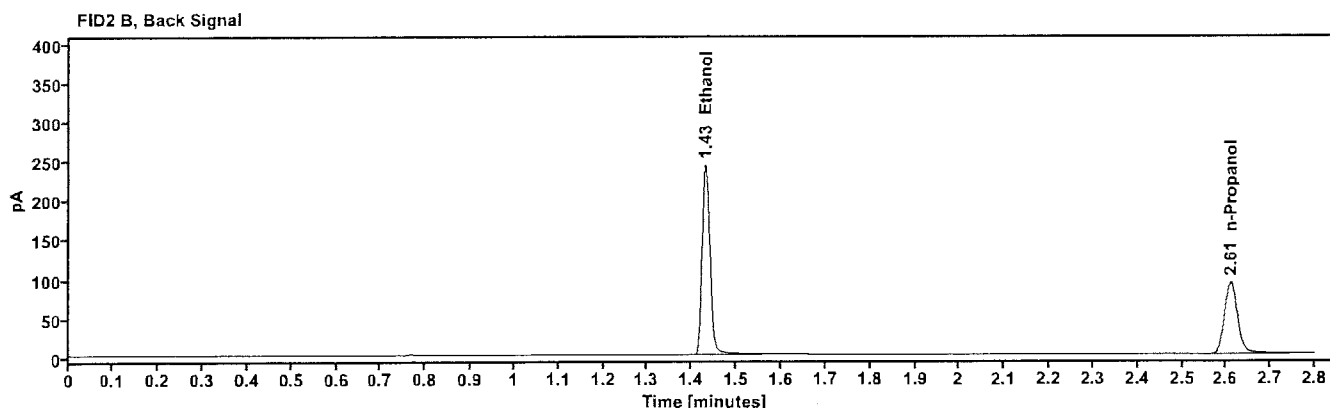
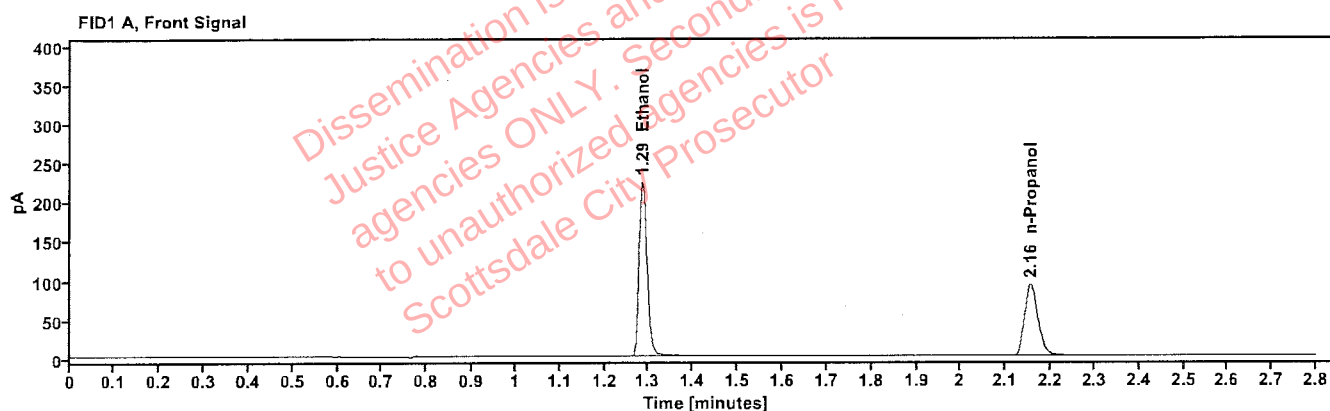


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.4069	1.286	285.902
n-Propanol	-----	2.155	190.846

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.431	295.032
n-Propanol	2.610	194.177

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	529862	Vial:	54
Description:	0.040	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 6:04:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\54.D	Analyst:	Abbott

Abbott

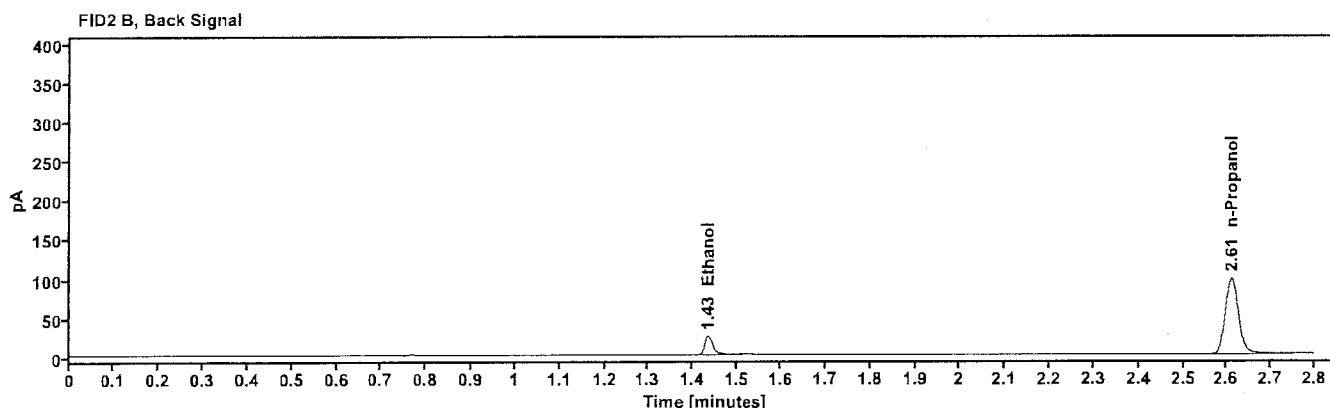
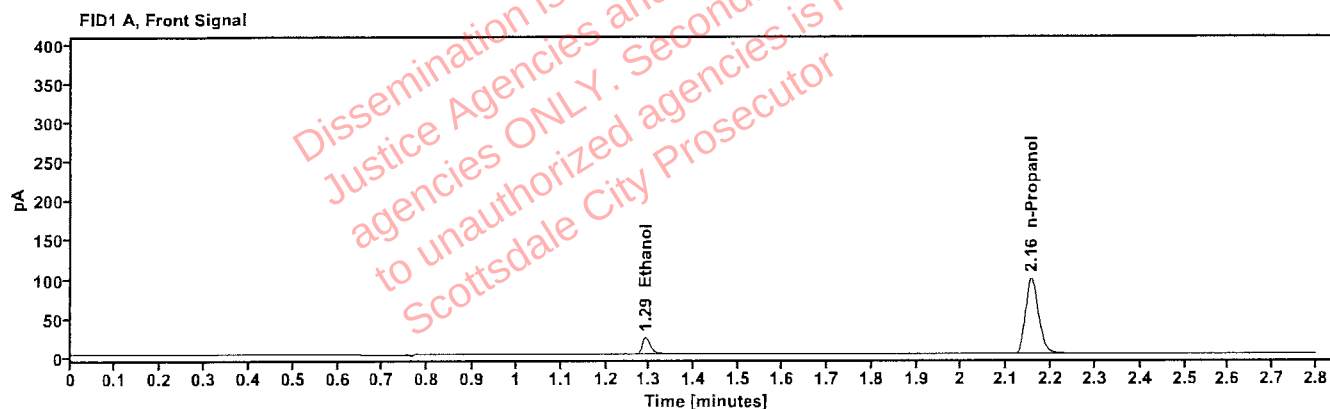


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0410	1.291	29.292
n-Propanol	-----	2.155	201.438

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.434	30.176
n-Propanol	2.610	205.346

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	4190425129/3	Vial:	55
Description:	0.203	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 6:09:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\55.D	Analyst:	Abbott

[Handwritten signature]

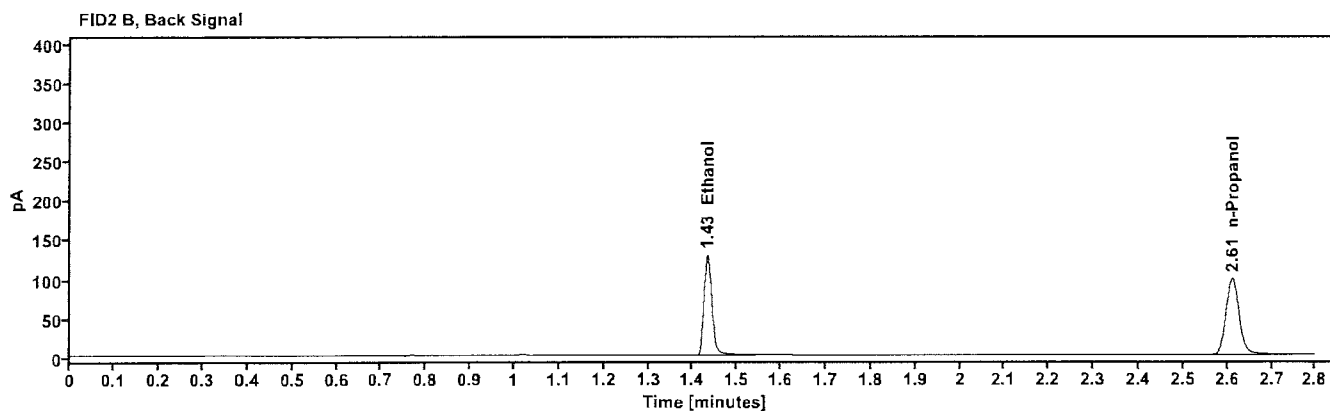
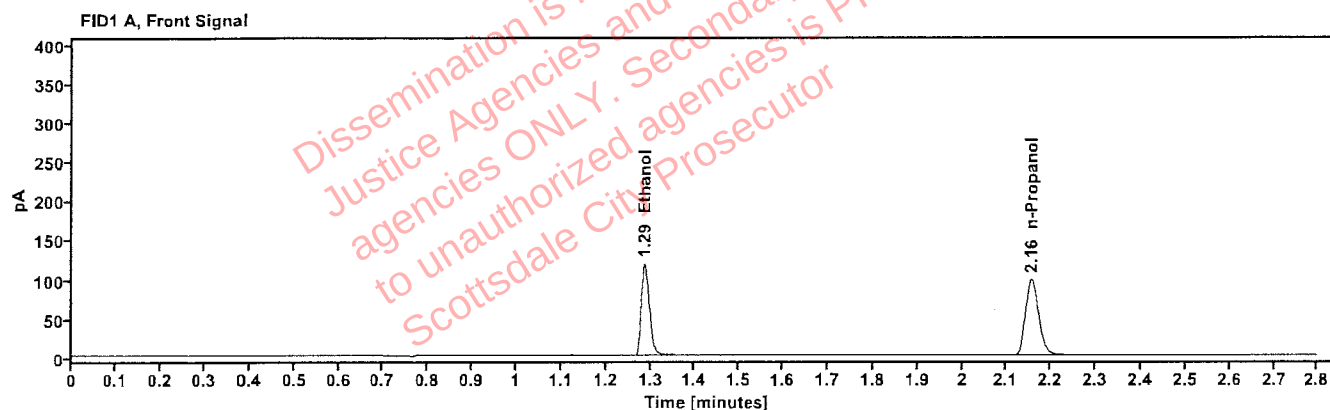


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2043	1.287	151.921
n-Propanol	—	2.155	202.820

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	156.789
n-Propanol	2.610	206.738

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	060326AQ	Vial:	56
Description:	negative	LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 6:12:55 PM
Data file:	C:\Chem32\1\Data\06Aug26\56.D	Analyst:	Abbott <i>JA</i>

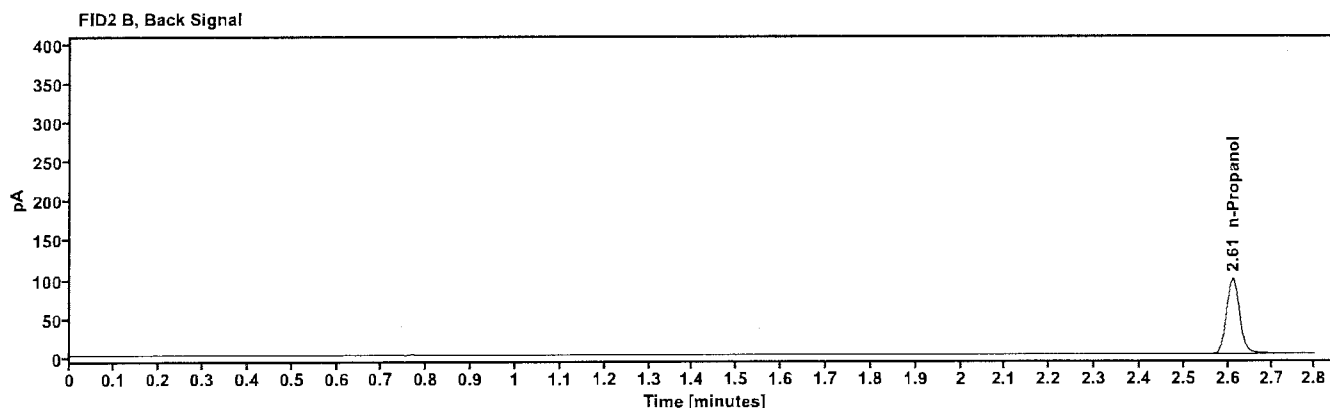
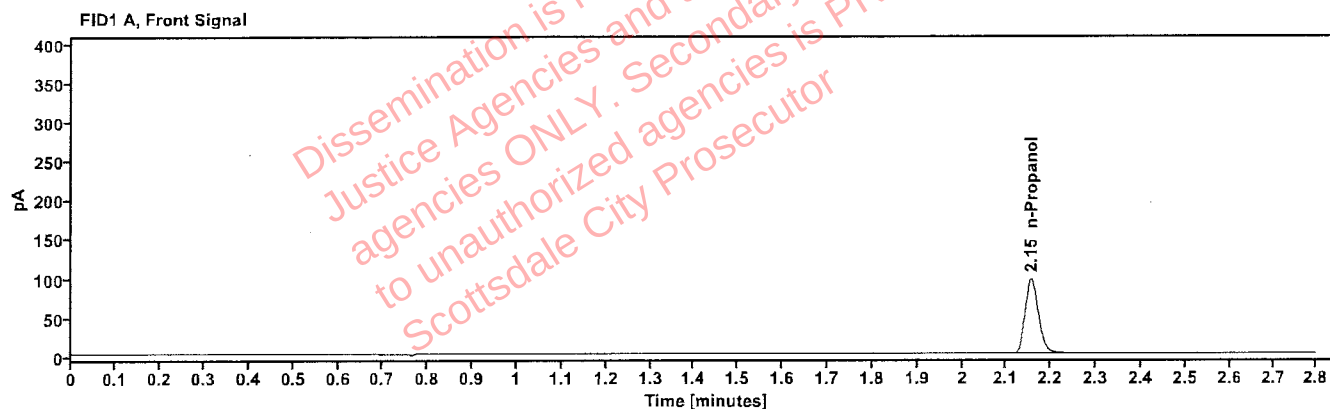


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
n-Propanol	-----	2.155	200.292

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
n-Propanol	2.610	204.317

Sequence Summary

Page 1 of 2

Sequence name: 06Aug26

Instrument: US2320A024 US23180007

Analyst: Abbott

Vial	Sample	Description	Type	LIMS ID	Method
1	FN03062417	0.020 Calibrator	Calibration		ethanol quant.M
2	FN04042408	0.100 Calibrator	Calibration		ethanol quant.M
3	FN07012401	0.200 Calibrator	Calibration		ethanol quant.M
4	FN04092401	0.400 Calibrator	Calibration		ethanol quant.M
5	061026WB	negative	Control		ethanol quant.M
6	031726MIX	Volatile Mix	Control		ethanol quant.M
7	C1475825	0.400	Control		ethanol quant.M
8	529862	0.040	Control		ethanol quant.M
9	4190425129/3	0.203	Control		ethanol quant.M
10	456-01O.P.		Sample		ethanol quant.M
11	456-01O.P.		Sample		ethanol quant.M
12	2613847-3		Sample		ethanol quant.M
13	2613847-3		Sample		ethanol quant.M
14	2614986-9		Sample		ethanol quant.M
15	2614986-9		Sample		ethanol quant.M
16	404-01G.S.		Sample		ethanol quant.M
17	404-01G.S.		Sample		ethanol quant.M
18	446-01J.C.		Sample		ethanol quant.M
19	446-01J.C.		Sample		ethanol quant.M
20	C1455120	0.100	Control		ethanol quant.M
21	444-01J.C.		Sample		ethanol quant.M
22	444-01J.C.		Sample		ethanol quant.M
23	421-01A.P.		Sample		ethanol quant.M
24	421-01A.P.		Sample		ethanol quant.M
25	456-01J.R.		Sample		ethanol quant.M
26	456-01J.R.		Sample		ethanol quant.M
27	434-01C.Z.		Sample		ethanol quant.M
28	434-01C.Z.		Sample		ethanol quant.M
29	445-01M.T.		Sample		ethanol quant.M
30	445-01M.T.		Sample		ethanol quant.M
31	4190425129/3	0.203	Control		ethanol quant.M
32	2615463-2		Sample		ethanol quant.M
33	2615463-2		Sample		ethanol quant.M
34	2615461-2		Sample		ethanol quant.M
35	2615461-2		Sample		ethanol quant.M
36	2615387-1		Sample		ethanol quant.M
37	2615387-1		Sample		ethanol quant.M
38	2615382-2		Sample		ethanol quant.M
39	2615382-2		Sample		ethanol quant.M
40	2615606-1		Sample		ethanol quant.M
41	2615606-1		Sample		ethanol quant.M
42	C1455120	0.100	Control		ethanol quant.M
43	2615599-1		Sample		ethanol quant.M
44	2615599-1		Sample		ethanol quant.M
45	2615540-1		Sample		ethanol quant.M
46	2615540-1		Sample		ethanol quant.M
47	2615059-5		Sample		ethanol quant.M

Sequence Summary

Page 2 of 2



48	2615059-5		Sample		ethanol quant.M
49	2615520-2		Sample		ethanol quant.M
50	2615520-2		Sample		ethanol quant.M
51	2615522-1		Sample		ethanol quant.M
52	2615522-1		Sample		ethanol quant.M
53	C1475825	0.400	Control		ethanol quant.M
54	529862	0.040	Control		ethanol quant.M
55	4190425129/3	0.203	Control		ethanol quant.M
56	060326AQ	negative	Control		ethanol quant.M

Dissemination is restricted to Criminal
Justice Agencies and authorized
agencies ONLY. Secondary dissemination
to unauthorized agencies is PROHIBITED.
Scottsdale City Prosecutor

Scottsdale Police Department Crime Laboratory

Summary of Cases

SEQUENCE NAME: 06Aug26

ANALYST: Abbott

Vials	Test 1 (g/dL)	Test 2 (g/dL)	Mean (g/dL)	Percent Difference*	Absolute Difference(g/dL)*
10 11	0.0878	0.0876	0.08770	0.11	0.00010
12 13	0.2156	0.2162	0.21590	0.14	0.00030
14 15	0.1588	0.1594	0.15910	0.19	0.00030
16 17	0.3057	0.3059	0.30580	0.03	0.00010
18 19	0.2028	0.2034	0.20310	0.15	0.00030
21 22	0.2085	0.2093	0.20890	0.19	0.00040
23 24	0.0000	0.0000	0.00000	0.00	0.00000
25 26	0.2259	0.2264	0.22615	0.11	0.00025
27 28	0.1816	0.1802	0.18090	0.39	0.00070
29 30	0.1537	0.1531	0.15340	0.20	0.00030
32 33	0.1760	0.1731	0.17455	0.83	0.00145
34 35	0.0748	0.0748	0.07480	0.00	0.00000
36 37	0.0853	0.0845	0.08490	0.47	0.00040
38 39	0.2685	0.2656	0.26705	0.54	0.00145
40 41	0.1525	0.1521	0.15230	0.13	0.00020
43 44	0.2413	0.2426	0.24195	0.27	0.00065
45 46	0.0674	0.0682	0.06780	0.59	0.00040
47 48	0.2154	0.2159	0.21565	0.12	0.00025
49 50	0.1426	0.1406	0.14160	0.71	0.00100
51 52	0.1322	0.1331	0.13265	0.34	0.00045

*Calculated differences are differences from the mean of the two results.

Scottsdale Forensic Lab

Blood Alcohol Pipetting Log

ANALYST: Abbott *gh*

SEQUENCE: 06Aug26

[illegible]

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	456-010.P.	Vial:	10
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:08:32 PM
Data file:	C:\Chem32\1\Data\06Aug26\10.D	Analyst:	Abbott

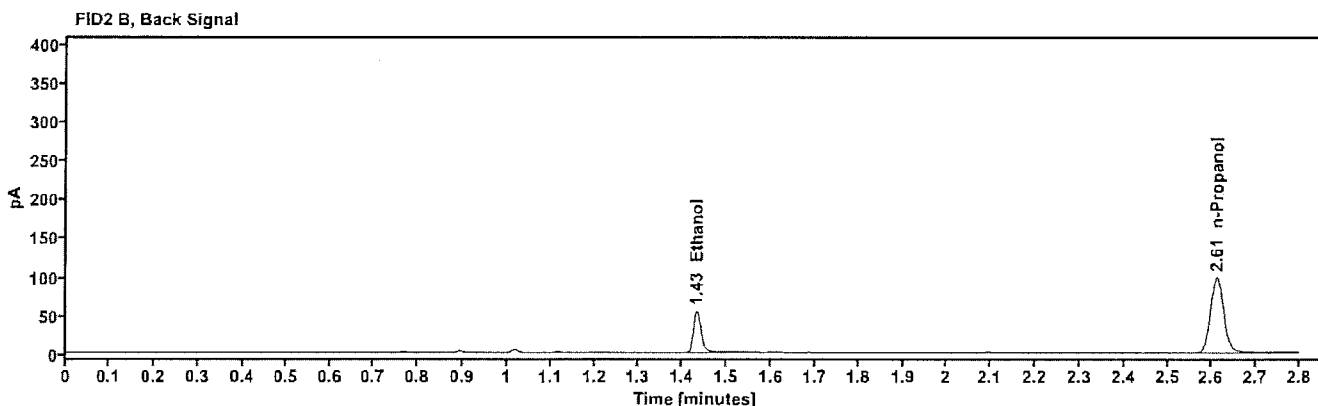
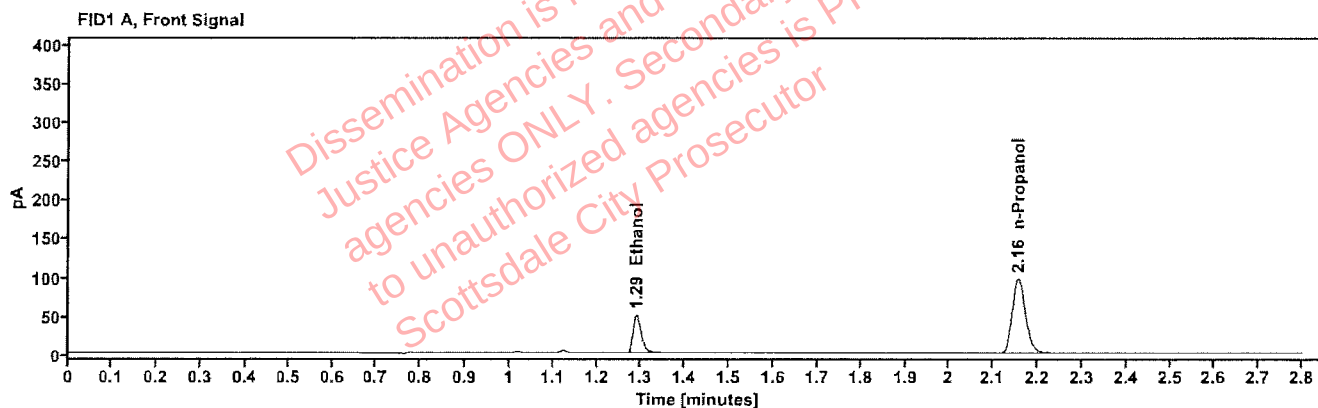


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0878	1.289	63.633
n-Propanol	-----	2.156	199.806

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.433	65.217
n-Propanol	2.611	204.633

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	456-010.P.	Vial:	11
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:12:20 PM
Data file:	C:\Chem32\1\Data\06Aug26\11.D	Analyst:	Abbott

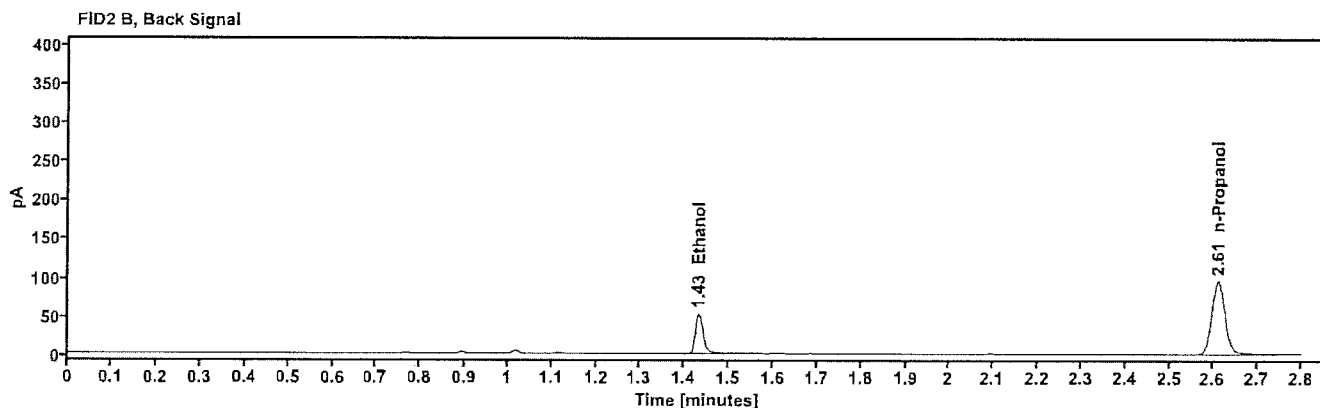
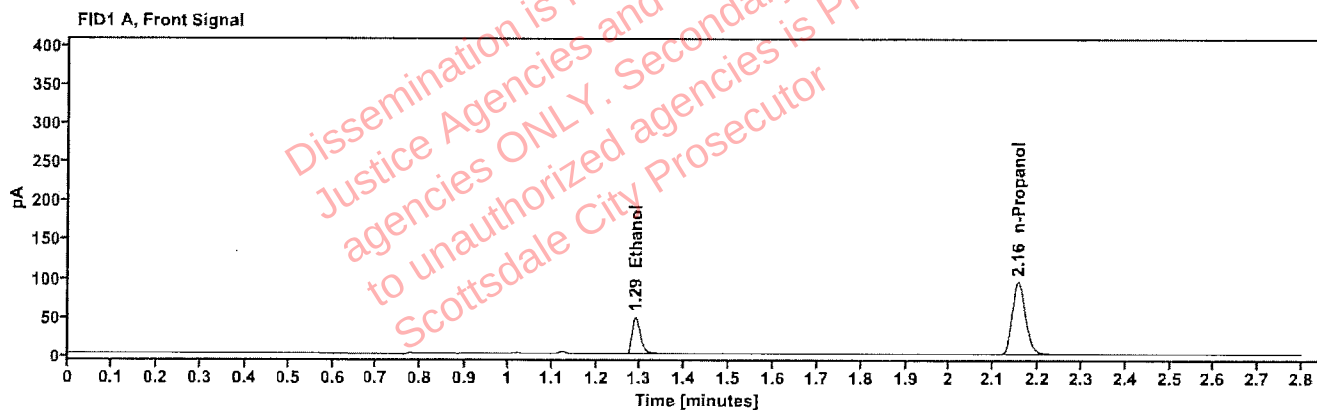


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0876	1.289	61.753
n-Propanol	—	2.155	194.350

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	62.997
n-Propanol	2.610	198.882

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2613847-3	Vial:	12
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:16:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\12.D	Analyst:	Abbott

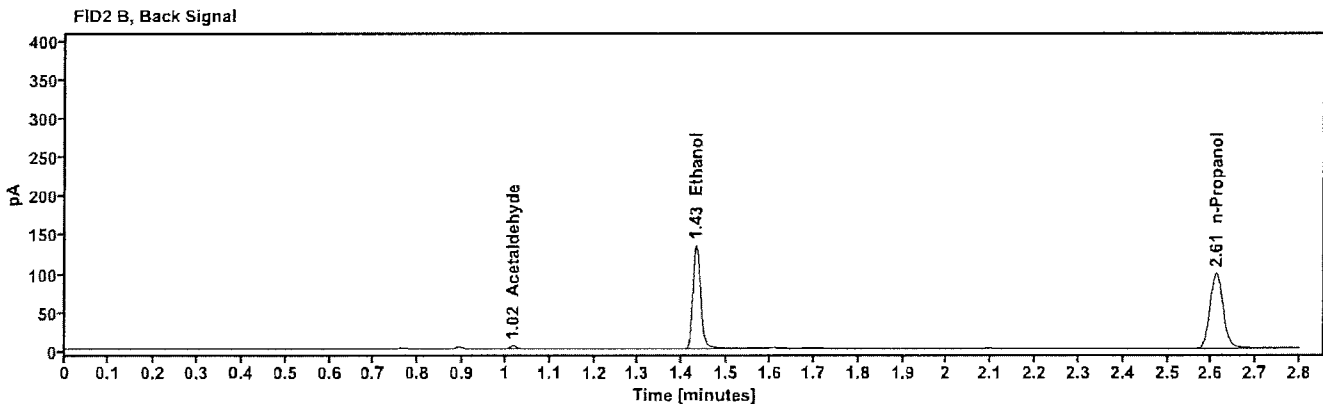
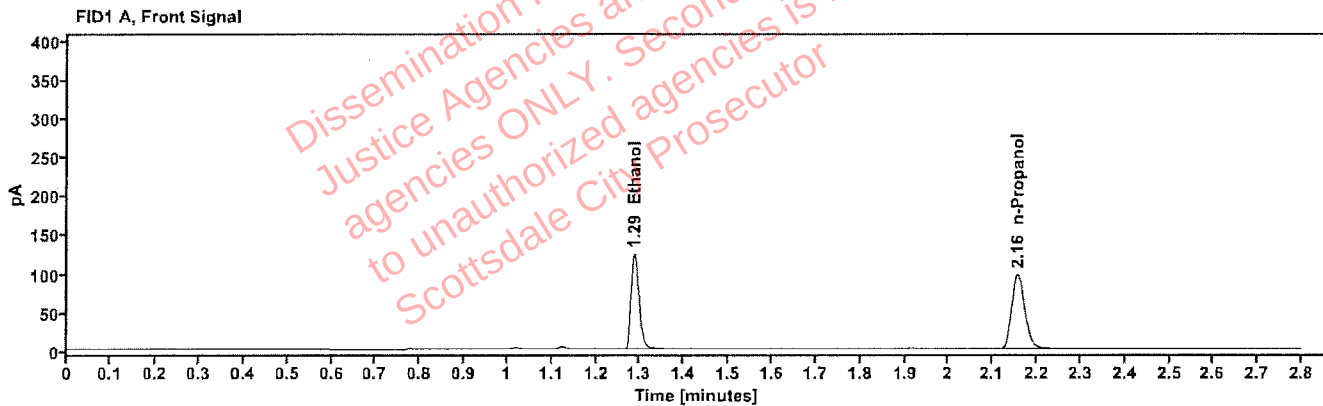


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2156	1.287	157.999
n-Propanol	----	2.155	199.776

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Acetaldehyde	1.017	3.873
Ethanol	1.432	163.043
n-Propanol	2.611	204.766

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2613847-3	Vial:	13
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:21:07 PM
Data file:	C:\Chem32\1\Data\06Aug26\13.D	Analyst:	Abbott

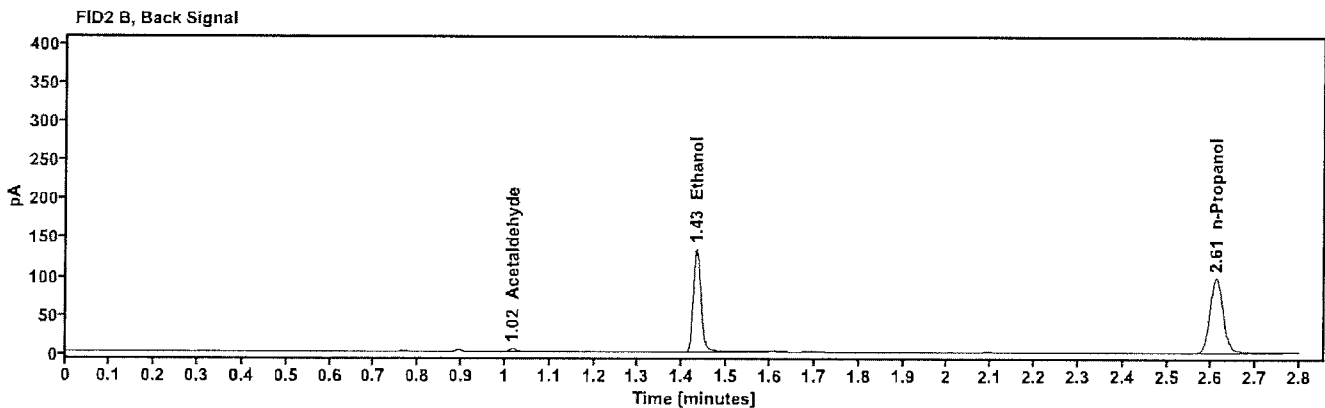
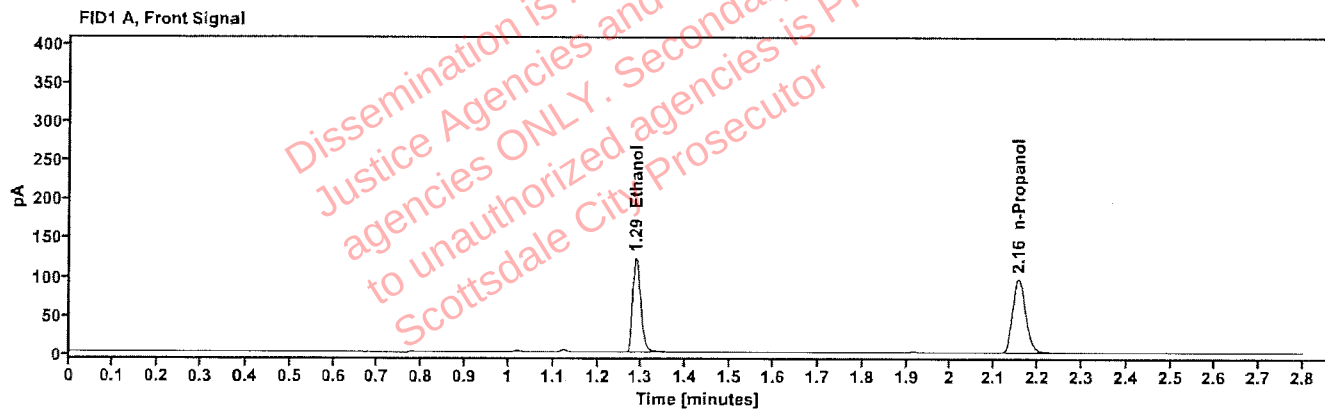


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2162	1.287	157.769
n-Propanol	-----	2.155	198.958

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Acetaldehyde	1.017	3.860
Ethanol	1.431	162.541
n-Propanol	2.611	203.823

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2614986-9	Vial:	14
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:24:56 PM
Data file:	C:\Chem32\1\Data\06Aug26\14.D	Analyst:	Abbott

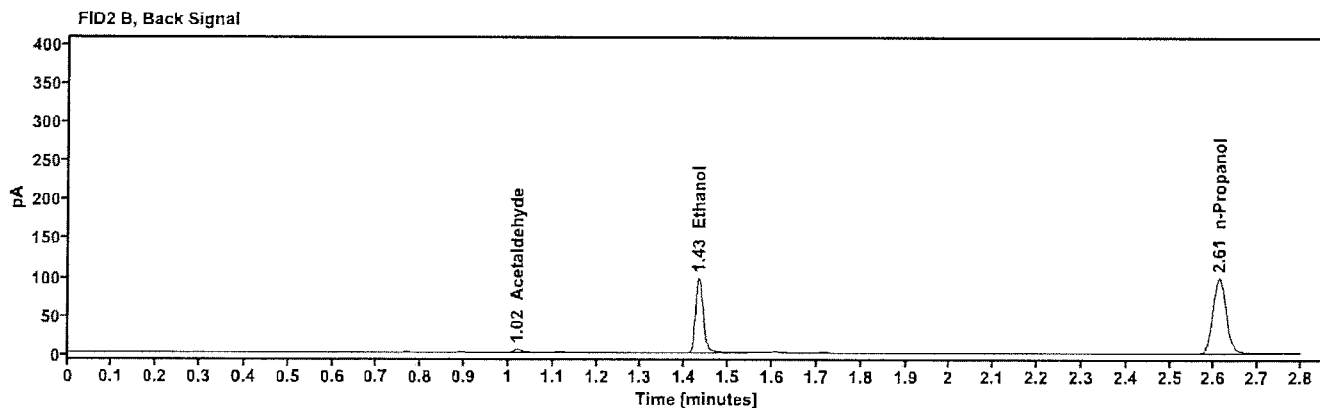
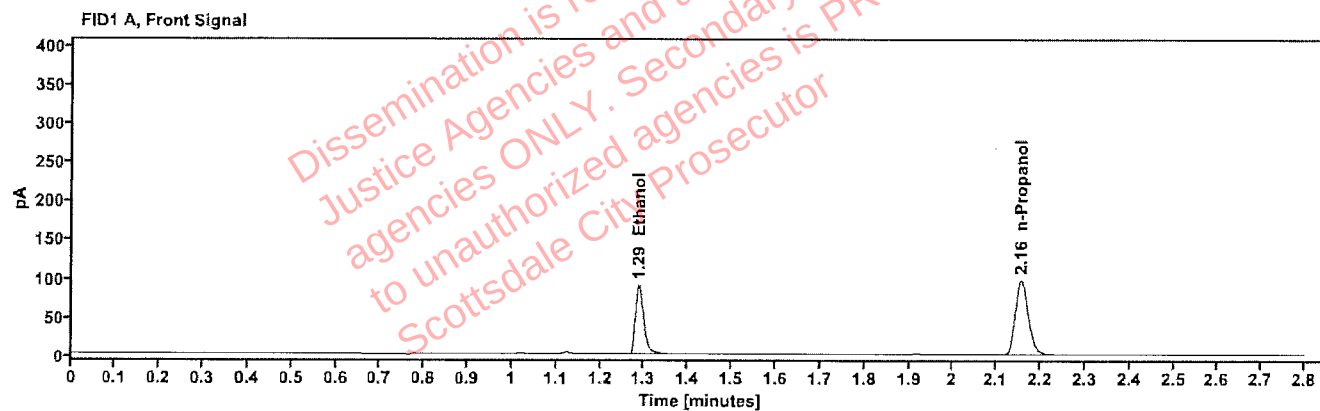


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1588	1.288	116.516
n-Propanol	-----	2.156	200.636

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Acetaldehyde	1.019	4.865
Ethanol	1.432	118.797
n-Propanol	2.611	204.904

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2614986-9	Vial:	15
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:28:44 PM
Data file:	C:\Chem32\1\Data\06Aug26\15.D	Analyst:	Abbott

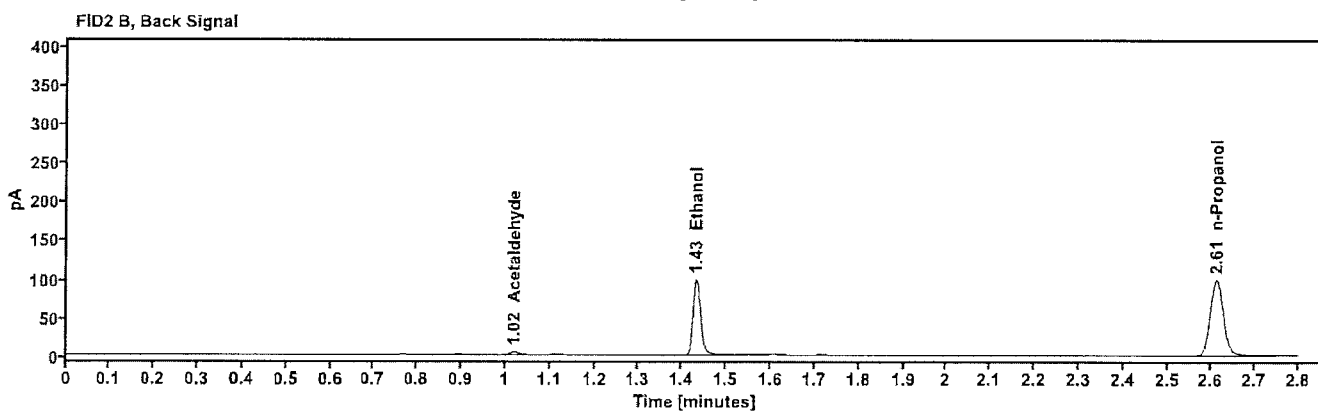
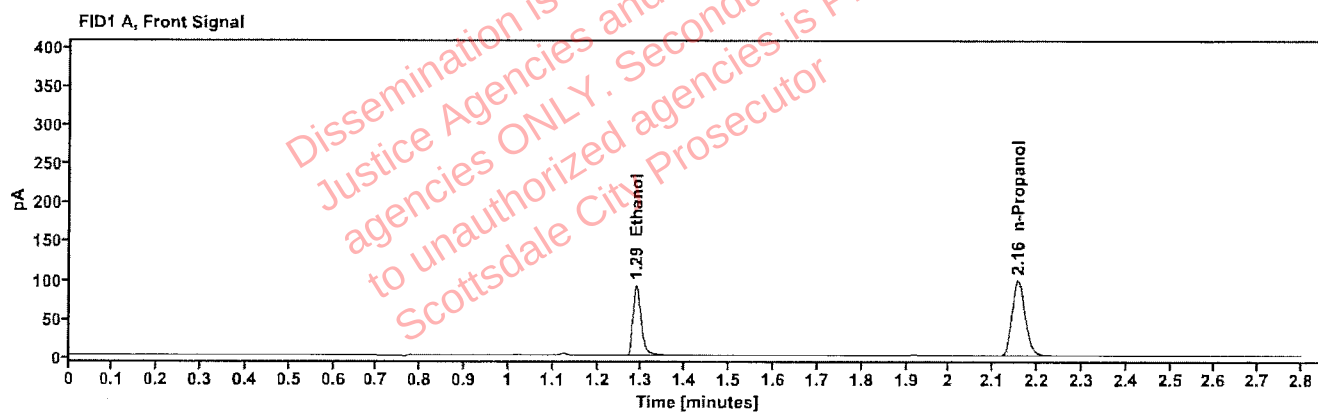


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1594	1.288	117.963
n-Propanol	—	2.156	202.285

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Acetaldehyde	1.019	4.847
Ethanol	1.432	120.395
n-Propanol	2.611	206.385

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	404-01 G.S.	Vial:	16
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:32:32 PM
Data file:	C:\Chem32\1\Data\06Aug26\16.D	Analyst:	Abbott

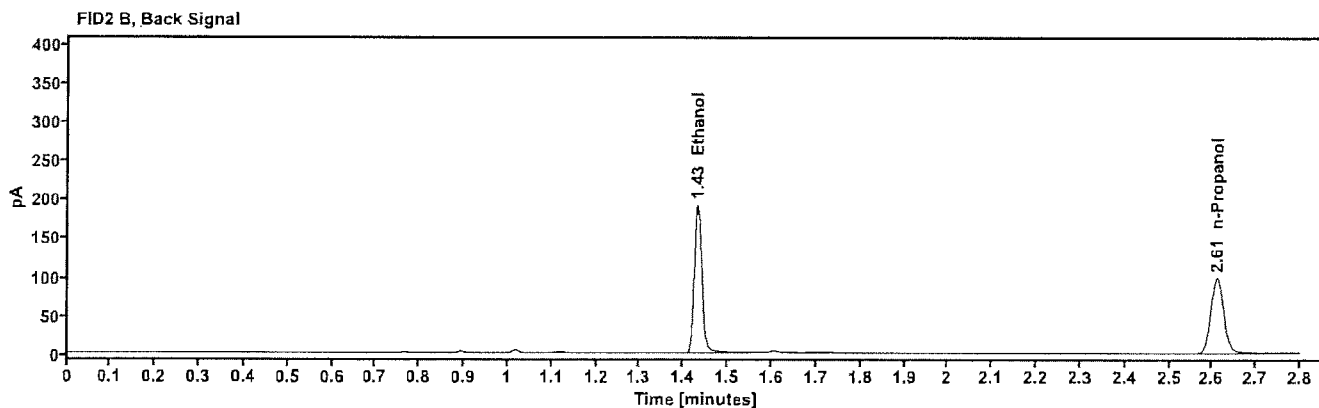
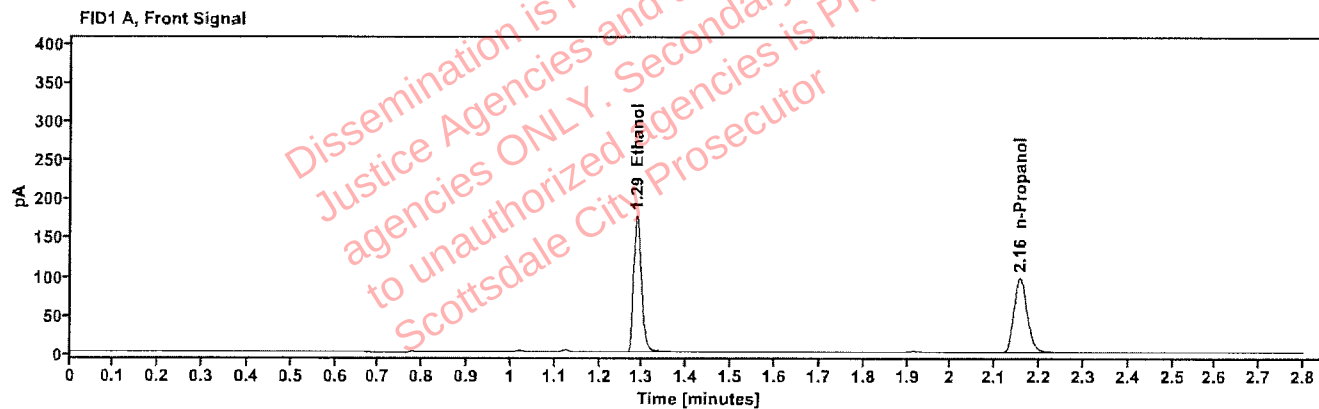


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.3057	1.286	225.111
n-Propanol	—	2.156	200.274

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.431	232.257
n-Propanol	2.611	202.534

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	404-01G.S.	Vial:	17
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:36:20 PM
Data file:	C:\Chem32\11\Data\06Aug26\17.D	Analyst:	Abbott

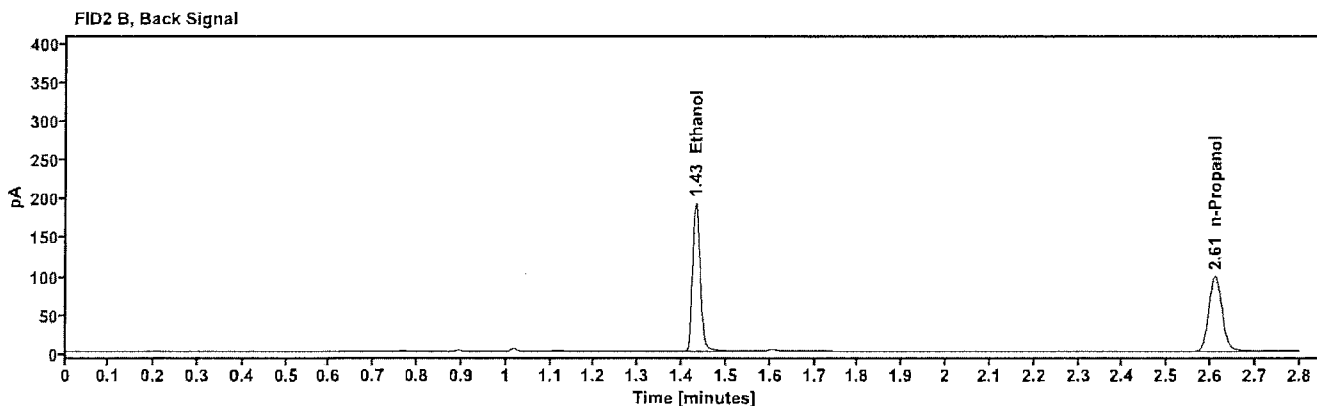
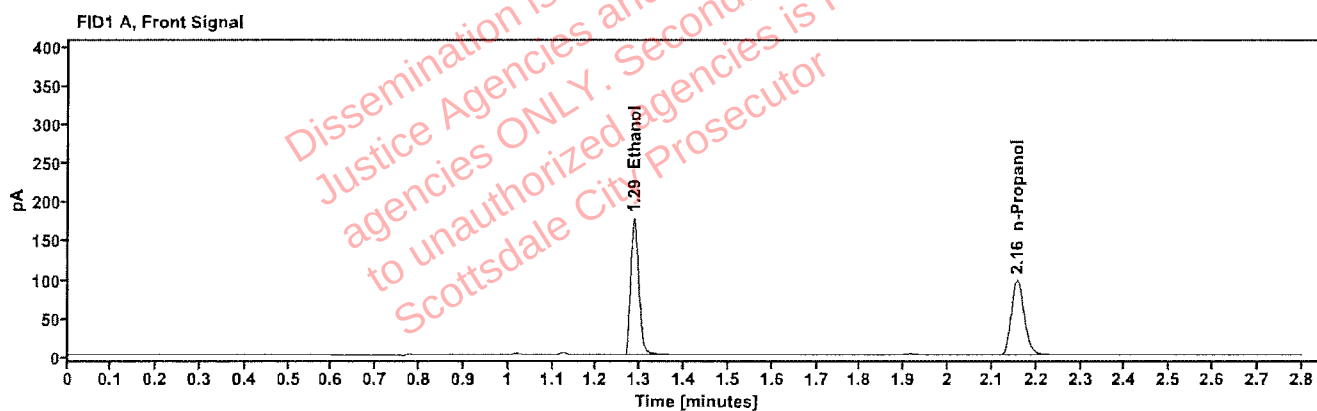


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.3059	1.286	225.869
n-Propanol	-----	2.155	200.851

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.431	233.086
n-Propanol	2.610	203.229

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	446-01J.C.	Vial:	18
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:40:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\18.D	Analyst:	Abbott

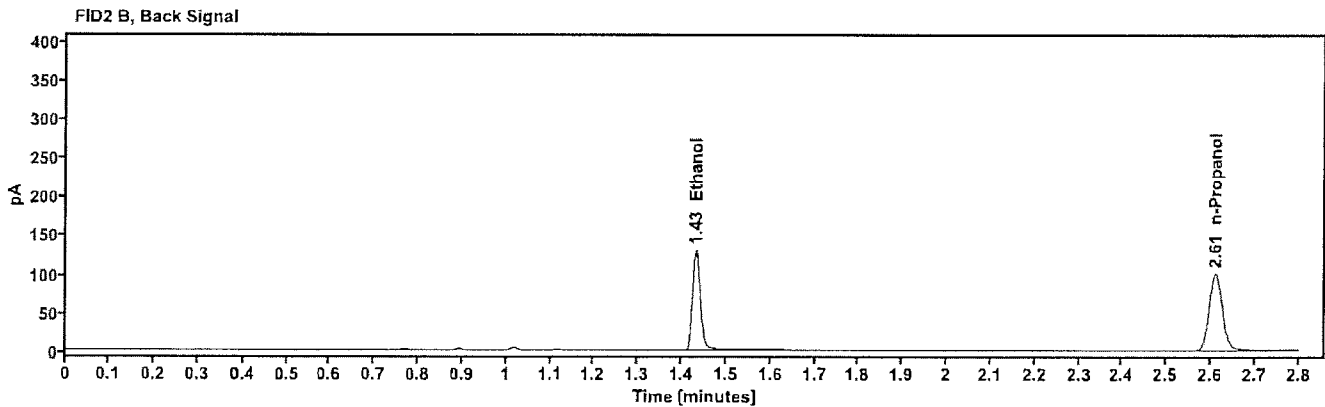
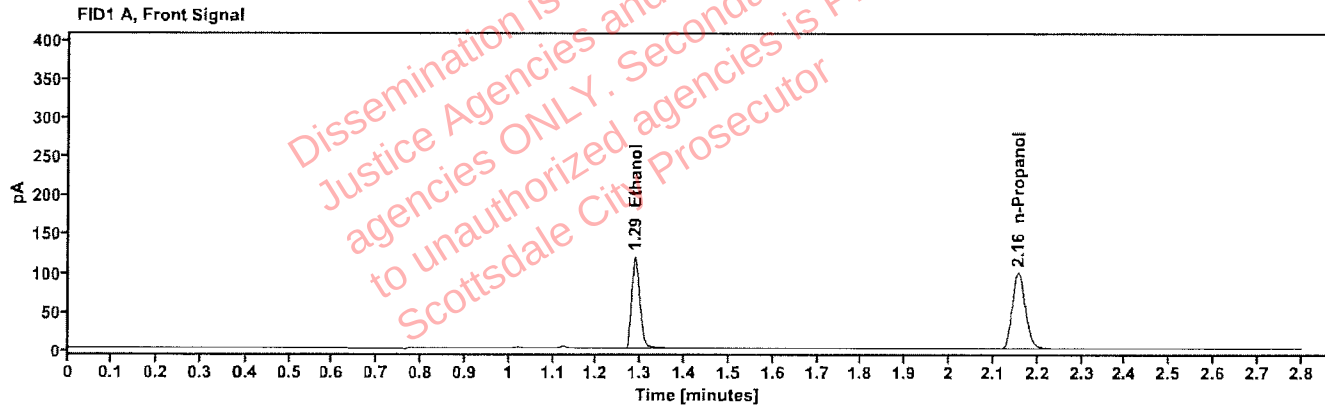


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2028	1.287	152.453
n-Propanol	----	2.156	205.030

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	157.478
n-Propanol	2.611	209.523

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	446-01J.C.	Vial:	19
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:45:07 PM
Data file:	C:\Chem32\1\Data\06Aug26\19.D	Analyst:	Abbott

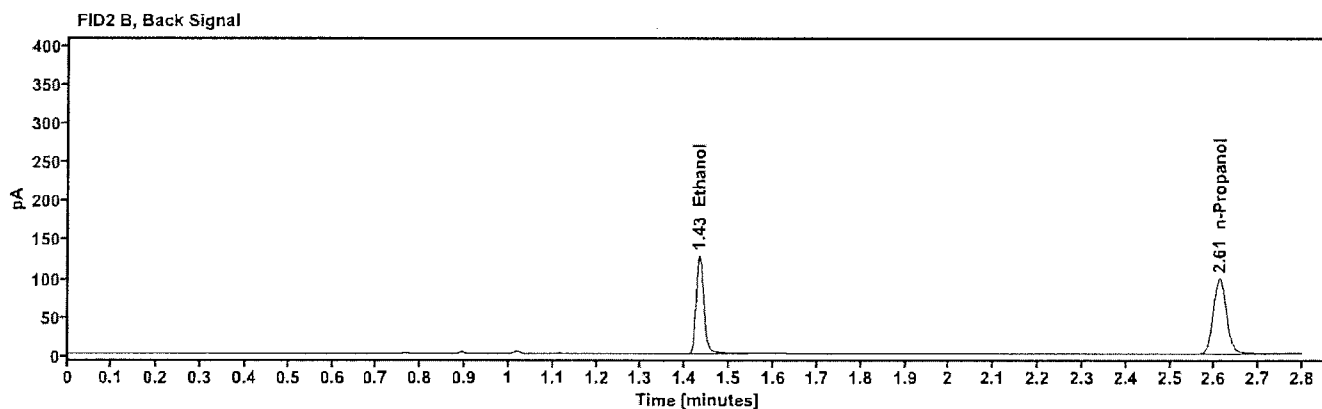
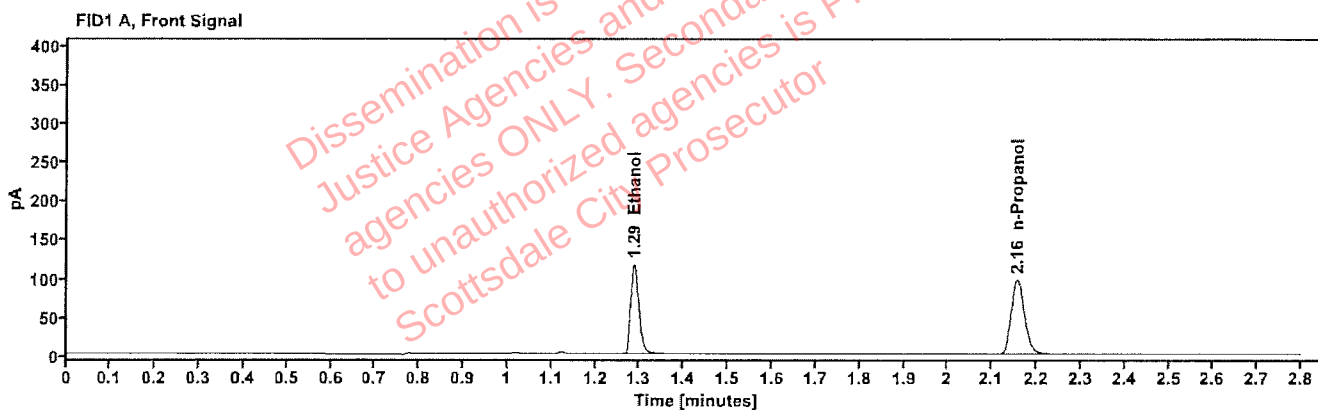


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2034	1.287	149.619
n-Propanol	-----	2.156	200.585

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	154.556
n-Propanol	2.611	205.240

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	444-01J.C.	Vial:	21
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:52:44 PM
Data file:	C:\Chem32\1\Data\06Aug26\21.D	Analyst:	Abbott

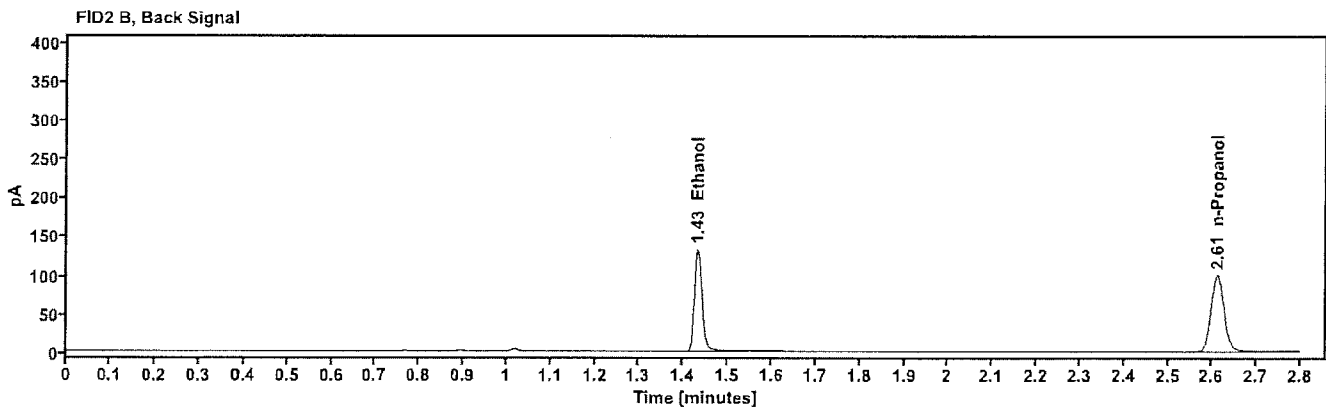
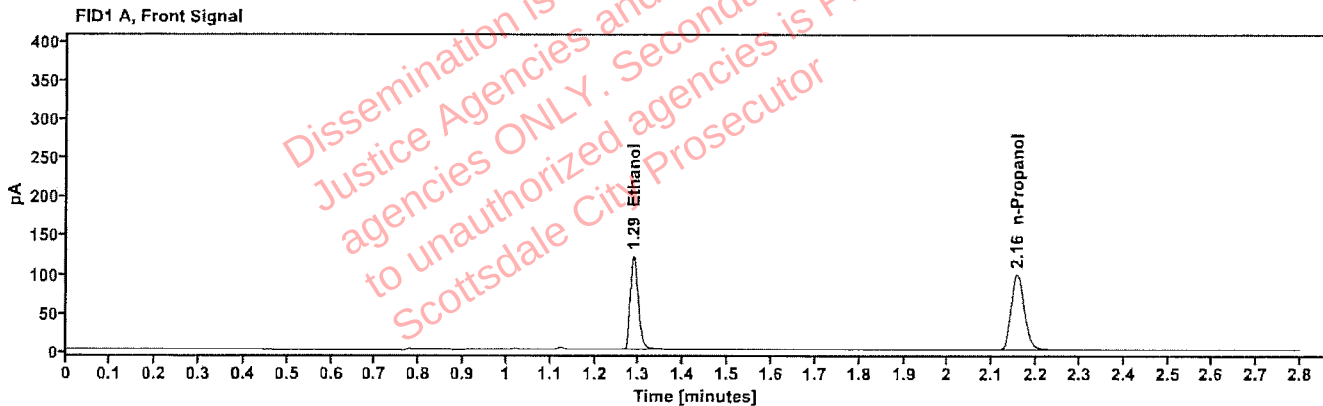


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2085	1.287	155.751
n-Propanol	----	2.156	203.735

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	160.774
n-Propanol	2.611	208.550

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	444-01J.C.	Vial:	22
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 3:56:32 PM
Data file:	C:\Chem32\1\Data\06Aug26\22.D	Analyst:	Abbott

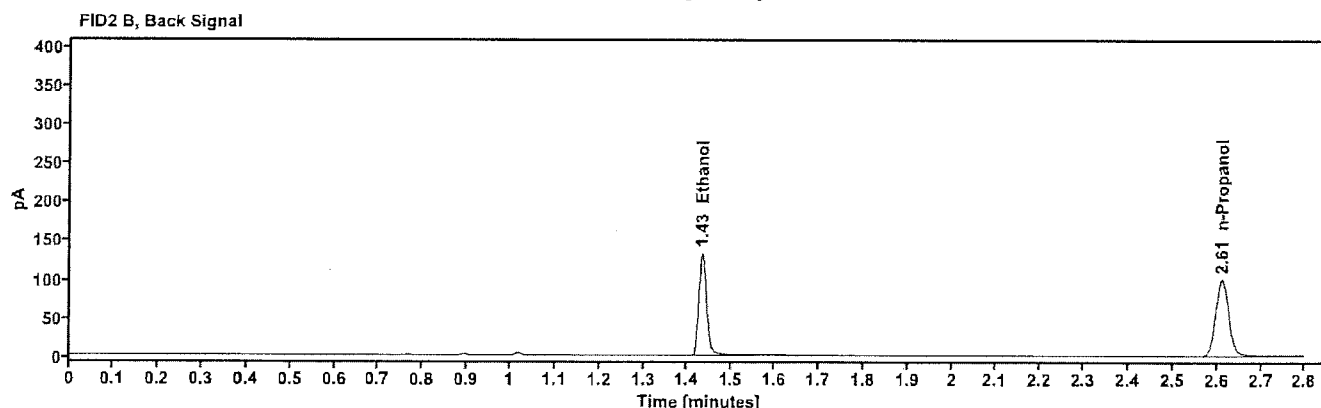
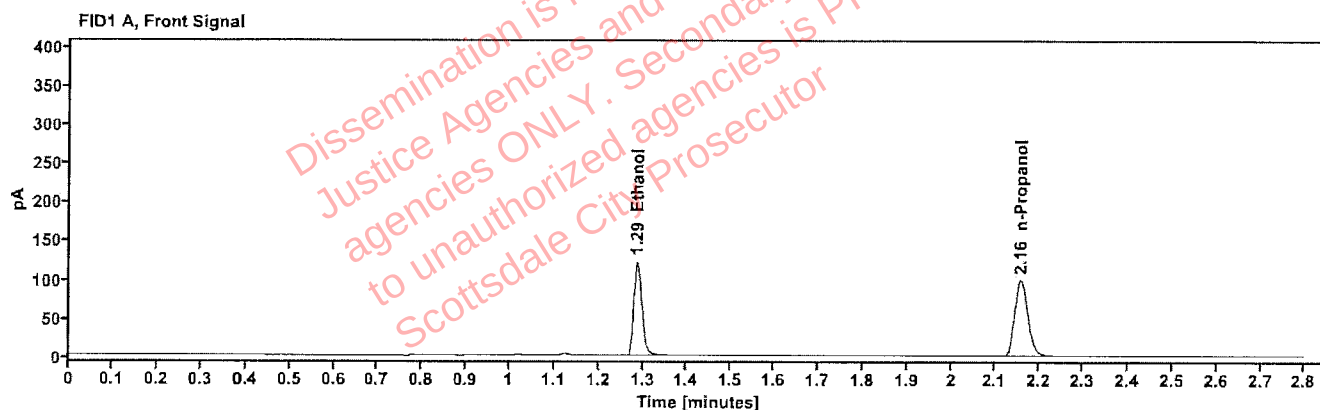


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2093	1.287	156.373
n-Propanol	-----	2.156	203.774

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	161.571
n-Propanol	2.611	208.459

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	421-01A.P.	Vial:	23
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:00:20 PM
Data file:	C:\Chem32\1\Data\06Aug26\23.D	Analyst:	Abbott

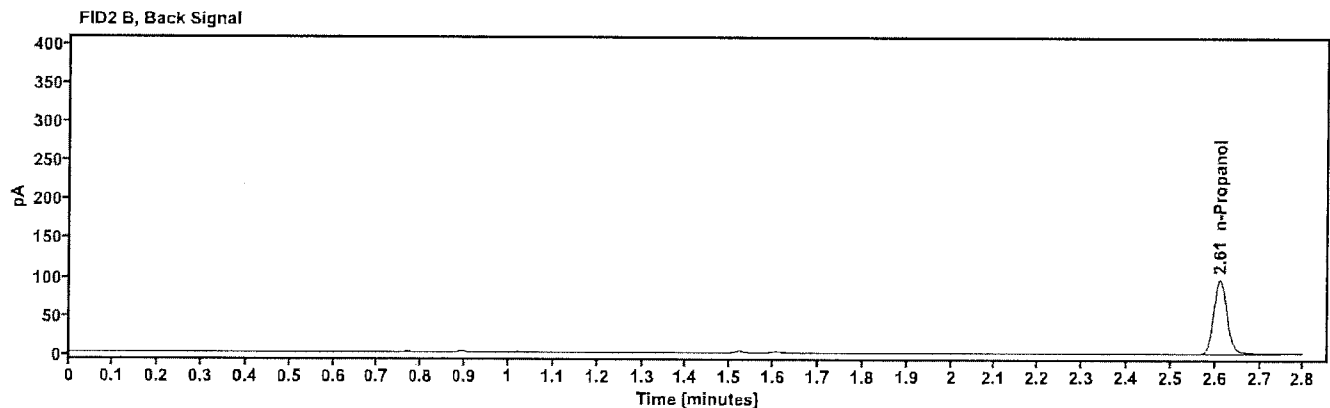
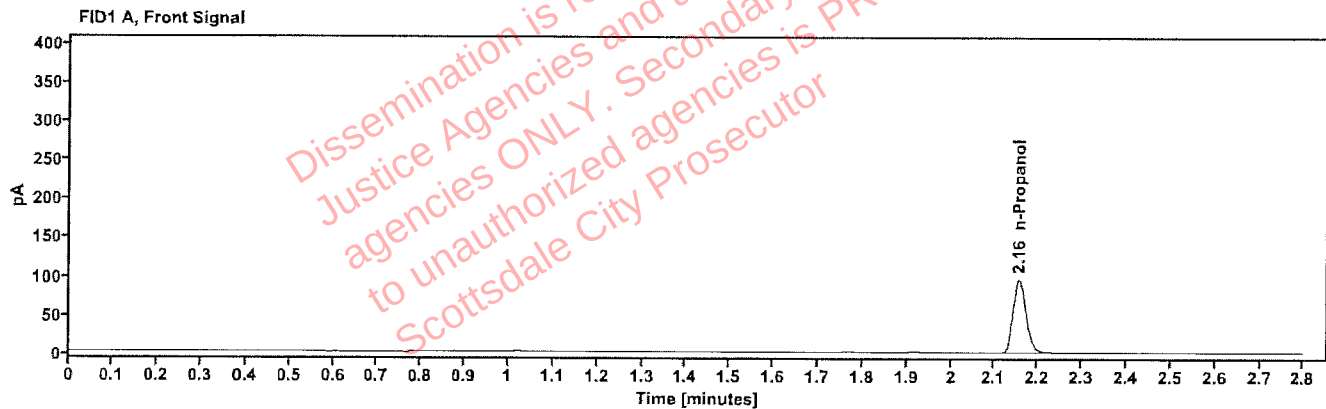


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
n-Propanol	----	2.156	197.444

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
n-Propanol	2.611	201.571

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	421-01A.P.	Vial:	24
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:04:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\24.D	Analyst:	Abbott

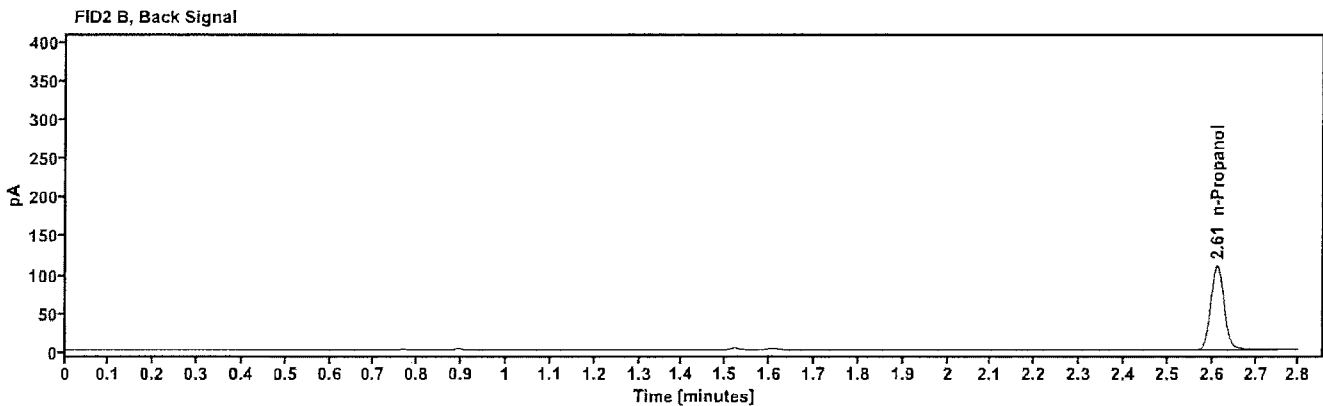
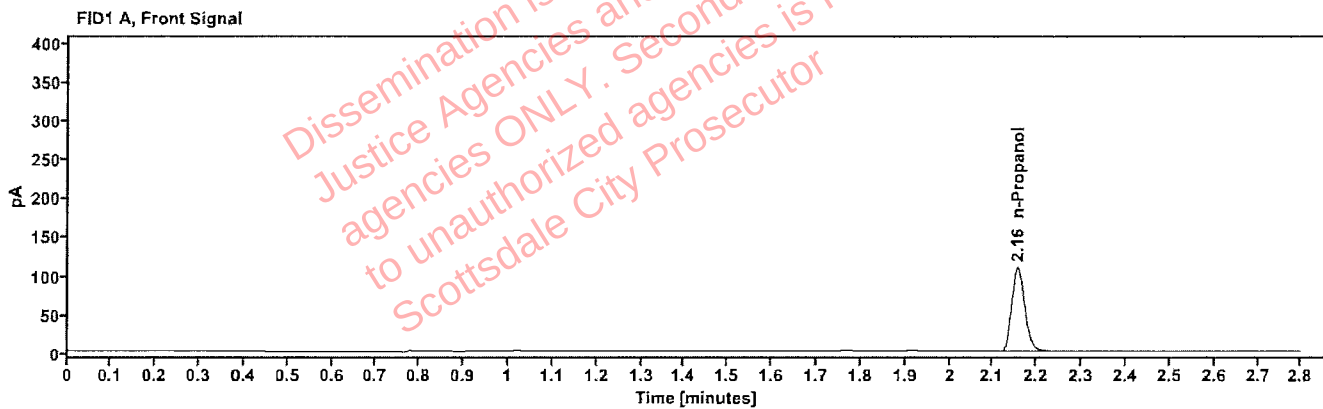


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
n-Propanol	---	2.156	224.415

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
n-Propanol	2.611	229.267

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	456-01J.R.	Vial:	25
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:09:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\25.D	Analyst:	Abbott

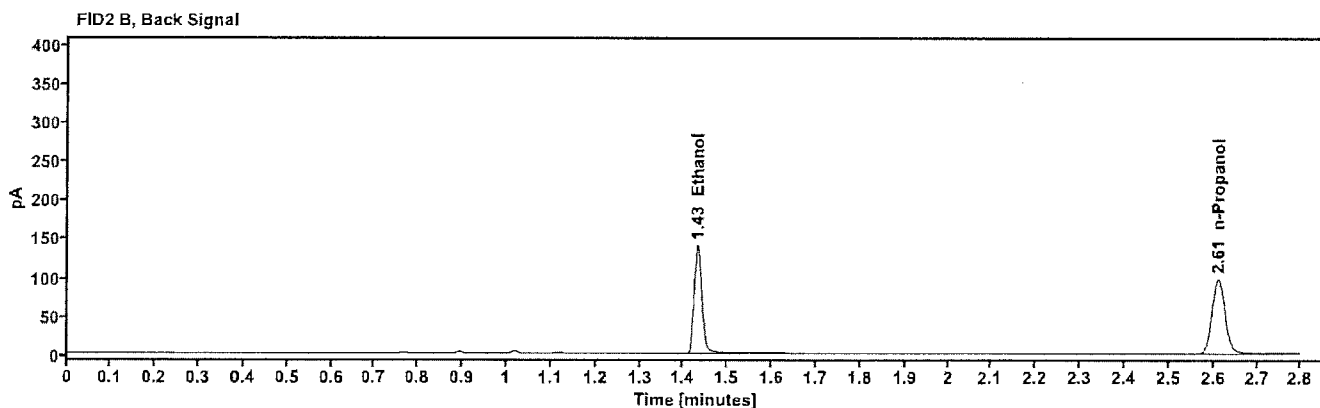
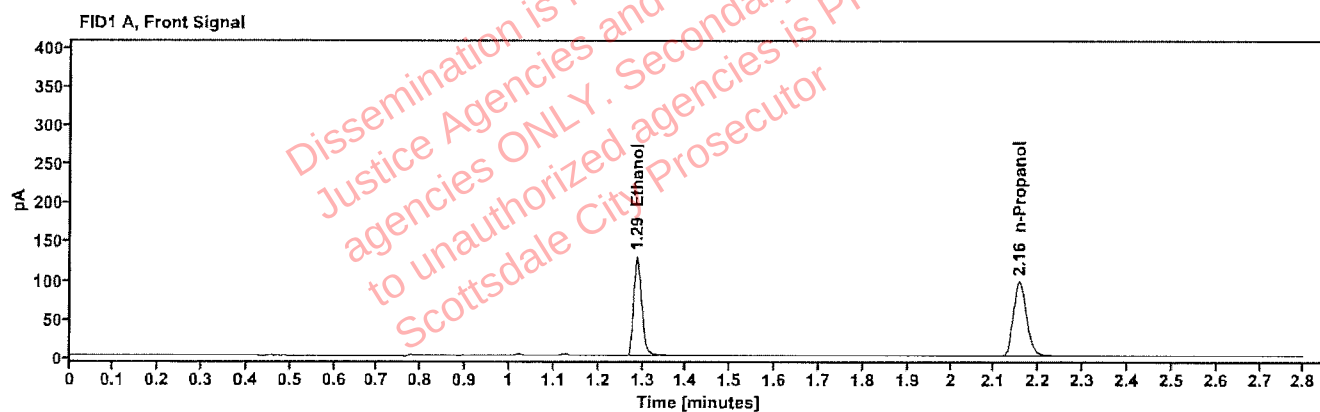


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2259	1.287	165.487
n-Propanol	-----	2.155	199.618

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	170.827
n-Propanol	2.611	204.209

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	456-01J.R.	Vial:	26
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:12:56 PM
Data file:	C:\Chem32\1\Data\06Aug26\26.D	Analyst:	Abbott

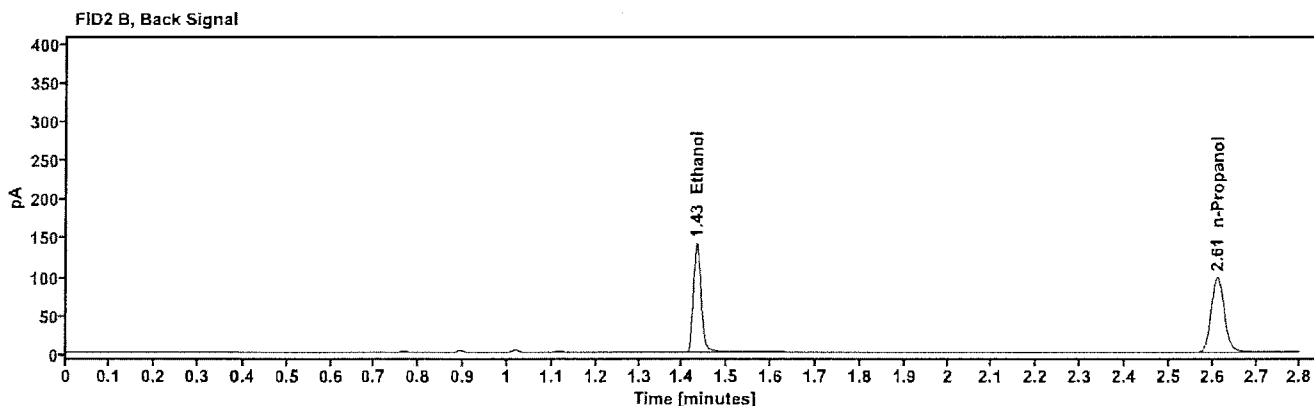
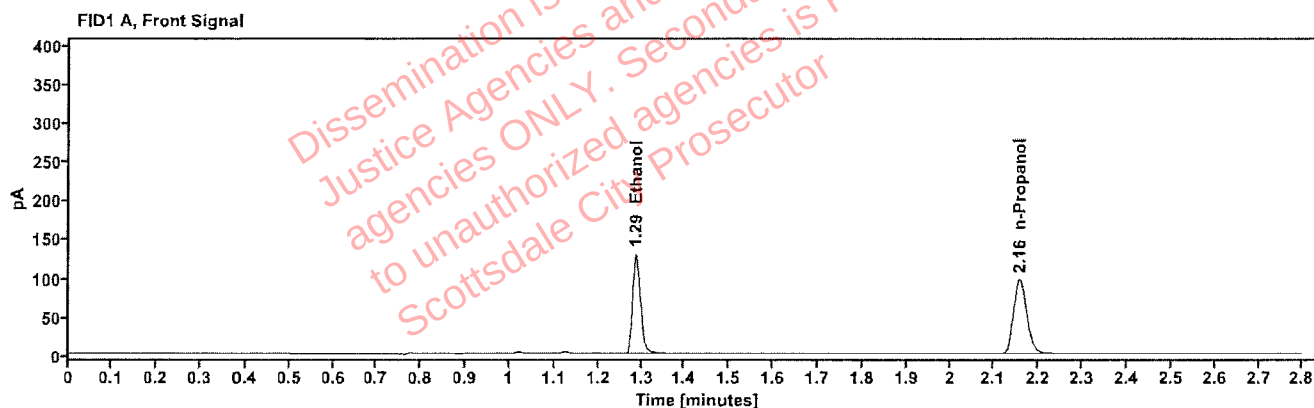


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2264	1.287	166.460
n-Propanol	-----	2.156	200.412

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	172.164
n-Propanol	2.611	204.869

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	434-01C.Z.	Vial:	27
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:16:44 PM
Data file:	C:\Chem32\1\Data\06Aug26\27.D	Analyst:	Abbott

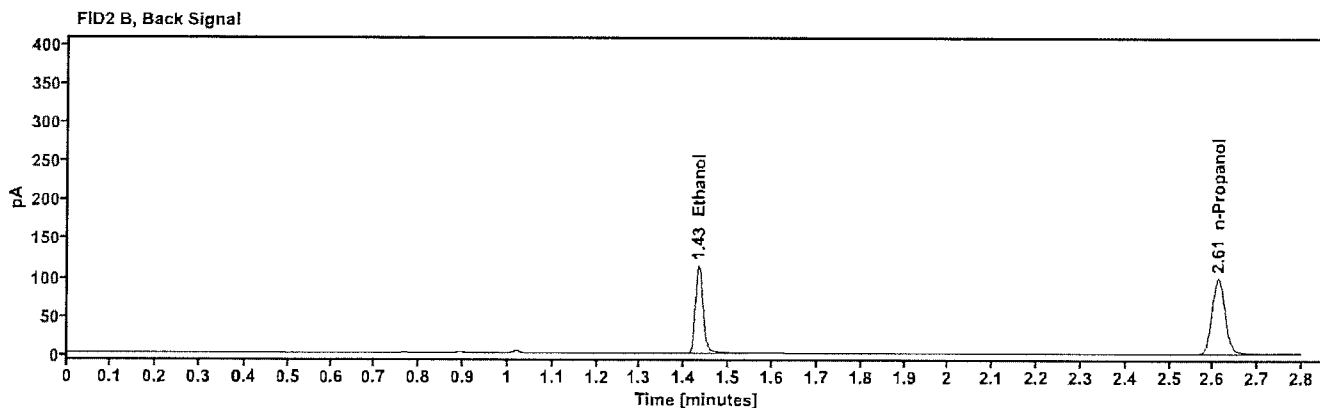
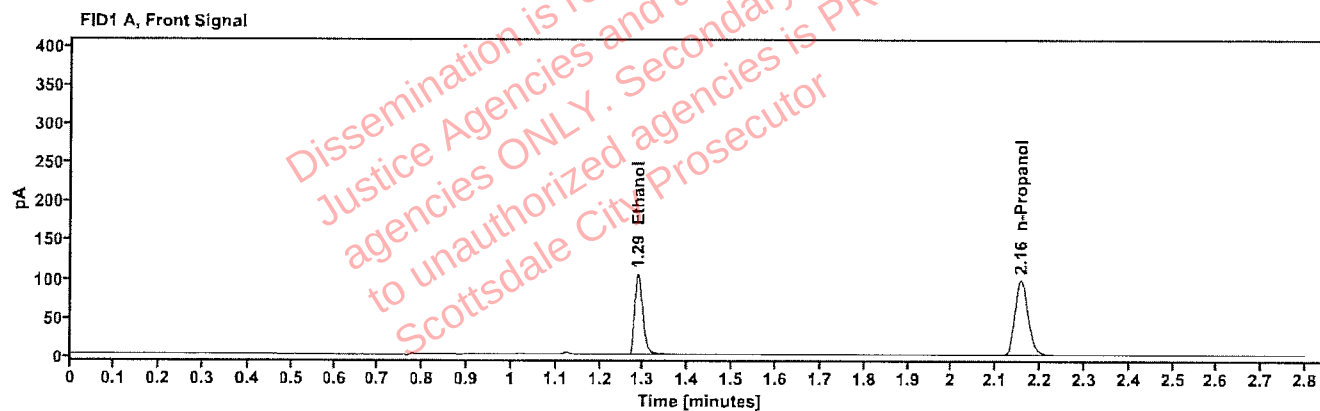


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1816	1.288	134.109
n-Propanol	----	2.156	201.632

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	138.360
n-Propanol	2.611	205.977

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	434-01C.Z.	Vial:	28
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:20:32 PM
Data file:	C:\Chem32\1\Data\06Aug26\28.D	Analyst:	Abbott

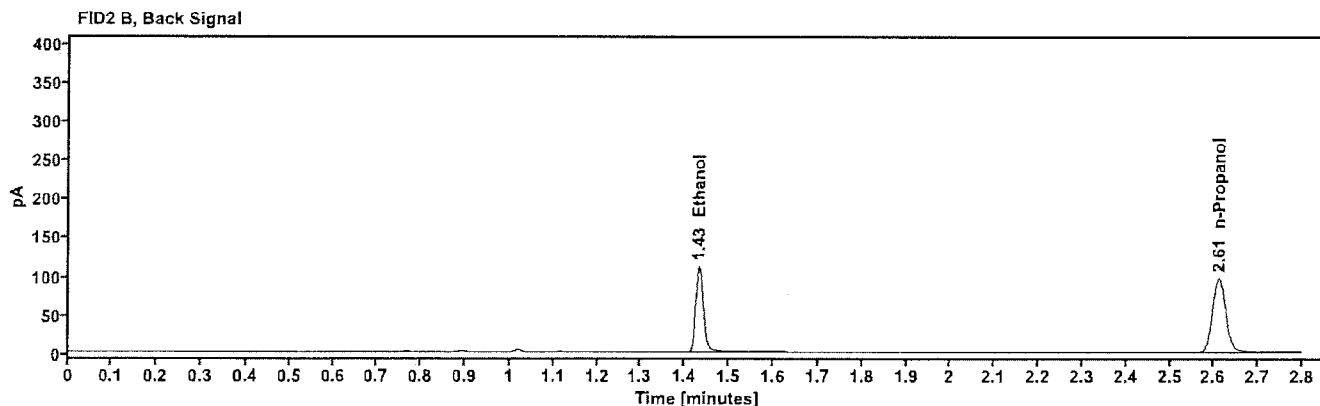
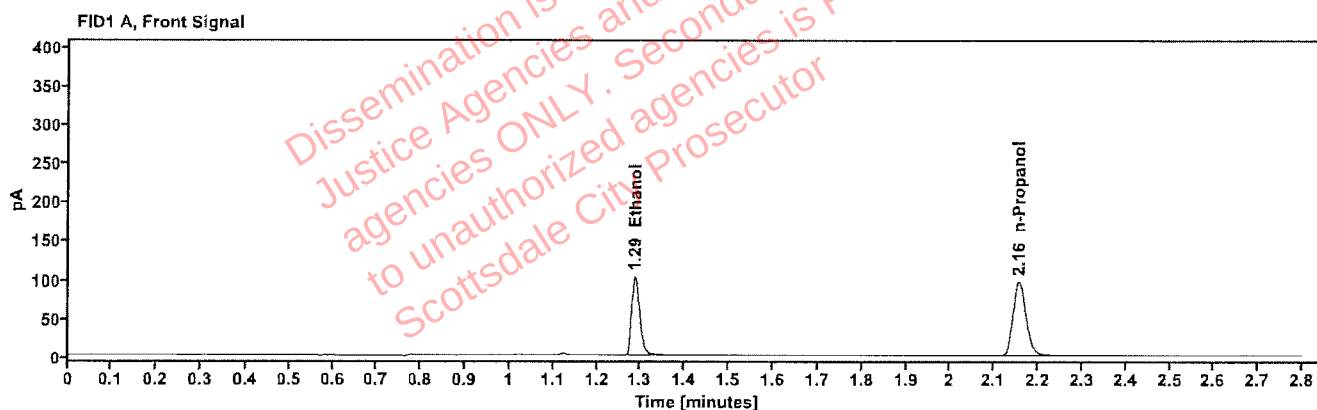


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1802	1.287	131.134
n-Propanol	----	2.155	198.642

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	135.806
n-Propanol	2.611	202.841

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	445-01M.T.	Vial:	29
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:24:20 PM
Data file:	C:\Chem32\1\Data\06Aug26\29.D	Analyst:	Abbott

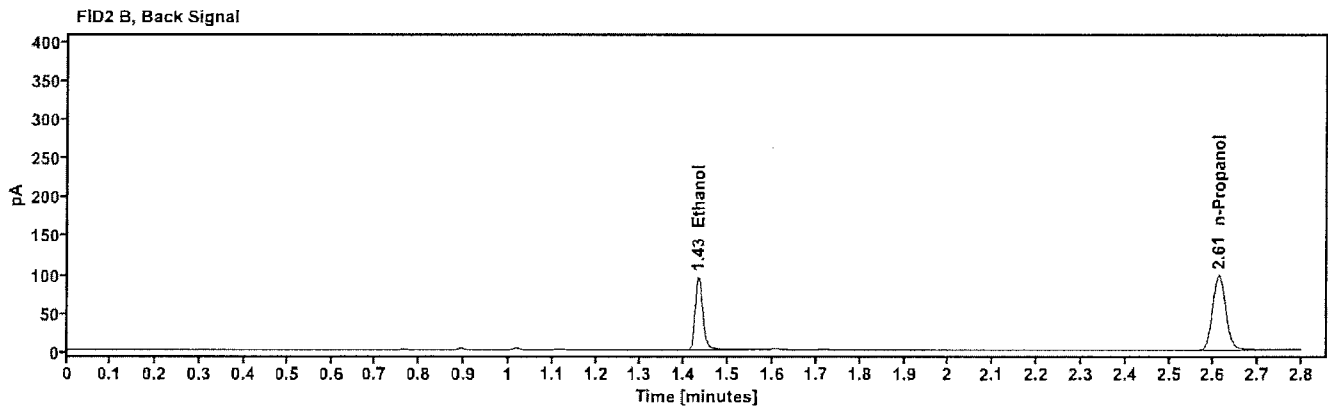
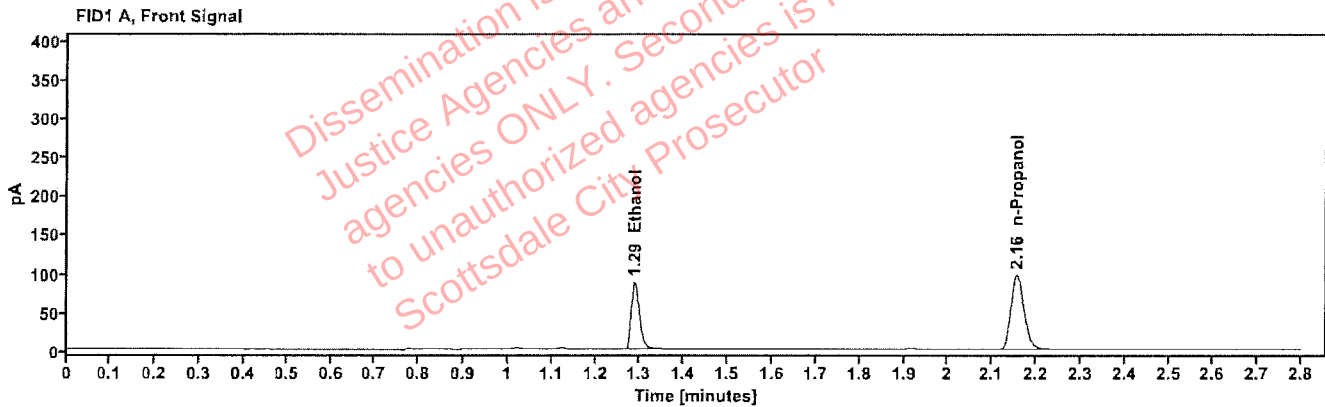


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1537	1.288	112.028
n-Propanol	-----	2.156	199.368

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	115.820
n-Propanol	2.611	203.924

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	445-01M.T.	Vial:	30
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:28:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\30.D	Analyst:	Abbott

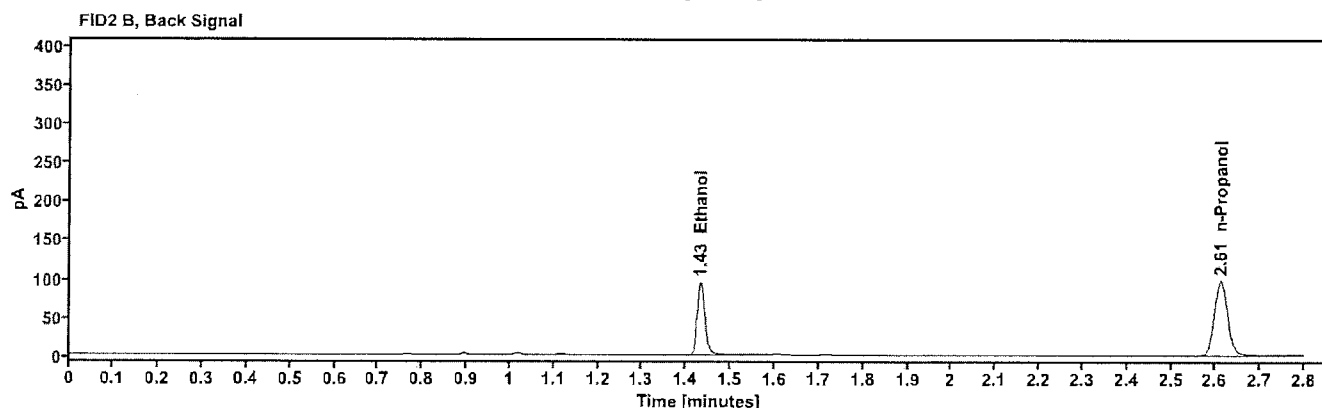
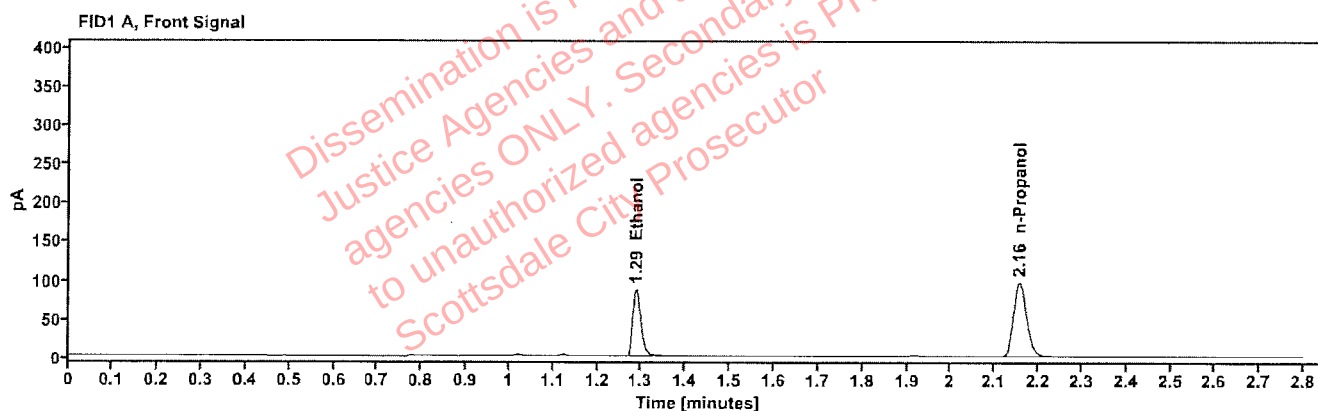


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1531	1.288	111.541
n-Propanol	-----	2.156	199.263

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	115.277
n-Propanol	2.611	203.772

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615463-2	Vial:	32
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:36:56 PM
Data file:	C:\Chem32\1\Data\06Aug26\32.D	Analyst:	Abbott

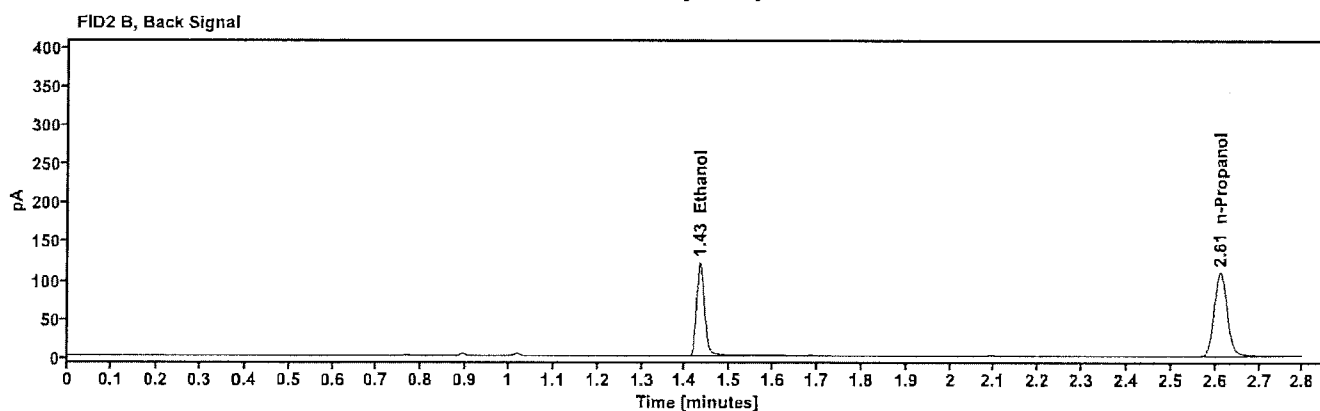
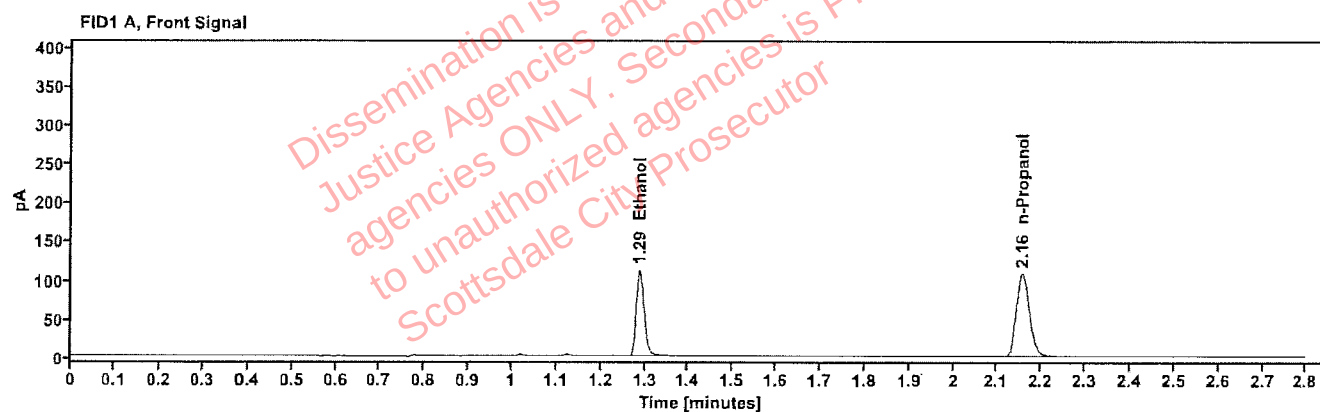


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1760	1.288	142.993
n-Propanol	-----	2.156	221.869

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	147.615
n-Propanol	2.611	227.107

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615463-2	Vial:	33
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:40:44 PM
Data file:	C:\Chem32\1\Data\06Aug26\33.D	Analyst:	Abbott

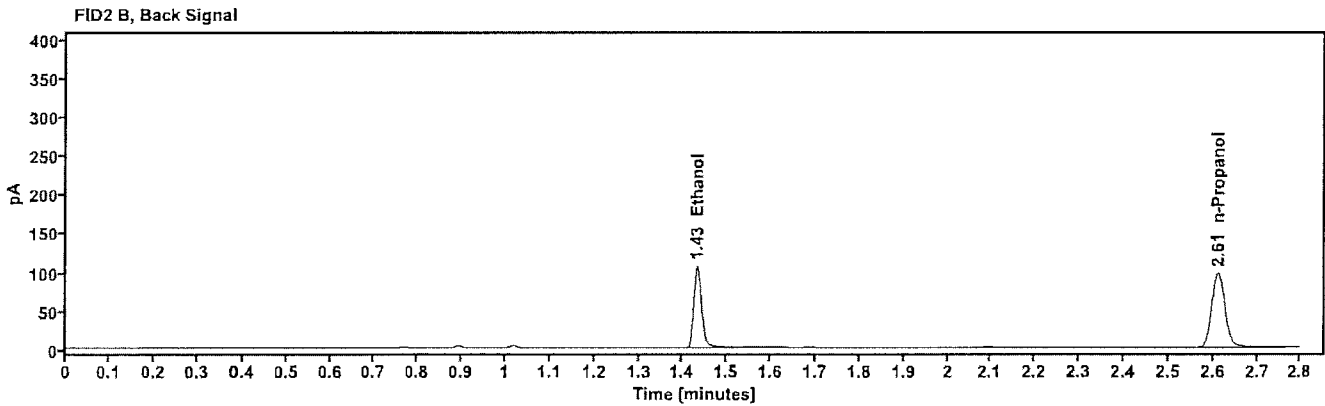
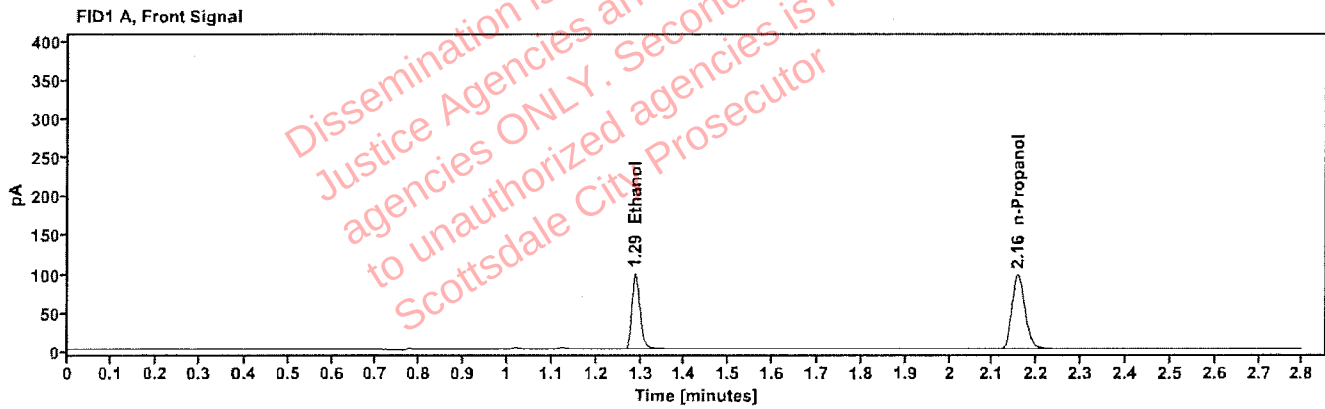


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1731	1.288	125.991
n-Propanol	---	2.156	198.815

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	130.052
n-Propanol	2.611	203.681

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615461-2	Vial:	34
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:44:31 PM
Data file:	C:\Chem32\1\Data\06Aug26\34.D	Analyst:	Abbott

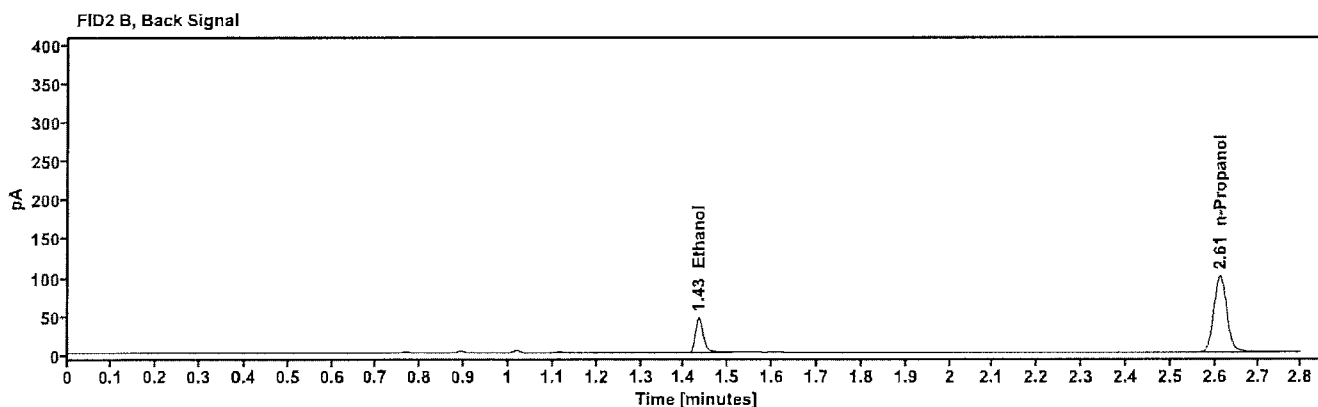
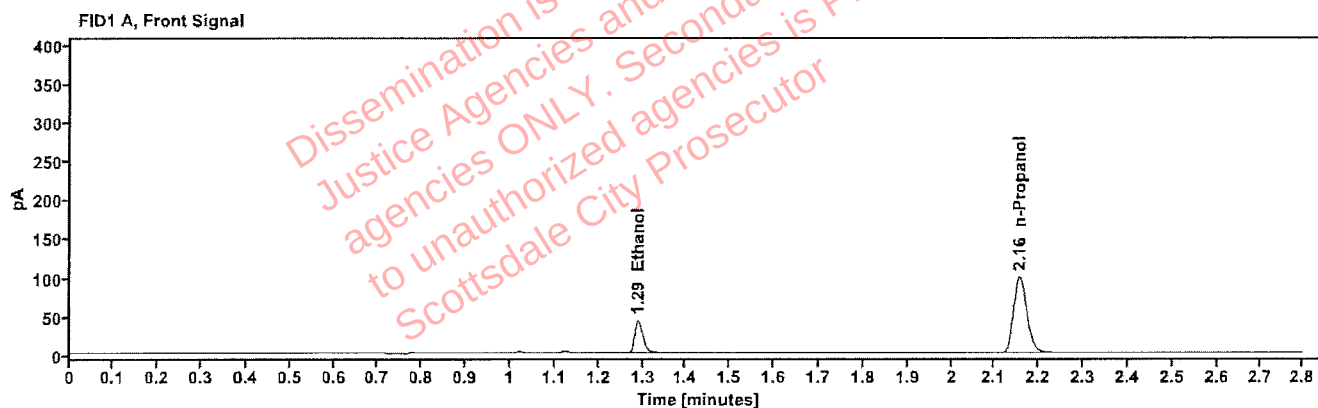


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0748	1.289	55.109
n-Propanol	-----	2.156	203.859

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.433	56.313
n-Propanol	2.611	208.433

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615461-2	Vial:	35
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:48:20 PM
Data file:	C:\Chem32\1\Data\06Aug26\35.D	Analyst:	Abbott

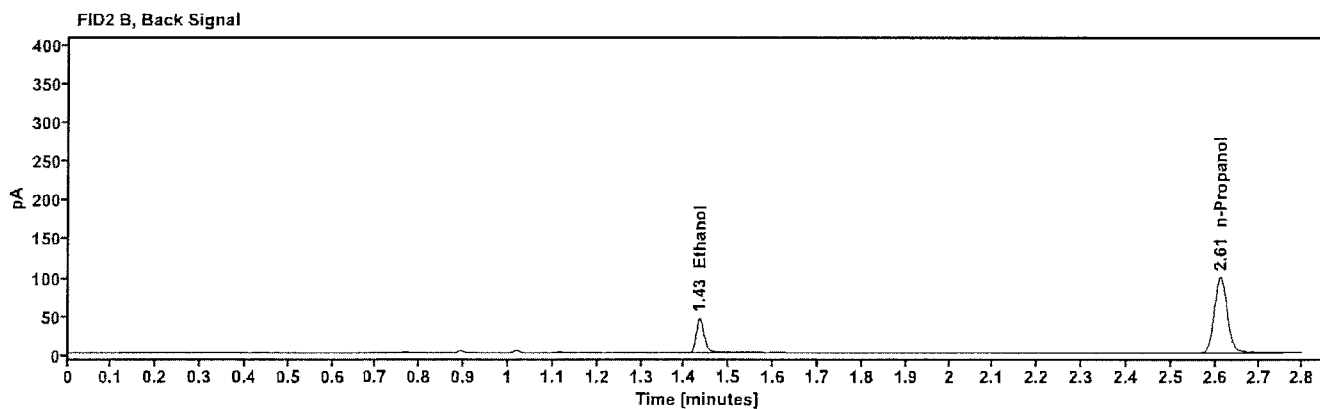
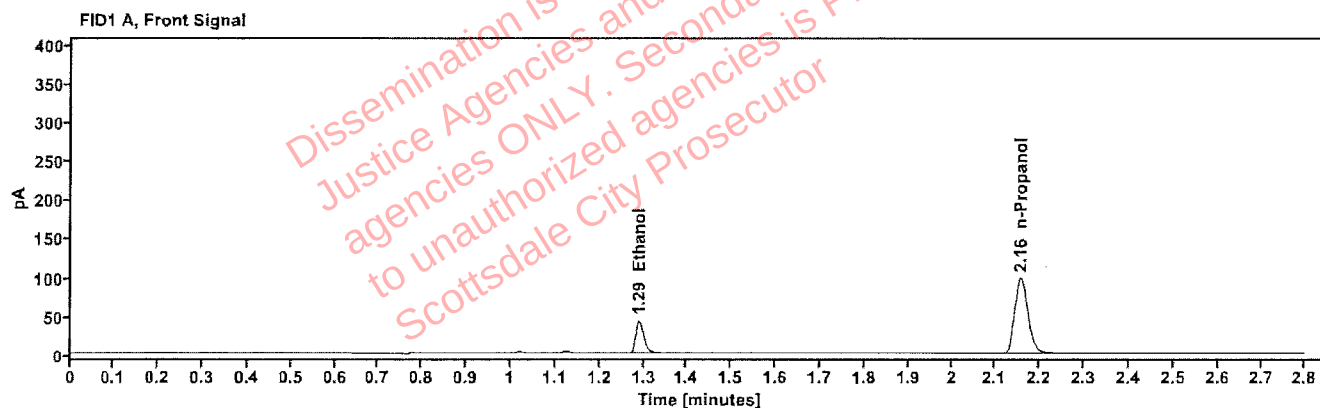


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0748	1.290	54.693
n-Propanol	-----	2.156	202.280

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.433	55.670
n-Propanol	2.611	207.003

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615387-1	Vial:	36
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:52:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\36.D	Analyst:	Abbott

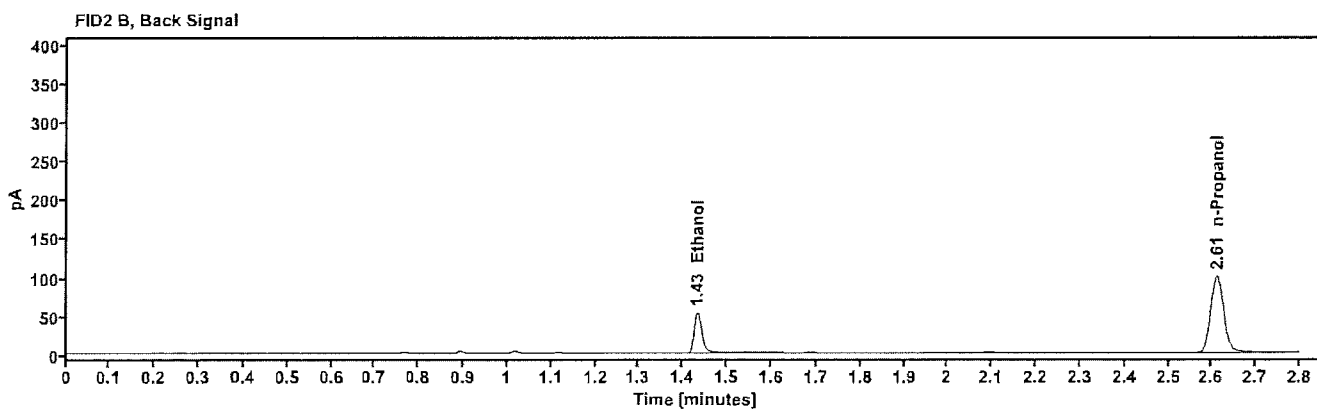
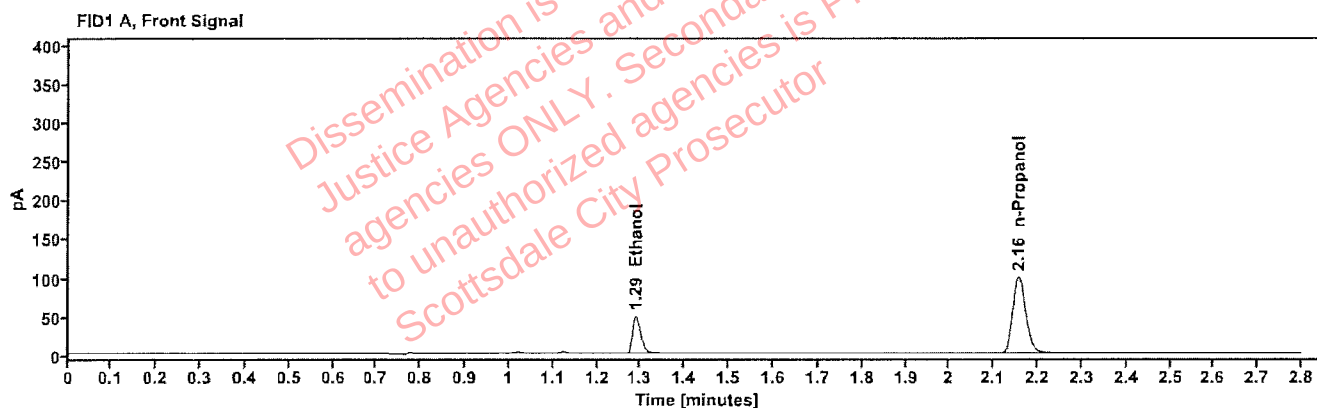


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0853	1.289	63.207
n-Propanol	---	2.156	204.498

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.433	64.676
n-Propanol	2.611	209.792

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615387-1	Vial:	37
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 4:57:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\37.D	Analyst:	Abbott

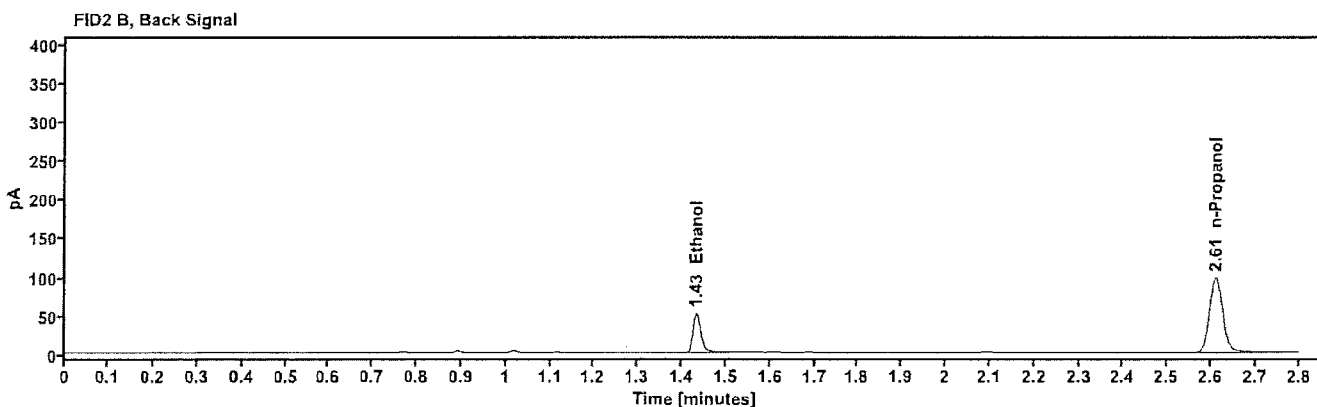
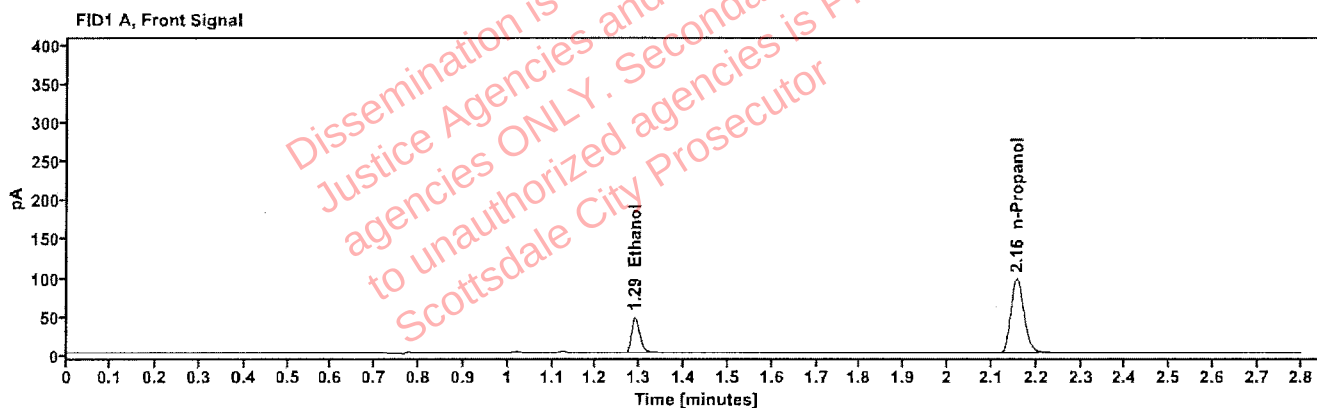


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0845	1.289	61.097
n-Propanol	-----	2.156	199.414

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.433	62.466
n-Propanol	2.611	204.556

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615382-2	Vial:	38
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:00:56 PM
Data file:	C:\Chem32\1\Data\06Aug26\38.D	Analyst:	Abbott

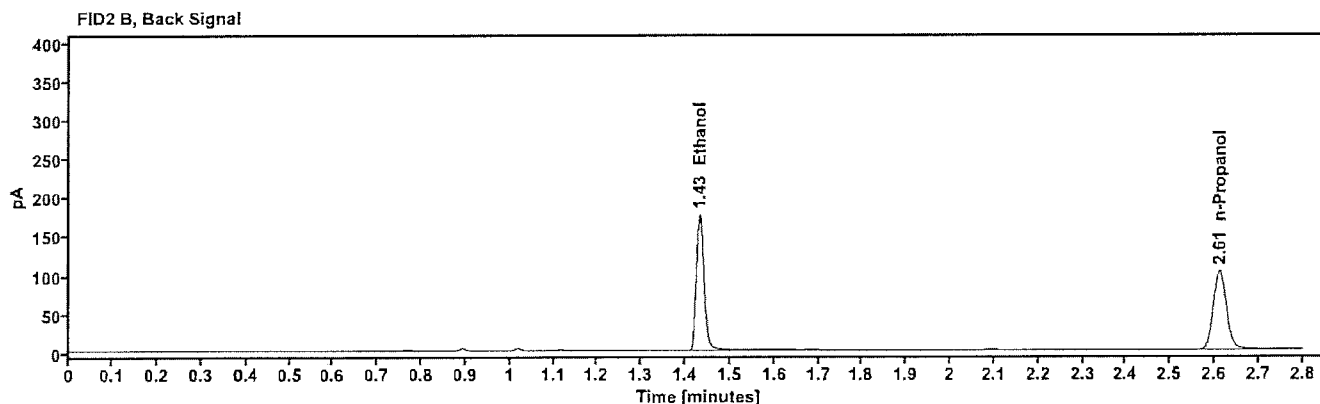
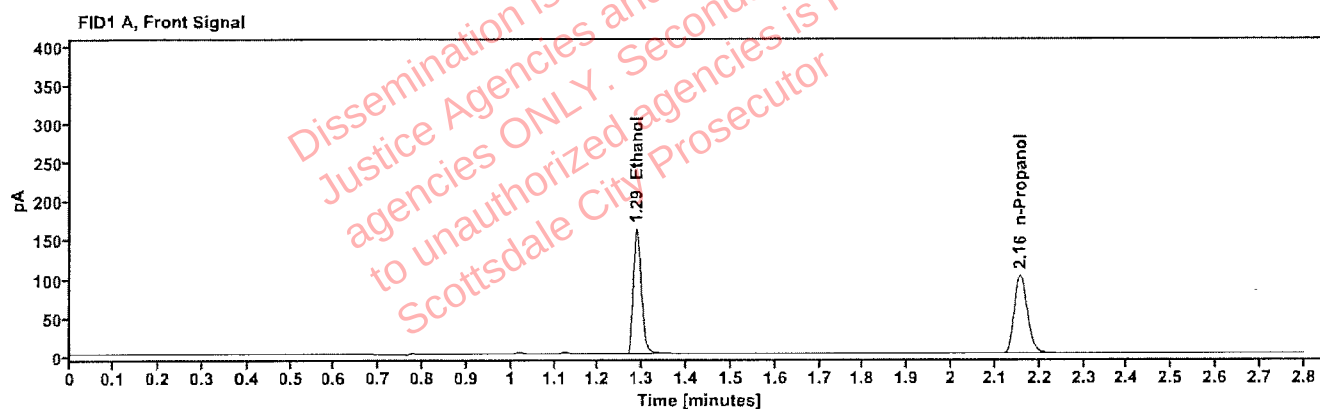


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2685	1.287	206.554
n-Propanol	-----	2.156	209.391

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	213.041
n-Propanol	2.611	214.582

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615382-2	Vial:	39
Description:	[REDACTED]	LIMS ID:	[REDACTED]
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:04:44 PM
Data file:	C:\Chem32\11\Data\06Aug26\39.D	Analyst:	Abbott

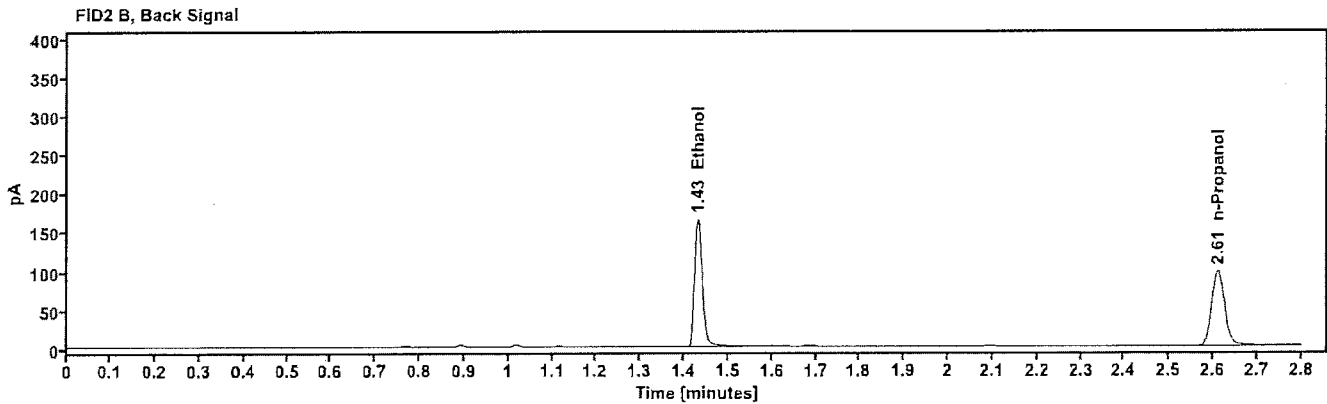
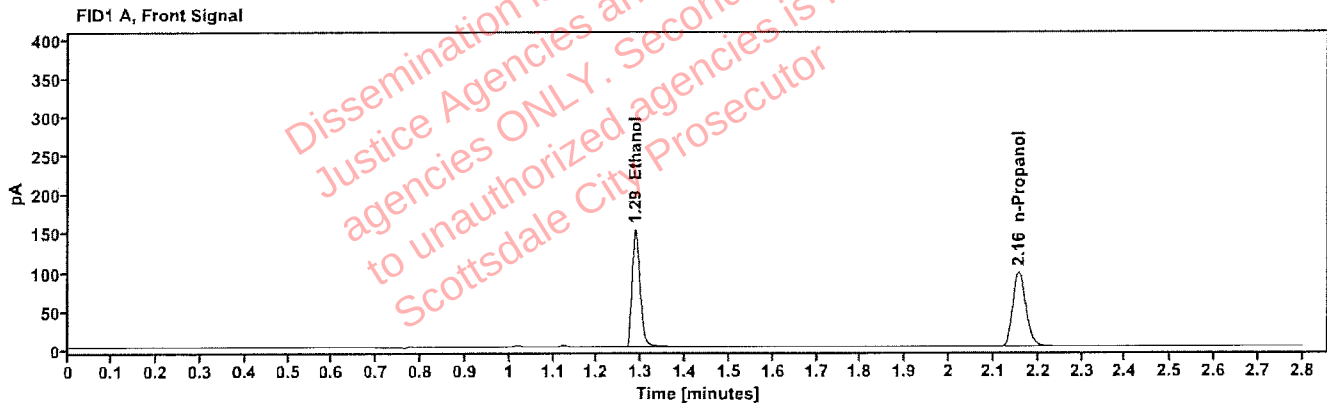


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2656	1.287	195.172
n-Propanol	-----	2.156	200.061

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.431	201.073
n-Propanol	2.611	204.910

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615606-1	Vial:	40
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:08:31 PM
Data file:	C:\Chem32\1\Data\06Aug26\40.D	Analyst:	Abbott

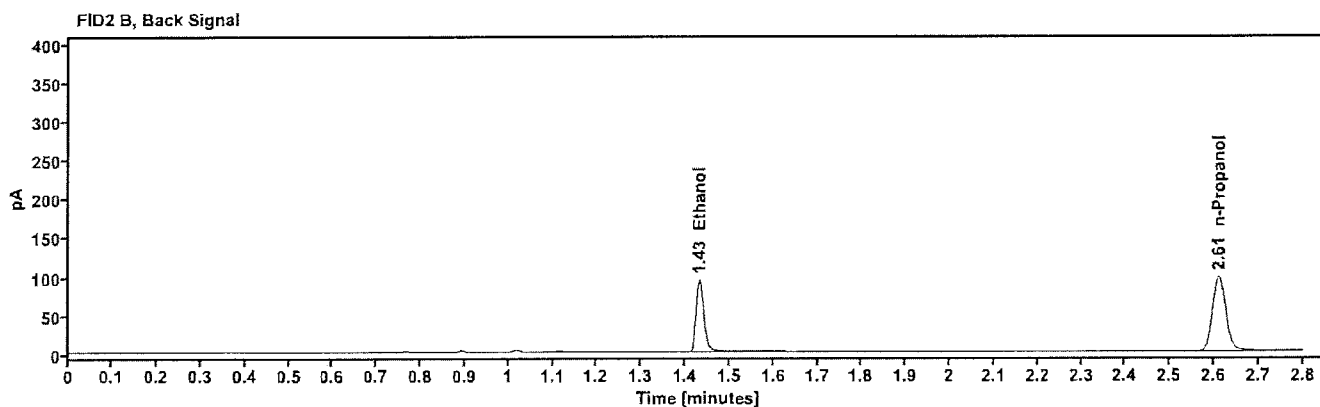
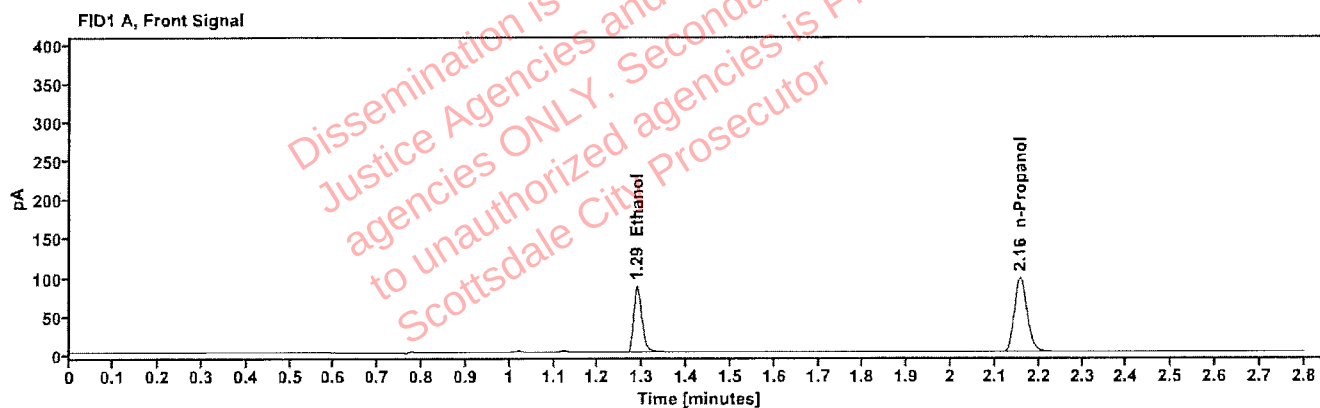


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1525	1.288	111.060
n-Propanol	-----	2.156	199.147

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	114.507
n-Propanol	2.611	203.941

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615606-1	Vial:	41
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:12:20 PM
Data file:	C:\Chem32\1\Data\06Aug26\41.D	Analyst:	Abbott

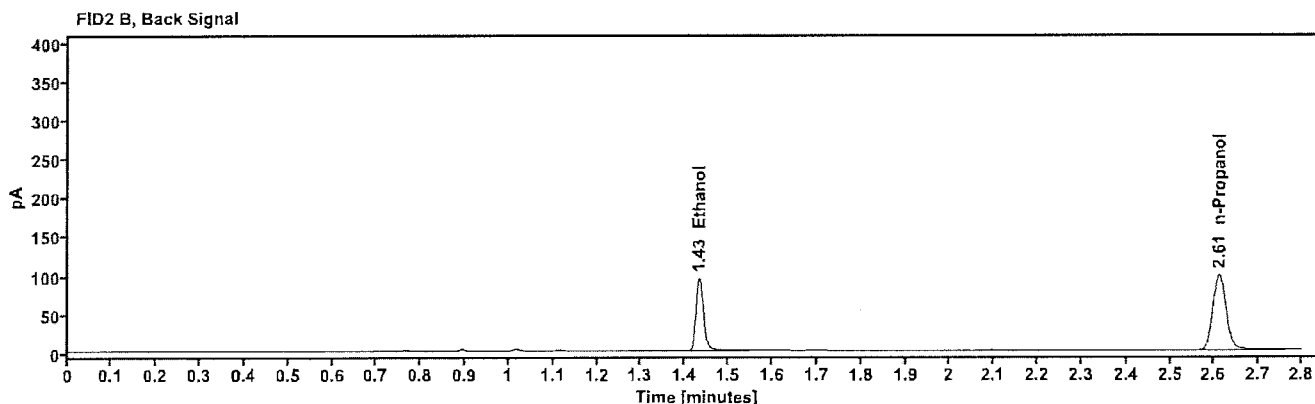
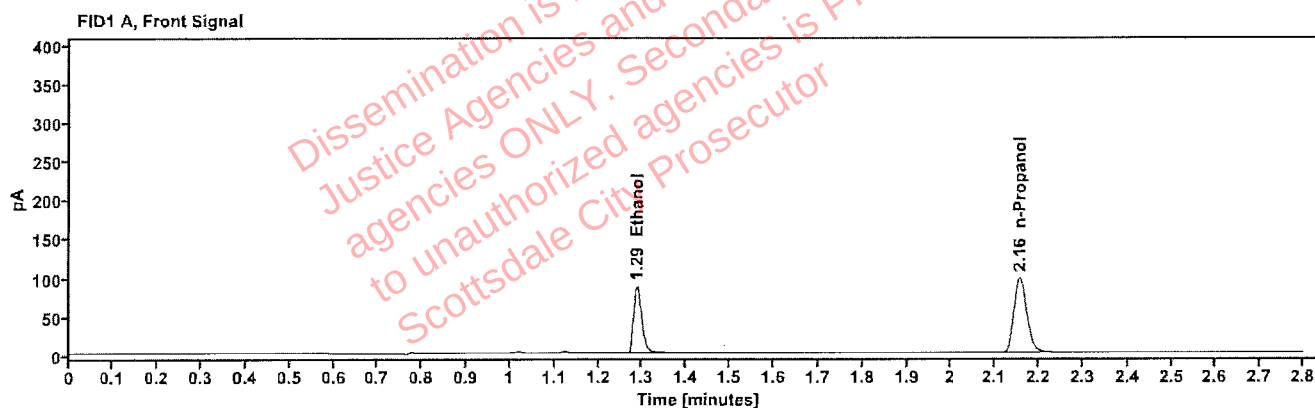


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1521	1.288	111.989
n-Propanol	-----	2.156	201.332

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	115.644
n-Propanol	2.611	206.002

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615599-1	Vial:	43
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:21:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\43.D	Analyst:	Abbott

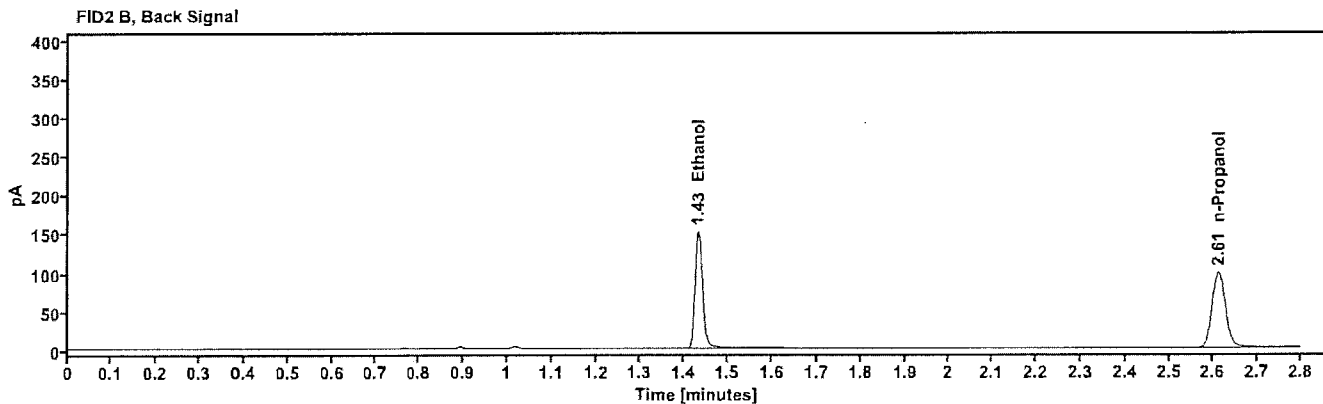
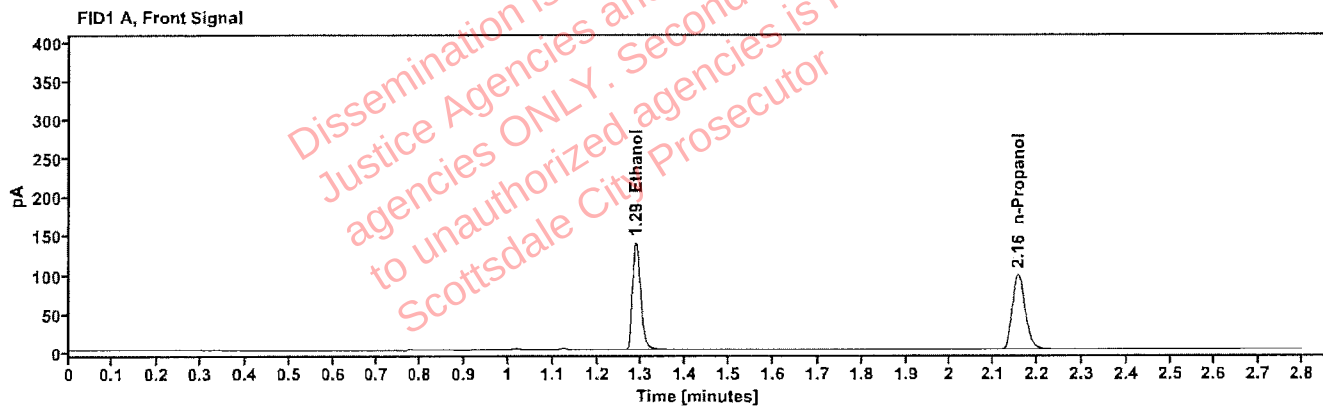


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2413	1.287	178.655
n-Propanol	----	2.156	201.716

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	184.255
n-Propanol	2.611	206.673

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615599-1	Vial:	44
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:24:56 PM
Data file:	C:\Chem32\1\Data\06Aug26\44.D	Analyst:	Abbott

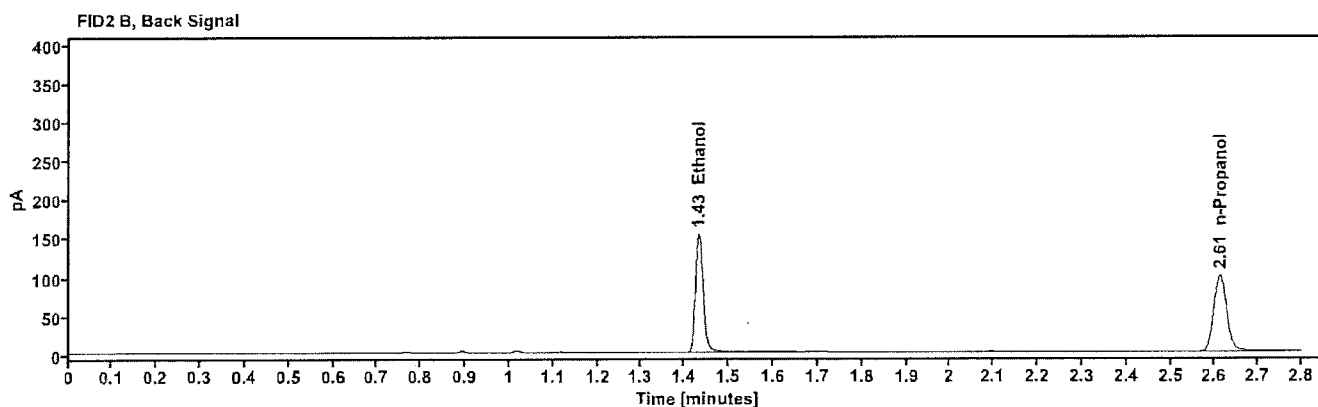
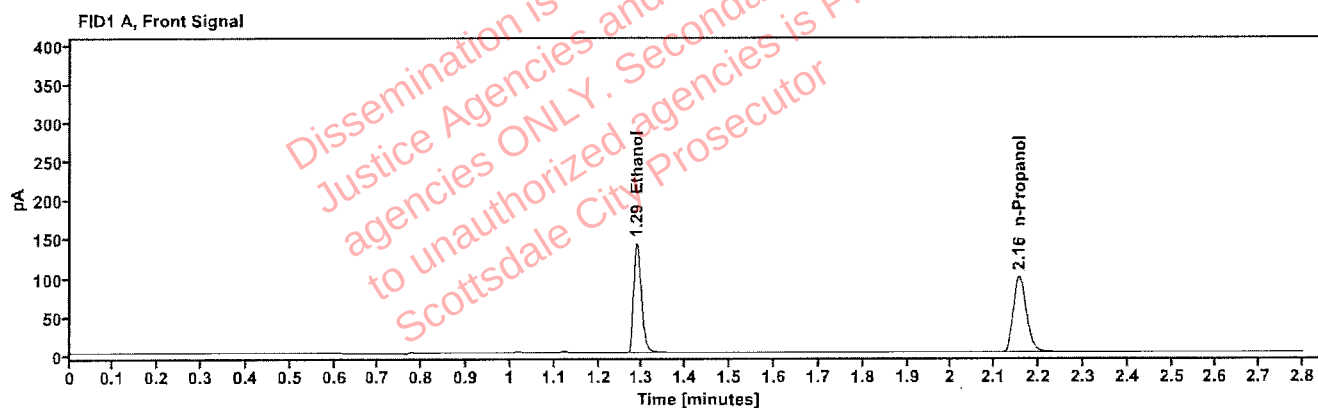


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2426	1.287	181.120
n-Propanol	---	2.156	203.397

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	186.948
n-Propanol	2.611	208.138

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615540-1	Vial:	45
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:28:44 PM
Data file:	C:\Chem32\1\Data\06Aug26\45.D	Analyst:	Abbott

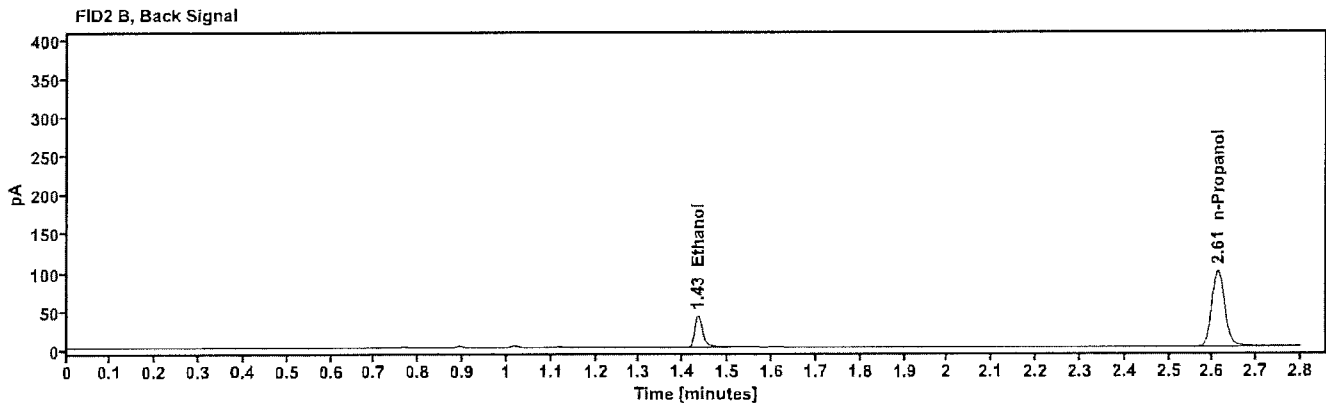
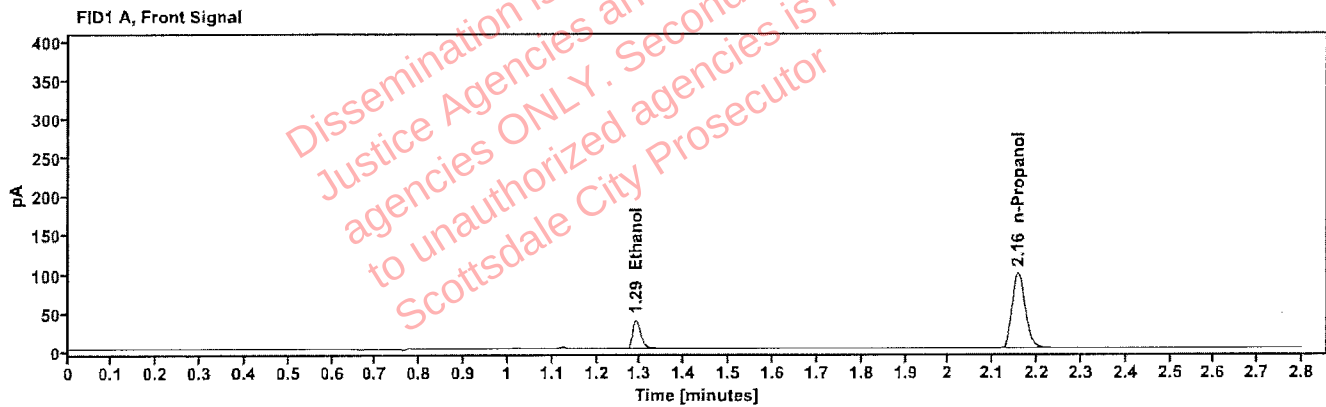


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0674	1.290	49.277
n-Propanol	-----	2.156	202.819

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.434	50.156
n-Propanol	2.611	207.435

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615540-1	Vial:	46
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:32:31 PM
Data file:	C:\Chem32\1\Data\06Aug26\46.D	Analyst:	Abbott

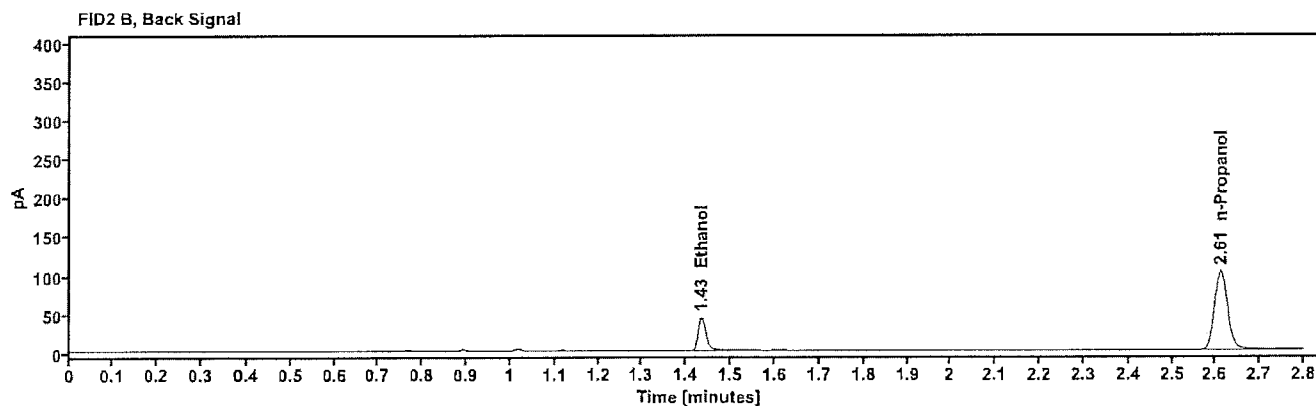
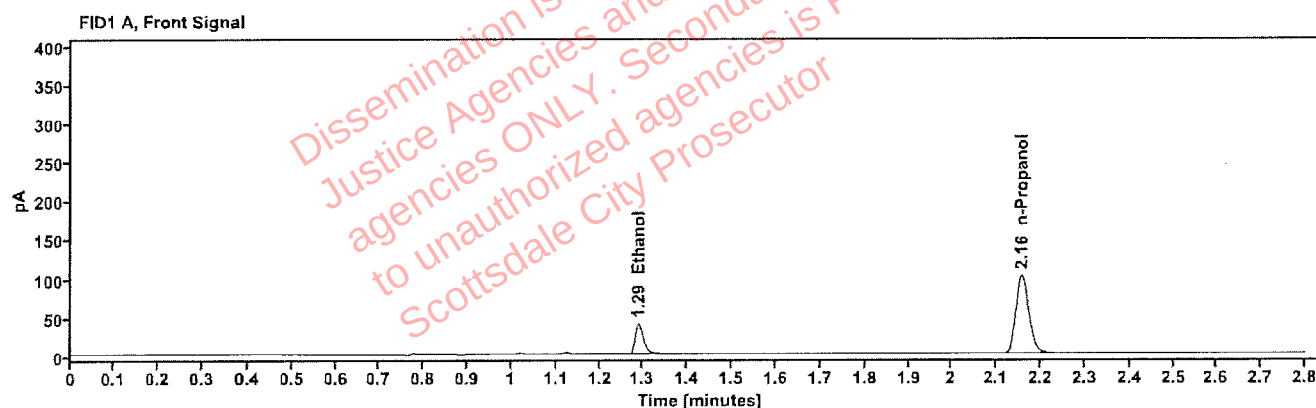


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.0682	1.290	51.673
n-Propanol	-----	2.156	209.922

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.434	52.512
n-Propanol	2.611	214.339

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615059-5	Vial:	47
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:36:19 PM
Data file:	C:\Chem32\1\Data\06Aug26\47.D	Analyst:	Abbott

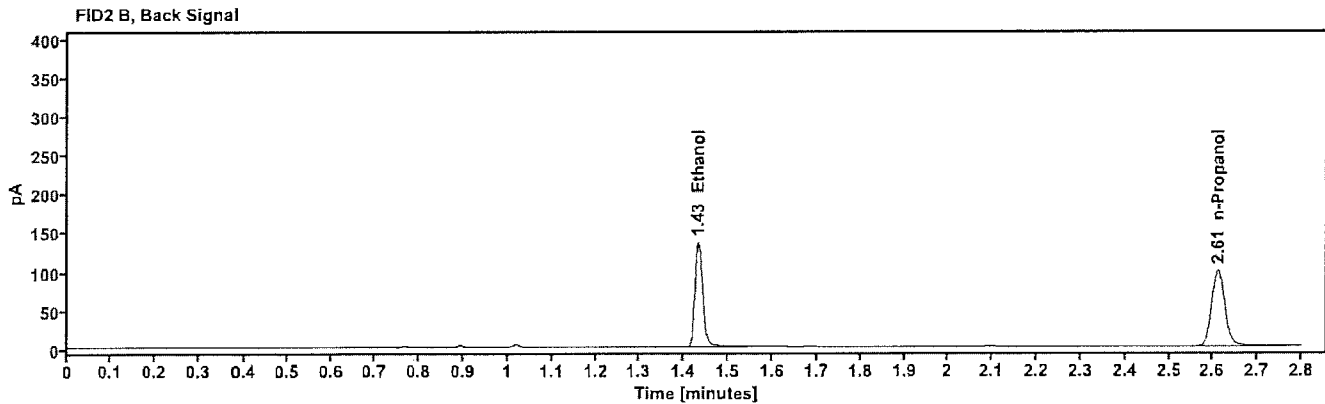
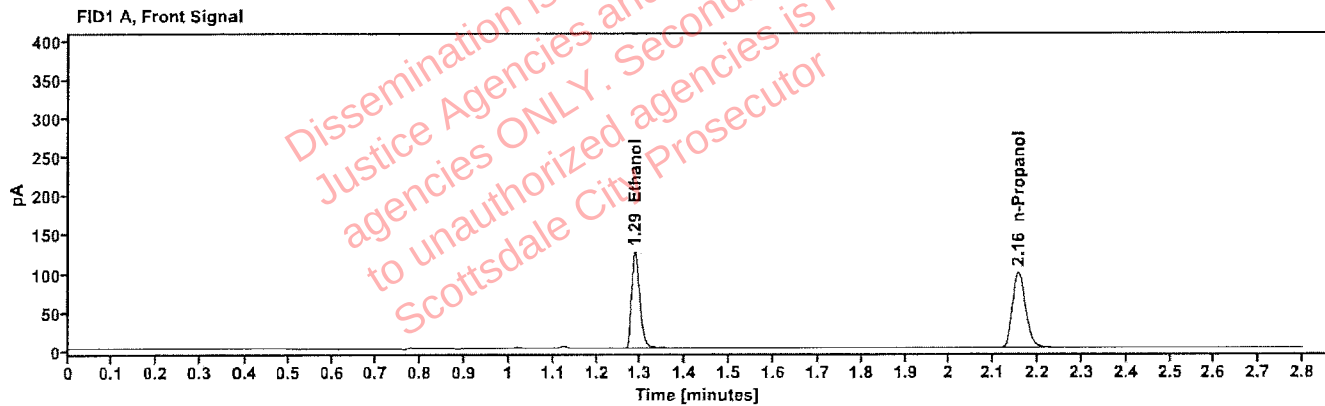


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2154	1.287	160.274
n-Propanol	-----	2.156	202.890

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	165.181
n-Propanol	2.611	208.074

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615059-5	Vial:	48
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:40:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\48.D	Analyst:	Abbott

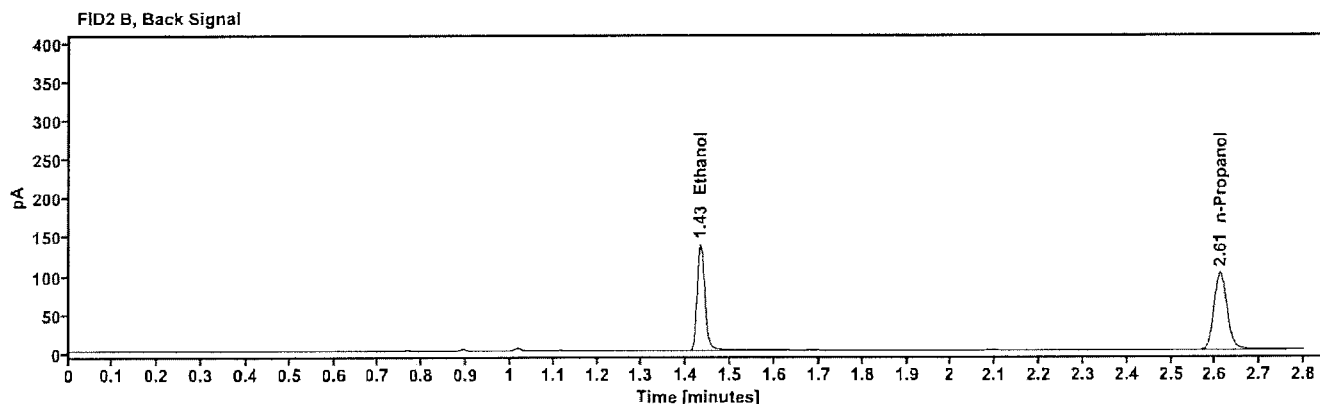
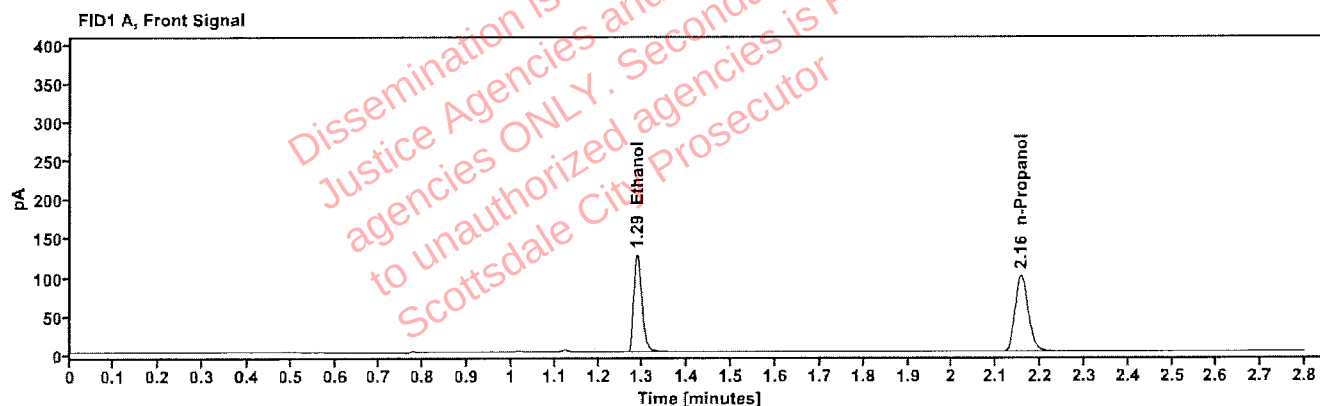


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.2159	1.287	162.222
n-Propanol		2.156	204.855

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	167.226
n-Propanol	2.611	209.786

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615520-2	Vial:	49
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:45:08 PM
Data file:	C:\Chem32\1\Data\06Aug26\49.D	Analyst:	Abbott

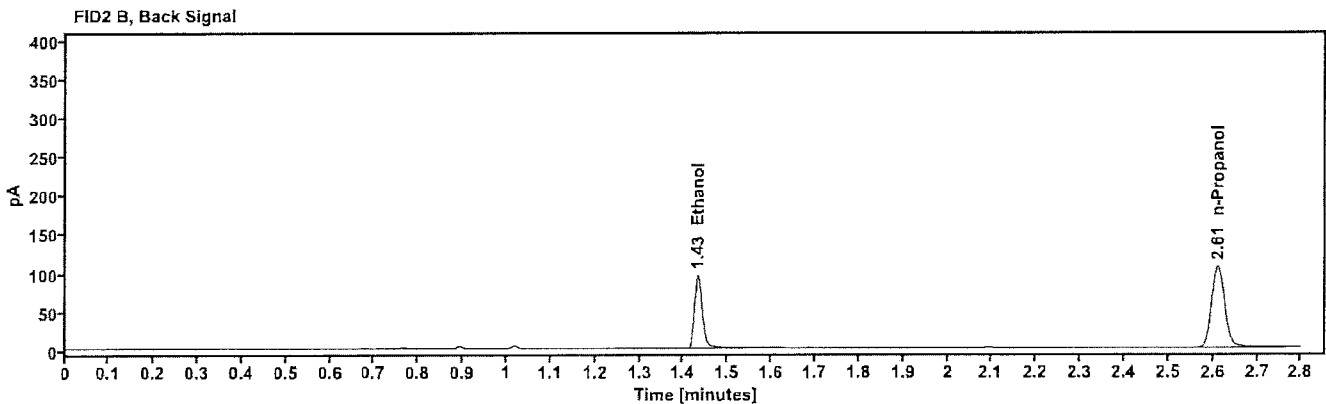
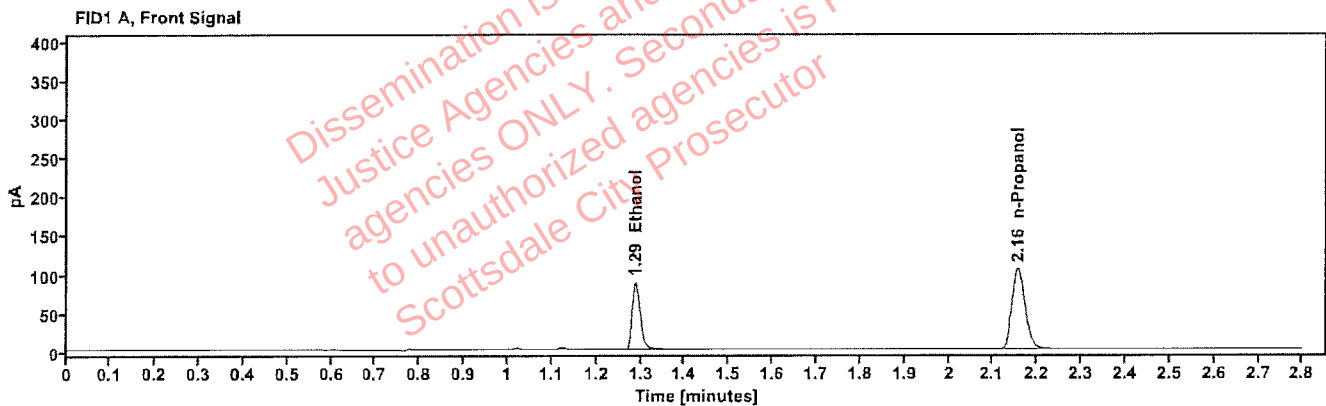


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1426	1.288	112.711
n-Propanol	-----	2.155	216.320

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	116.103
n-Propanol	2.611	221.743

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615520-2	Vial:	50
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:48:55 PM
Data file:	C:\Chem32\1\Data\06Aug26\50.D	Analyst:	Abbott

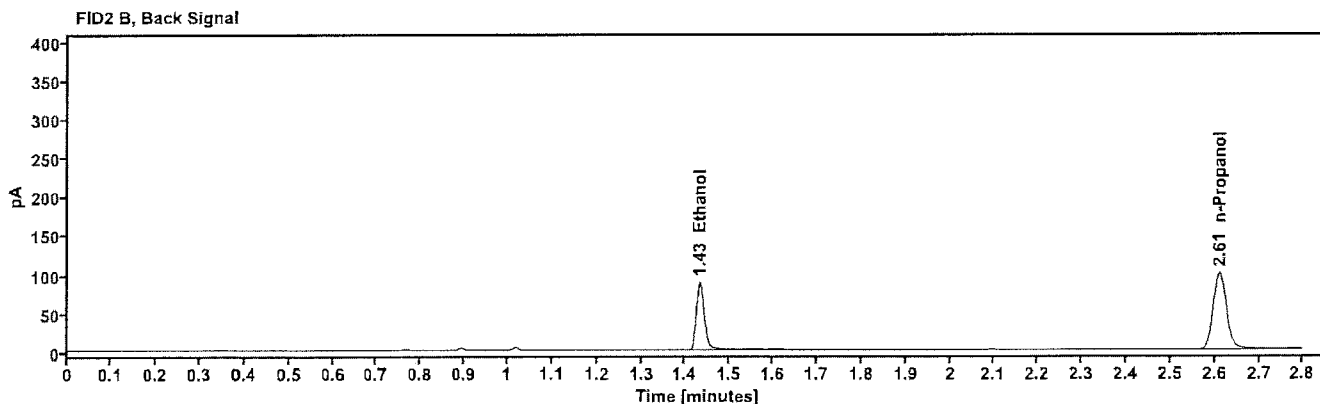
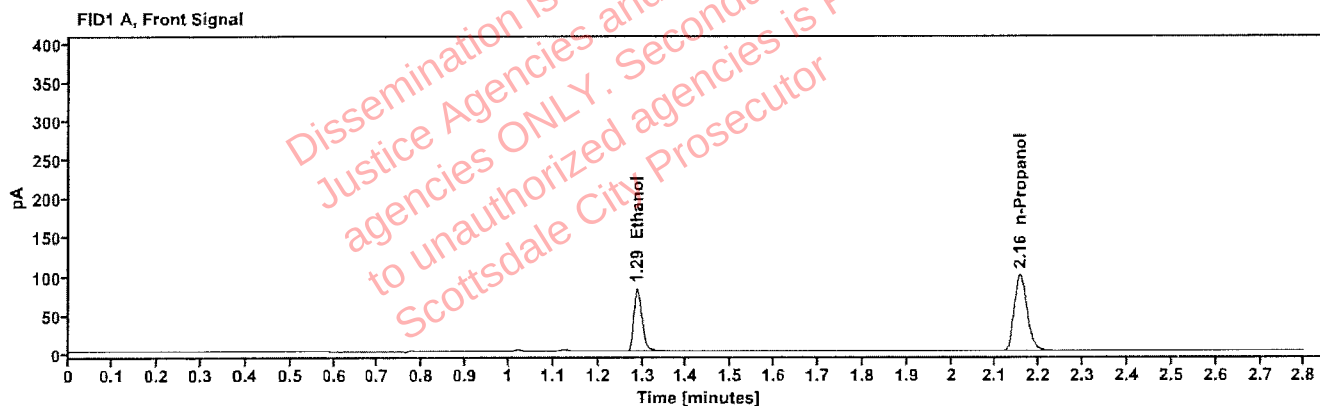


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1406	1.288	105.023
n-Propanol	---	2.155	204.507

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	108.326
n-Propanol	2.611	209.477

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615522-1	Vial:	51
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:52:44 PM
Data file:	C:\Chem32\1\Data\06Aug26\51.D	Analyst:	Abbott

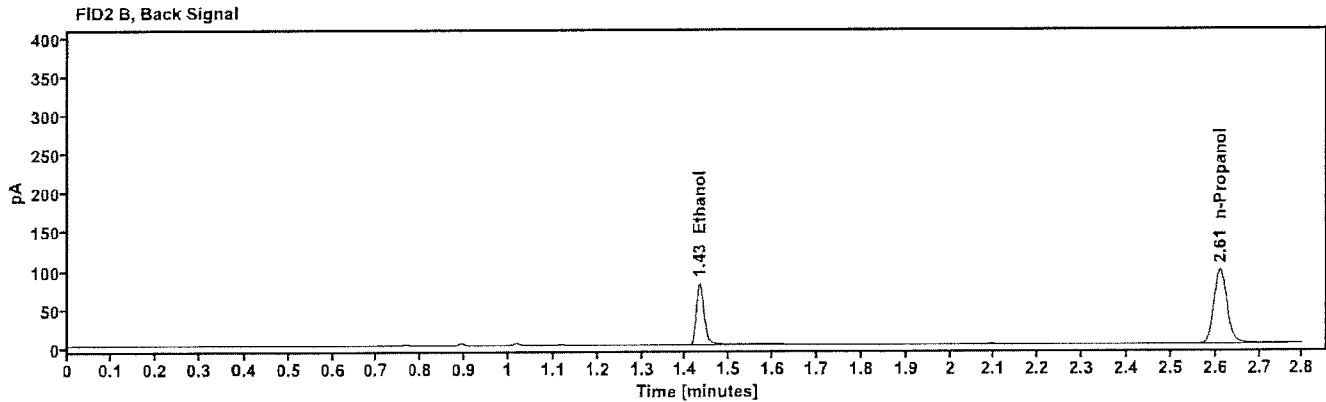
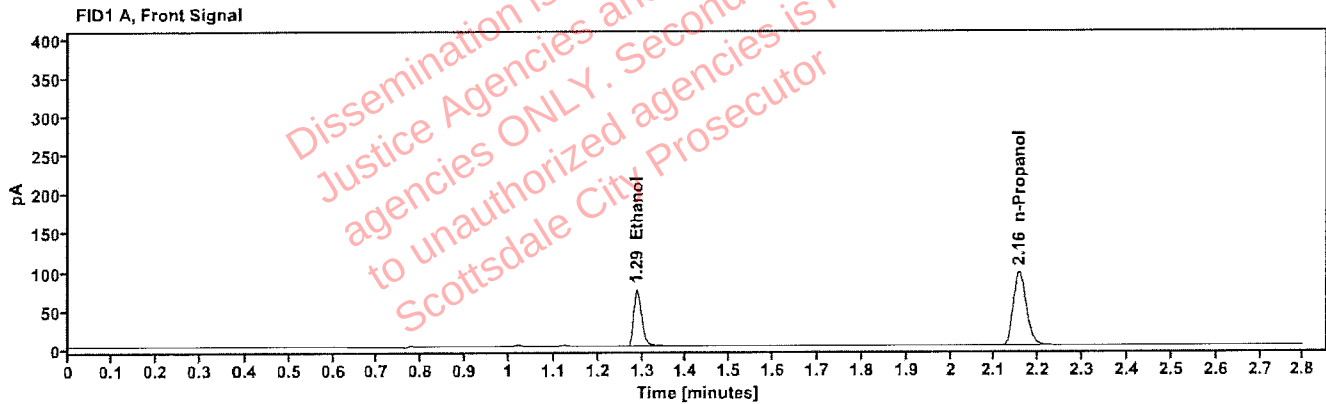


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1322	1.288	95.488
n-Propanol	---	2.155	197.839

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	98.489
n-Propanol	2.610	202.441

Scottsdale Police Department Crime Lab Volatiles Analysis

Sample:	2615522-1	Vial:	52
Description:		LIMS ID:	
Method:	ethanol quant.M	Sequence:	06Aug26
Instrument:	US2320A024 US23180007	Injection date:	8/6/2026 5:56:32 PM
Data file:	C:\Chem32\1\Data\06Aug26\52.D	Analyst:	Abbott

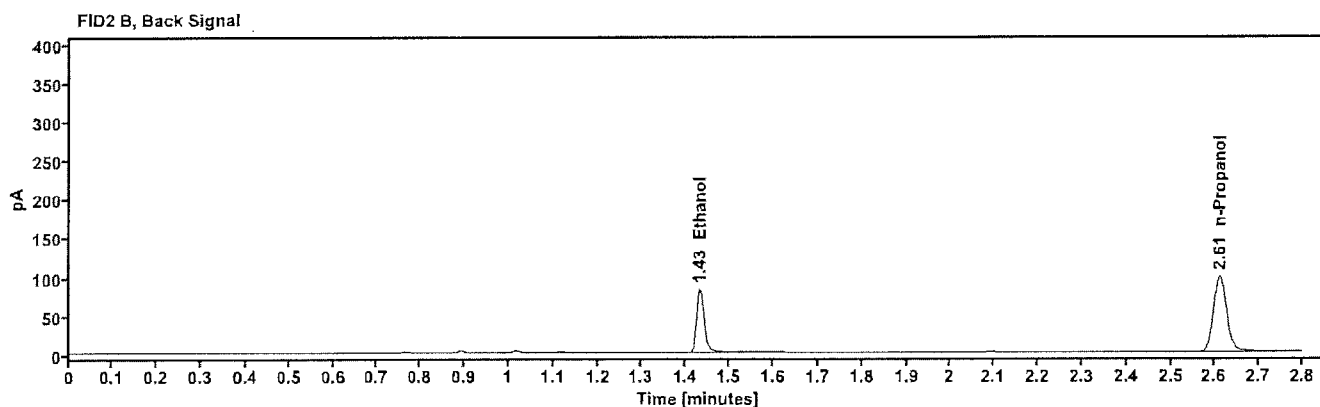
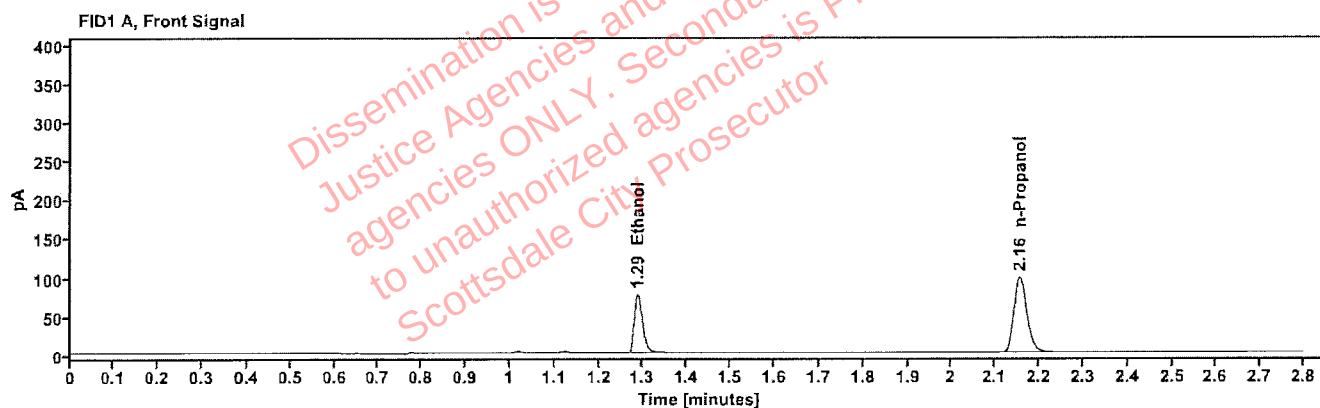


Table 1: FID 1 A (column DB-ALC1)

Compound	Amount (g/100mL)	Time (min)	Peak Area
quantified Ethanol	0.1331	1.288	97.928
n-Propanol	—	2.155	201.584

Table 2: FID 2 B (column DB-ALC2)

Compound	Time (min)	Peak Area
Ethanol	1.432	101.013
n-Propanol	2.611	206.437